

Example of Use of Alternate Carrier Gas in GC Analysis: Change from He Gas to N₂ Gas

Helium (He) has mainly been used as the carrier gas for gas chromatography (GC), but limited availability and a sharp rise in the price of He have become issues. In recent years, hydrogen (H₂) and nitrogen (N₂) have been actively used as alternate carrier gases for He. As an example of an analysis in which the carrier gas was changed from He to N₂, this article introduces an example of an analysis by the vegetable oil (olive oil) immersion method for plastic materials using GC-FID. In this analysis, we used the overall migration test referred to the European Committee for Standardization, EN1186 "Materials and articles in contact with foodstuffs - Plastics (No 10/2011)"

T. Ishii

Points to Note when Changing Carrier Gases

Fig. 1 shows an example of the characteristics of HETP (height equivalent of one theoretical plate), which shows the separation efficiency of a column when using He, N₂, or H₂ as the carrier gas.

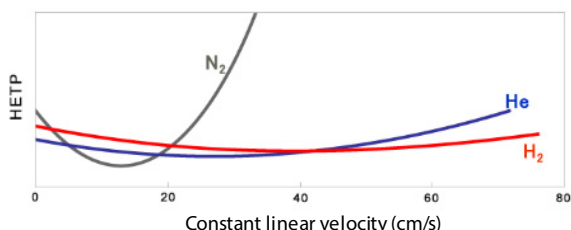


Fig. 1 Example of Relationship of Constant Linear Velocity and HETP of Carrier Gases

Although the HETP value changes depending on the constant linear velocity, in order to maximize the separation performance of column, it is important to conduct the analysis at the constant linear velocity which minimizes HETP as far as possible. For example, the optimum condition with He is in the vicinity of 30 to 40 cm/s. On the other hand, when analyzed using N₂ under He - optimized conditions, separation performance of the column will decrease due to the difference in HETP. The left side of Fig. 2 shows part of a chromatogram when using He, but as shown on the right, when the same constant linear velocity is used with N₂, it is no longer possible to separate C 18:0 and C 18:1.

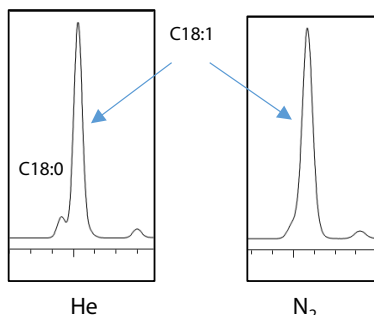


Fig. 2 Difference of Separation with He and N₂ at Linear Velocity of 40 cm/s

Calculation of Optimum Linear Velocity by EZGC®

In this experiment, the optimum conditions when using different carrier gases were studied using the method conversion program EZGC® Method Translator, which is available from Restek Corporation via the internet. Fig. 3 shows the screen of EZGC® Method Translator. When the conditions for He are input in the "Original" column and the alternate carrier gas type is selected under "Translation" column, the optimum parameters for the alternate carrier gas are displayed. It is also possible to select the priority of the analysis by using the "Result" item. For example, select "Speed" when you wish to shorten the analysis time, and select "Efficiency" when separation efficiency is the priority. Because the focus in this experiment was separation of C 18:0 and C 18:1, the method was translated by selecting "Efficiency," which prioritizes separation efficiency, and the program calculated the optimum constant linear velocity of 18.11 cm/s.

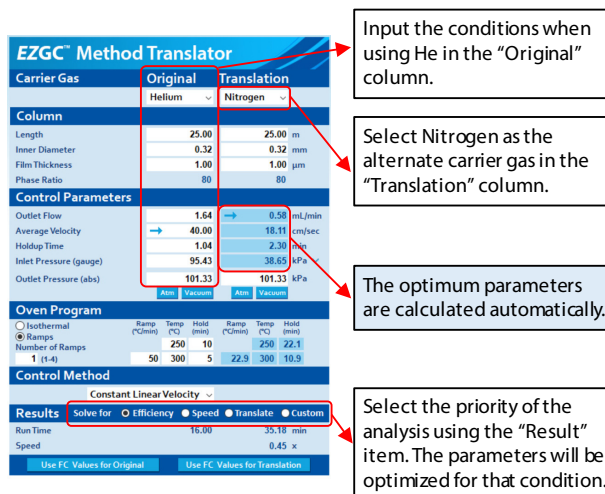


Fig. 3 EZGC® Method Translator Screen

Use of Gas Selector

Gas selector (P/N: S221-84916-41), which is a dedicated option of Nexis™ GC-2030, not only changes the gas flow lines, but also makes it possible to switch between two types of carrier gas from the software. Use of the gas selector in combination with EZGC® enables easy calculation and control of N₂ for low-cost analyses which do not require high resolution, and He or H₂ when high resolution, high speed analysis is necessary. For details concerning the gas selector, please refer to Application News No. G328.

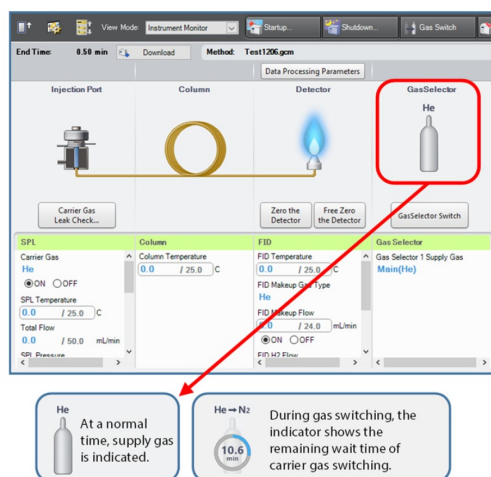


Fig. 4 Screen of Carrier Gas Switching Function of Gas Selector

Example of Analysis Using Optimal Linear Velocity for Carrier Gases

When an analysis was carried out with N₂ as the carrier gas using the linear velocity calculated by EZGC[®] (Table 1), the chromatogram shown in Fig. 5 was obtained. Although the analysis time was long in comparison with an analysis using He carrier gas (Fig. 6), an evaluation of the resolution R of C18:0 and C18:1 confirmed that resolution was almost the same, as R = 0.93 was obtained with the nitrogen carrier and R = 0.90 was obtained with the helium carrier.

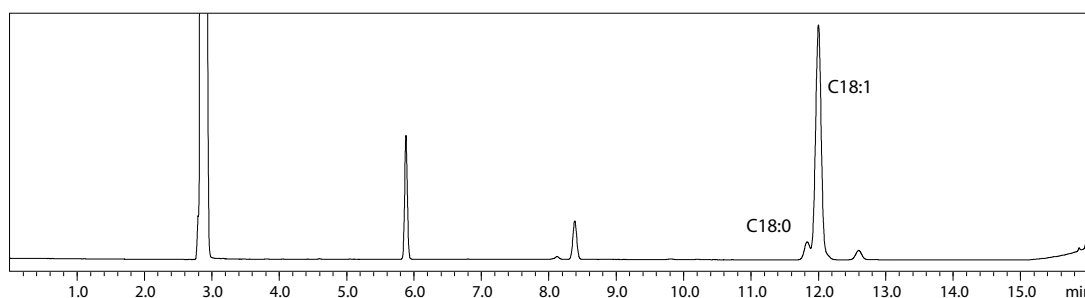


Fig. 5 Chromatogram when Using N₂

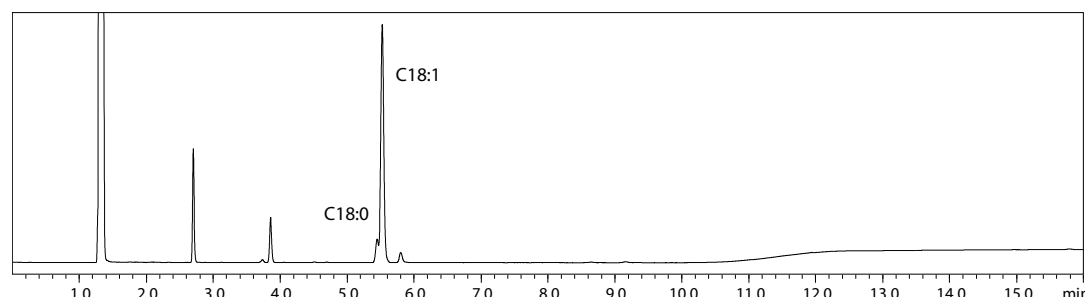


Fig. 6 Chromatogram when Using He

Acknowledgement: The data introduced in this Application News article were provided by SGS Japan Inc.

Nexis and LabSolutions are trademarks of Shimadzu Corporation in Japan and/or other countries. EZGC is either a trademark or a registered trademark of Restek Corporation in the United States and/or other countries. To use EZGC[®], please access the program at the Restek website. If the program does not run on your browser, please use the downloadable application at the bottom of the screen.

First Edition: Nov. 2020



Shimadzu Corporation

www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. Shimadzu disclaims any proprietary interest in trademarks and trade names used in this publication other than its own. See <http://www.shimadzu.com/about/trademarks/index.html> for details.

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

© Shimadzu Corporation, 2020

Table 1 Analysis Conditions

Injection	: Split (1 : 50), 320 °C
Carrier gas	: He, N ₂
Carrier gas control	: Column linear velocity He 40.0 cm/s N ₂ 18.1 cm/s (calculated by EZGC [®])
Column	: BP-1 (25 m × 0.32 mm I.D., 1.0 μm)
Column temp.	: He 250 °C (10 min) – 50 °C/min – 300 °C (5 min) N ₂ 250 °C (15 min) – 50 °C/min – 300 °C (5 min)
Detector	: Flame ionization detector (FID), 320 °C

* Calculation of the optimal linear velocity by EZGC[®] is applicable to the GC-2010 series, GC-2025, and GC-2014 (capillary models) in addition to the Nexis™ GC-2030.

* The optional gas selector can only be used with the Nexis™ GC-2030.

Conclusion

An analysis was conducted using the method conversion program EZGC[®] Method Translator, assuming that the carrier gas is changed from He to N₂. This experiment confirmed that the optimum linear velocity condition for N₂ carrier gas can be obtained easily by only inputting the conditions for He currently being used.

If the optional gas separator is used, simple switching of carrier gases is possible in LabSolutions™. Automatic switching of carrier gases is also possible during continuous analysis, and multiple analyses with different separation requirements using the optimal carrier gas for each analysis can also be conducted automatically in a single sequence.