

Application News

No. Q122

Powder Property Analysis

Characterization of Particles in Suspension Eye Drops: Detection of Maximum Length of Irregular Particles

Eye drops can be classified as aqueous (nonaqueous) eye drops in which the ingredients are dissolved in water (oil), and aqueous (nonaqueous) suspension eye drops in which the ingredients are not dissolved but the particles are suspended. The Japanese Pharmacopoeia states that the maximum particle size of particles in suspension eye drops is generally 75 μm or less.

The laser diffraction method is widely used as a technique for measurement of particle size distribution because measurement time is short and the measurement range is wide. However, this technique also has certain problems in cases where it is necessary to obtain the maximum length, as it is difficult to detect coarse particles that exist in very small quantities relative to the total amount of particles, and it is not possible to measure the maximum length of aspherical particles because the particle size is calculated as the sphere equivalent diameter.

This article introduces an example of characterization of the maximum length of particles in a suspension eye drop product and a mixed sample containing both spherical and acicular particles by acquiring particle images and analyzing the particle shape, size distribution, and concentration with the iSpect™ DIA-10 dynamic particle image analysis system (Fig. 1).

H. Maeda

Samples


Two types of samples were analyzed, the first being a suspension eye drop product which is available in the market, and the second, a dispersion containing a mixed sample of polystyrene particles and glass fibers as an example of an admixture containing aspherical coarse particles. The dispersion medium of both samples was water.

Measurement Results


Fig. 2 shows the measurement results for the suspension eye drop product. The thumbnail image shows that the particle shape is nearly circular, but each particle has a different shape. From the particle size distribution, no particles with sizes of approximately 15 μm or larger were detected. The concentration of particles with sizes of 5 μm or larger was 2,652 particles/mL.



Fig. 1 iSpect™ DIA-10 Dynamic Particle Image Analysis System

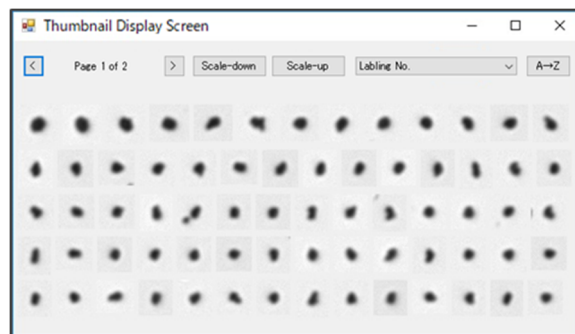

a) Thumbnail images

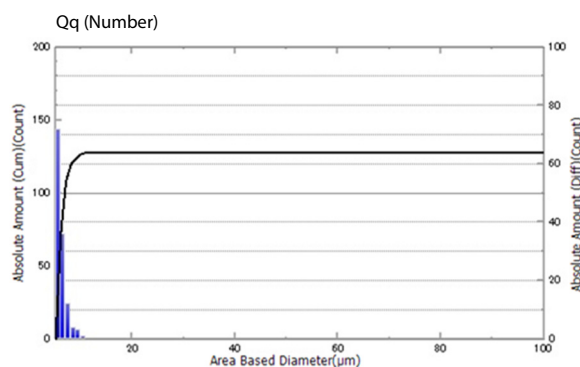

b) Particle size distribution

Fig. 2 Measurement Results of Suspension Eye Drop Product

Some particle size distribution analysis systems can obtain the measurement results of the particle diameter in the form of the equivalent diameter of a sphere of the same volume (volume based diameter) or a circle of the same area (area based diameter). When considering the maximum length, values near the area based diameter and maximum length can be obtained if the particle shape is spherical. However, the maximum length of acicular (needle-shaped) particles and other non-spherical shapes is generally larger than the area based diameter. Fig. 3 shows the definitions of particle size.

As a concrete example, Fig. 4 shows the results of measurement of a mixed sample containing polystyrene particles and glass fibers. From the thumbnail images, it can be understood that this is a mixture of spherical particles and acicular particles.

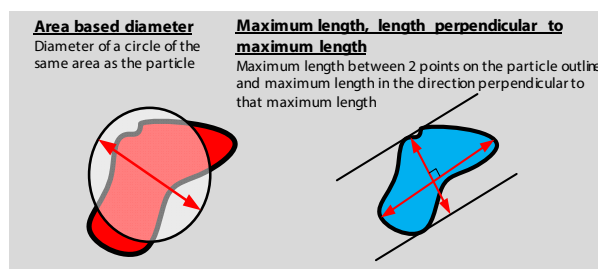


Fig. 3 Definitions of Particle Size

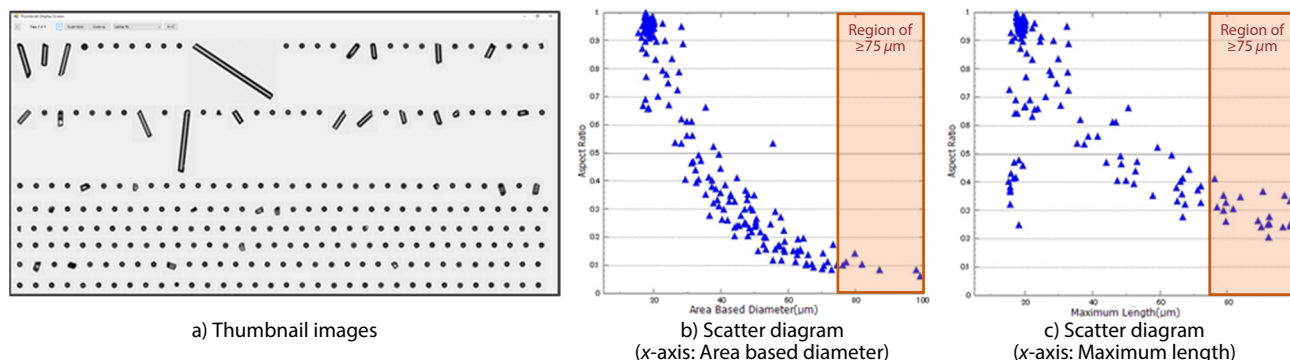


Fig. 4 Measurement Results of Mixed Sample of Polystyrene Particles and Glass Fibers

For example, the hatched regions on the right side of the scatter diagrams in Fig. 4b and c are the regions of particle sizes of 75 μm and larger. The concentration of particles included in this region is 47 particles/mL when calculated by the area based diameter (Fig. 4b), but 420 particles/mL when calculated by the maximum length (Fig. 4c). Thus, in comparison with the maximum length method, the number of particles with sizes of 75 μm and larger is smaller when the area based diameter method is used. Based on this, calculation of the maximum length from particle images, and not the area based diameter, is effective when it is necessary to evaluate the maximum length.

As shown in Fig. 5, classification of the particles in mixtures of particles with different shape parameters (in Fig. 5, Aspect ratio = Length perpendicular to maximum length ÷ Maximum length) is possible by using those parameters, and it is also possible to obtain the ratios of the respective particles. (Here, the relative contents are spherical particles 69.4% and acicular particles 30.6%.)

Conclusion

In an analysis of a suspension eye drop product with an iSpect DIA-10, it was possible to confirm that all detected particles had sizes of less than 75 μm (actual result with this sample <15 μm). In analyses of mixed samples containing particles with different shape parameters, such as polystyrene particles and glass fibers, it is also possible to obtain the ratios of the respective particles by classifying the particles by their shape parameter. As demonstrated in these experiments, characterization of the maximum particle length and classification of particles in particle groups that had been difficult to characterize by the laser diffraction method is possible by using the iSpect DIA-10.

iSpect is a trademark of Shimadzu Corporation in Japan and/or other countries.

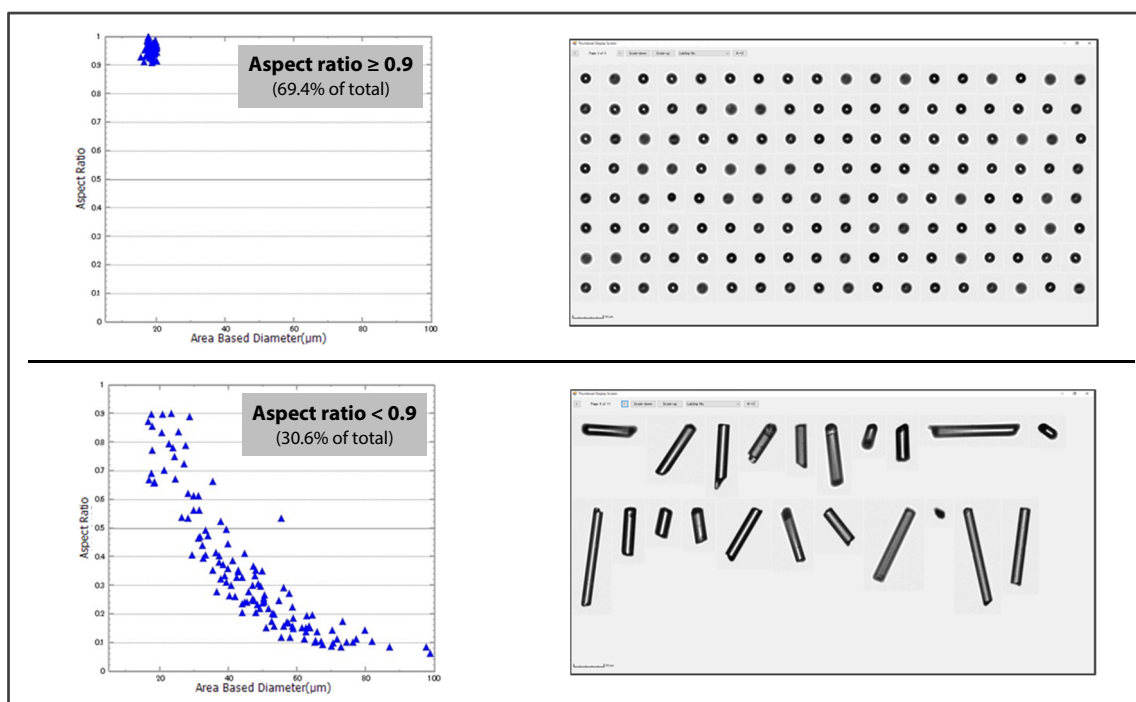


Fig. 5 Example of Particle Classification

First Edition: Dec. 2019



Shimadzu Corporation

www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. Shimadzu disclaims any proprietary interest in trademarks and trade names used in this publication other than its own. See <http://www.shimadzu.com/about/trademarks/index.html> for details.

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

© Shimadzu Corporation, 2019