

Supercritical Fluid Chromatograph

# SFC Basic Guide





CO<sub>2</sub> cylinders used for supercritical fluid chromatography (SFC) are filled with supercritical carbon dioxide generated by industrial plants that manufacture ammonia, refine petroleum, and so on. In other words, no new carbon dioxide is generated when operating SFC. On the other hand, we generate new carbon dioxide in the environment while treating the organic solvent waste from HPLC. Therefore, switching the separation method from HPLC to SFC reduces the quantity of CO<sub>2</sub> greenhouse gas emitted, achieving one of the sustainable development goals (SDGs) proposed by the United Nations, "13. Take urgent action to combat climate change and its impacts."

This brochure provides a collection of basic information about SFC.

# I N D E X

<b>Chapter 1 Supercritical Fluid Chromatography</b>	2
1-1 Properties of Supercritical Fluids	2
1-1-1 Physical Properties of Supercritical Fluids	2
1-1-2 Supercritical Carbon Dioxide	3
1-2 Supercritical Fluid Chromatography	4
1-2-1 Chromatography in General	4
1-2-2 Supercritical Fluid Chromatography	4
1-2-3 Instrumentation for Supercritical Fluid Chromatography	5
1-3 Analytical Advantages of SFC	6
1-3-1 Low Polarity — Better Structural Identification	7
1-3-2 Miscibility with Organic Solvents — Separations Covering Larger Chemical Space	7
1-3-3 Low Viscosity — Achieving Ultra-Fast Analysis	8
1-3-4 Vaporization — Improved LC-MS Detection Sensitivity	10
1-4 SFC at Preparative Range	11
<b>Chapter 2 Nexera UC Supercritical Fluid Chromatograph</b>	13
2-1 Components of SFC System	13
2-2 Nexera UC/s: SFC/UHPLC Switching System	14
2-3 Preparative SFC Systems	15
2-3-1 Stacked Fraction System	15
2-3-2 Multi-Fraction System	15
2-3-3 Analytical Fraction System	16
2-3-4 LotusStream Gas-Liquid Separator	16
2-4 Online SFE-SFC System	17
2-5 Offline SFE-SFC System	19
<b>Chapter 3 SFC Analytical Condition Parameters</b>	20
3-1 SFC Columns	21
3-1-1 Column Selection	23
3-1-2 Shim-pack UC Series	22
3-2 SFC Mobile Phases	24
3-2-1 Modifiers	24
3-2-2 Using Additives to Adjust Separation and/or Improve Peak Shape	25
3-2-3 Other Factors that Affect Separation	26
3-3 Support for SFC Method Development: Method Scouting Systems	27
<b>Chapter 4 Examples of SFC Solutions</b>	28
4-1 Analyzing Triglycerides by SFC	28
4-2 Separating Chiral API Molecule using SFC	29
4-3 Using Offline SFE to Quickly Extract Compounds from Natural Products	31
4-4 Fast Analysis of Ciclesonide Inhalant by Online SFE-SFC-QTOFMS	32
4-5 Preparative Purification of Volatile Aroma Components by SFC	33
<b>Chapter 5 Frequently Asked Questions about SFC</b>	34
5-1 SFC Installation	34
5-2 Samples and Target Components	35
5-3 Mobile Phase	35
5-4 Columns	36
5-5 Miscellaneous	36

# Chapter 1 Supercritical Fluid Chromatography

Supercritical fluid chromatography (SFC) is a separation method that uses supercritical carbon dioxide as the mobile phase. Though SFC technology has existed for a long time, recent advancements in SFC instrument and column technologies have made it a more common technique used for regular lab use and measurements.

By using supercritical carbon dioxide, which can be mixed with water-soluble solvents despite its low polarity, SFC offers potential advanced separation methods that are different from conventional GC and HPLC separation. This chapter describes what is supercritical carbon dioxide with an overview about new separation methods which are made possible using supercritical carbon dioxide.

## 1-1 Properties of Supercritical Fluids

### 1-1-1 Physical Properties of Supercritical Fluids

Substances can change their state between solid, liquid, and gas depending on the temperature and pressure. Vaporization and sublimation are phenomena that occur when heat activates the thermal motion of solid or liquid molecules in a substance, to an extent that the molecules break free from the intermolecular forces keeping the substance in its current state. For example, solid CO<sub>2</sub> (dry ice) in a sealed container sublimates to reach a supercritical state under high-temperature, high-pressure conditions (Fig. 1). A supercritical state is when substances have properties of both gas and liquid and they are thus called Supercritical Fluid. This point in the phase diagram with corresponding conditions of temperature and pressure is called Critical Point (Table 1). Supercritical fluids have a higher coefficient of diffusion with lower density and viscosity than liquids. Consequently, higher diffusivity results in better distribution into the mobile phase and better separation, and lower viscosity makes SFC a faster method than HPLC.

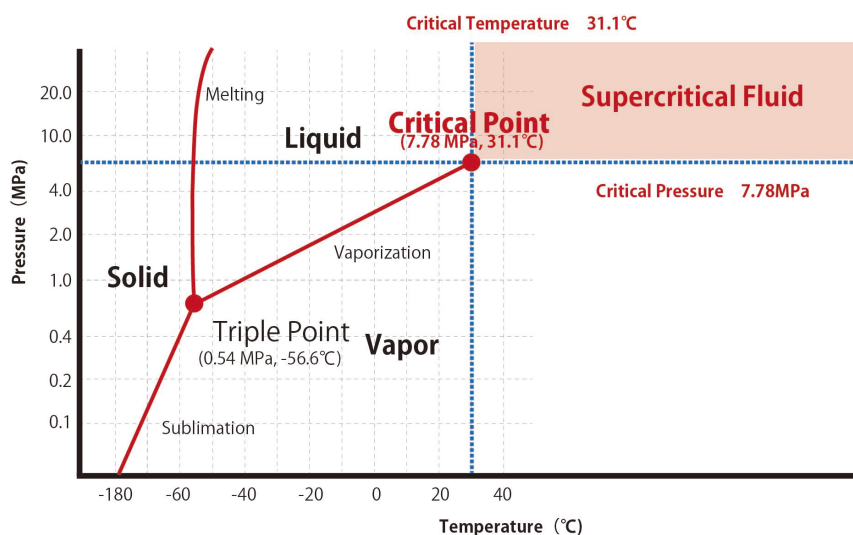


Fig. 1 Phase Diagram for Carbon Dioxide

Table 1 Physical Properties of Supercritical Fluids

	Diffusion Coefficient (m <sup>2</sup> /s)	Density (g/cm <sup>3</sup> )	Viscosity (g/cm·s)
Liquids	10 <sup>-6</sup>	1	10 <sup>-2</sup>
Supercritical Fluids	10 <sup>-3</sup>	0.2-0.8	10 <sup>-3</sup>
Gases	10 <sup>-1</sup>	10 <sup>-3</sup>	10 <sup>-4</sup>

## 1-1-2 Supercritical Carbon Dioxide

Critical pressures and critical temperatures for generating a variety of supercritical fluids are shown in Fig. 2. Water is commonly used for high-performance liquid chromatography, but due to its critical pressure of 22.1 MPa and critical temperature of 374 °C, it cannot be easily kept under conditions necessary for transitioning to a supercritical state. In contrast, supercritical carbon dioxide has a critical pressure of 7.38 MPa and a critical temperature of 31.1 °C, hence allowing transition to the supercritical state under relatively modest conditions.

Supercritical carbon dioxide offers the following advantages.

### (1) Handling

With critical points at 31.1 °C and 7.38 MPa, CO<sub>2</sub> becomes a supercritical fluid under modest conditions.

### (2) Safety

CO<sub>2</sub> is not flammable and has minimal risk of explosion even if used compressed. CO<sub>2</sub> offers intrinsically low toxicity when compared to other organic solvents.

### (3) Cost

CO<sub>2</sub> is mainly sold as a liquid stored in high-pressure cylinder and in most cases can be acquired at a low cost.

### (4) Environmental Friendliness

Other organic solvents are incinerated for disposal, which emits large amounts of energy and CO<sub>2</sub>. However, liquified CO<sub>2</sub> used in SFC is collected as a byproduct of fermentation processes or chemical plants. That means using liquified CO<sub>2</sub> causes minimal environmental impact.

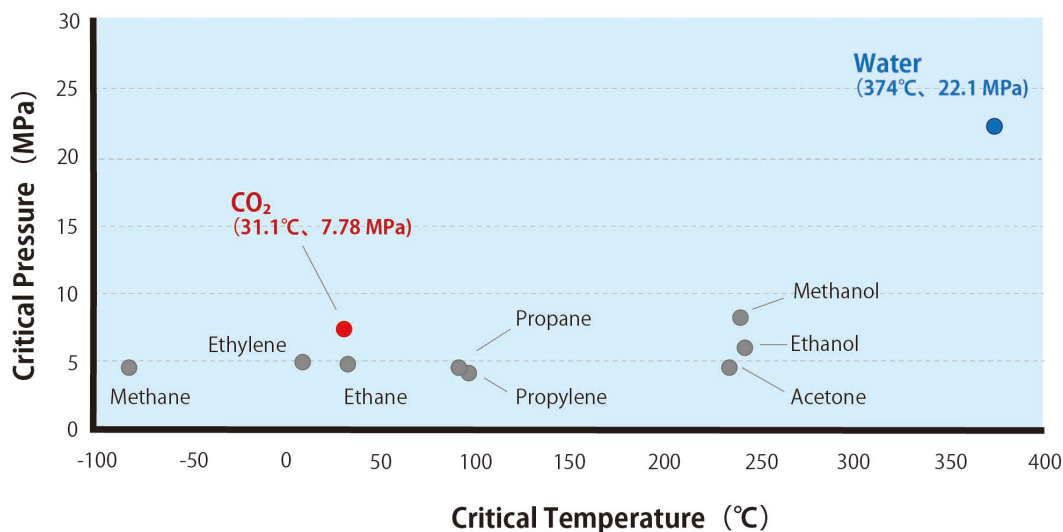


Fig. 2 Conditions for Generating Supercritical Fluids

## 1-2 Supercritical Fluid Chromatography

### 1-2-1 Chromatography in General

Chromatography is the collective term for various techniques used to separate target components in mixture samples, either chemically or physically, so that they can be qualitatively or quantitatively analyzed, collected, or for other purposes. In chromatography, components are separated within two phases, a mobile phase and a stationary phase. If a gas is used as the mobile phase, the technique is referred to as gas chromatography (GC). If a liquid is used, it is referred to as liquid chromatography (LC). Both GC and LC have become essential techniques for analyzing organic compounds.

LC systems that include a pump for pumping the mobile phase, an autosampler for automatically loading samples, and other devices for increasing performance are referred to as high-performance liquid chromatograph (HPLC) systems. In recent years, there has been increasingly widespread use of ultra high performance liquid chromatograph (UHPLC) systems. These systems can sustain very high back pressures generated by use of sub 3 $\mu$  particles as the stationary phase. The low particle size results in ultra fast analysis and sharper peak shapes.

### 1-2-2 Supercritical Fluid Chromatography

Research on supercritical fluid chromatography (SFC) as an analytical technique began in the 1960s. Typically, SFC refers to a chromatographic technique involving using supercritical carbon dioxide as the mobile phase. Since SFC uses CO<sub>2</sub> like a liquid, it is also referred to as "LC using CO<sub>2</sub>."

Before commercialization of SFC, measuring compounds with a wide range of properties required using multiple analytical techniques; however, using a combination of SFC and LC for "unified chromatography" enables simultaneous analysis of heterogeneous components.<sup>\*1</sup>

Fig. 3 shows the simultaneous analysis of 510 standard pesticides in a single injection using SFC. Table 2 shows the octanol-water partition coefficient, repeatability, and linearity results for a few target pesticides. The technique enables accurate simultaneous analysis of compounds with a wide range of polarities. For more details about the analysis in this example, refer to Application News No. L497 "Using the Nexera™ UC Online SFE-SFC-MS System to Analyze Residual Pesticides in Agricultural Products."

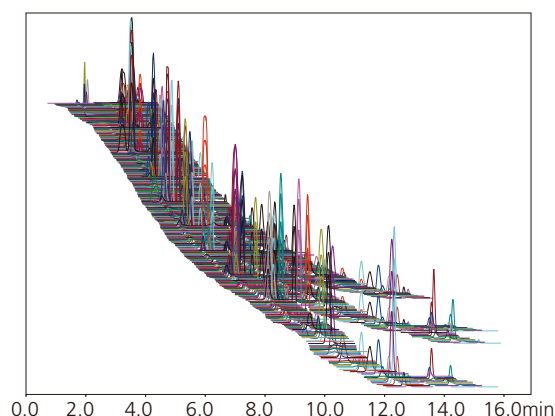


Fig. 3 Mass Chromatogram of Standard Pesticide Mixture Solution

Table 2 Repeatability and Linearity for Pesticide Analysis

Compounds	LogPow	Repeatability (%RSD, n=5)	Range (ng/g)	R <sup>2</sup>
Ethofenprox	6.9	6.1	1-100	0.9991
Hexaflumuron	5.68	6.8	1-100	0.9992
Benzofenap	4.69	1.4	2-200	0.9990
Mepronil	3.66	4.6	1-100	0.9993
Prometryn	3.34	2.7	1-100	0.9994
Fenamidone	2.8	3.0	2-200	0.9991
Ethylchlozate	2.5	3.0	1-100	0.9996
Imazosulfuron	1.6	6.2	1-100	0.9998
Bensulfuron methyl	0.79	8.1	1-100	0.9996
Primisulfuron methyl	0.2	5.5	1-100	0.9994
Halosulfuron methyl	-0.02	5.5	1-100	0.9996
Azimsulfuron	-1.4	4.2	1-100	0.9998

### Application News

#### Using the Nexera UC Online SFE-SFC-MS System to Analyze Residual Pesticides in Agricultural Products

Ever since the positive list system was enacted in 2006 in Japan for residual pesticides in foods, which applies to more than 800 types of pesticides, there has been an increasing demand for a system able to simultaneously analyze multiple pesticides with a wide range of properties, including pretreating samples. This article describes an example of using the Nexera UC online SFE-SFC-MS system to analyze residual pesticides in agricultural products.

[Application News >>](#)

## 1-2-3 Instrumentation for Supercritical Fluid Chromatography

SFC systems involve roughly the same configuration of instruments as HPLC systems. Fig. 4 shows a comparison of flow diagrams for SFC and HPLC systems. The four major differences between SFC and HPLC configurations are described below.

- **Pumps Designed Specifically for Pumping Supercritical Carbon Dioxide**

The carbon dioxide used for SFC is liquified by cooling. Therefore, the delivery pump must have built-in cooling functionality.

- **Back Pressure Regulator (BPR)**

A unit is required to keep the CO<sub>2</sub> in a solvent state (liquid or supercritical fluid) within the system and prevent it from vaporizing by maintaining the pressure level within flow channels. In an SFC system, the BPR unit is positioned downstream from a UV, PDA, or other detector or upstream from an evaporative light scattering detector (ELSD) or mass spectrometer (MS).<sup>\*2</sup> The BPR unit detects pressure in the flow channels and then rapidly opens and closes valves to maintain a constant pressure within the flow channels.

- **Modifier Delivery Pump**

For SFC analysis, an organic solvent such as methanol or acetonitrile (modifier) is pumped for mixing with the supercritical carbon dioxide. That means separate pumps are required for pumping the supercritical carbon dioxide and modifier.

- **Make-up Delivery Pump**

If an ELSD or MS is used for detection or if the system is used for preparative separation, then a solvent (make-up solvent) is pumped to prevent precipitation in the flow channels or to improve the recovery rate of components in separated fractions. Make-up solvent is also pumped to improve sensitivity during MS detection, because supercritical carbon dioxide does not promote ionization.

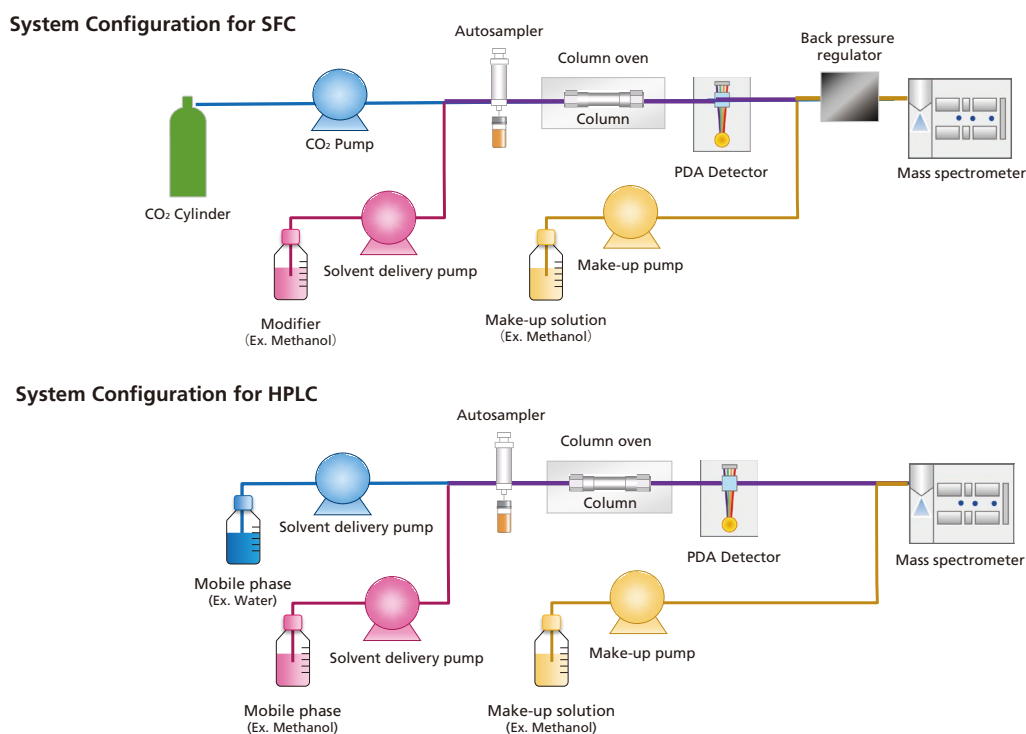


Fig. 4 Comparison of HPLC and SFC System Configurations

## 1-3 Analytical Advantages of SFC

The supercritical carbon dioxide used as a mobile phase for SFC offers unique advantages. Five of these advantages are described below.

### (1) Low Polarity

The polarity of supercritical carbon dioxide is similar to n-hexane.

### (2) Miscibility with Other Organic Solvents

Supercritical carbon dioxide can be mixed with methanol, which is a water-soluble solvent. In contrast, n-hexane, which offers low polarity similar to supercritical carbon dioxide, cannot be mixed with methanol, due to its relatively high polarity.

### (3) Low Viscosity

Supercritical carbon dioxide has lower viscosity than solvents such as water. That results in lower column back pressure, assuming the column has the same internal diameter and length, is packed with the same size particles, and supercritical carbon dioxide is pumped through the column at the same rate.

### (4) Vaporization

Supercritical carbon dioxide evaporates at a temperature and pressure lower than the critical point.

### (5) Solvent Cost

The similar polarity of supercritical CO<sub>2</sub> and n-hexane makes the former appropriate to be used for normal phase separations. Normal-phase chromatography consumes large amounts of solvent. However, because SFC consumes less solvent and involves treating less liquid waste, the overall solvent cost is lower than for HPLC (Fig. 5).

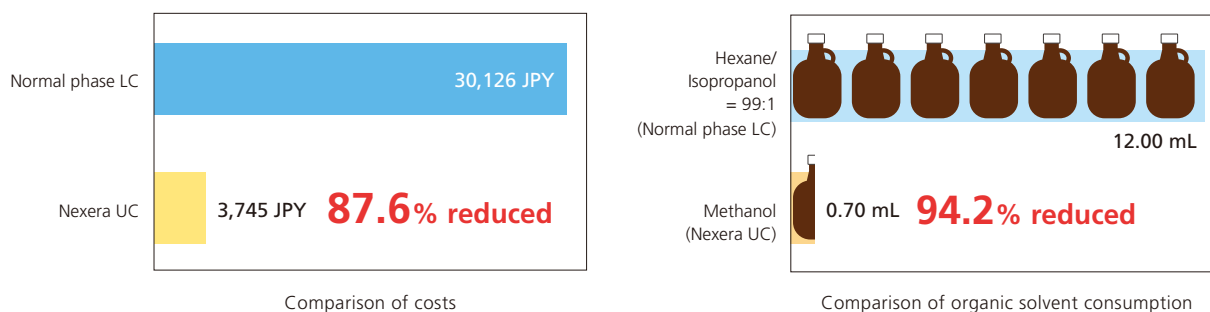


Fig. 5 Comparison of Analytical Costs and Solvent Consumption Quantities for HPLC and SFC

The following sections describe the analytical advantages of SFC provided by the supercritical carbon dioxide advantages (1) to (4) indicated above.



### 1-3-1 Low Polarity — Better Structural Identification

Normal-phase chromatography is a separation method that uses a solvent with low polarity, such as n-hexane, as the mobile phase, and a substance with high polarity, such as silica, as the stationary phase. One advantage of normal-phase chromatography is it offers better structure recognition ability than reversed-phase chromatography.

SFC uses supercritical carbon dioxide, which has similar properties as n-hexane, and also offers excellent structural recognition similar to normal-phase chromatography (Fig. 6). By using a column designed for chiral chromatography, which is packed with a stationary phase with optically active functional groups, SFC is suitable for separation of chiral compounds.

### 1-3-2 Miscibility with Organic Solvents — Separations Covering a Larger Chemical Space

Supercritical carbon dioxide can be mixed with certain organic solvents, such as methanol, ethanol, isopropyl alcohol (IPA), acetonitrile, and tetrahydrofuran (THF). By mixing the supercritical carbon dioxide with such modifiers, SFC can achieve a wide variety of separation patterns. This allows comprehensive analysis for a wide range of low-to-high polarity compounds not previously possible with HPLC. For example, earlier GC was required for analyzing fatty acids and HPLC for glycerides. With SFC, both fatty acids and glycerides can now be analyzed simultaneously due to the hexane-like properties of the supercritical carbon dioxide in combination with a column designed for reversed-phase chromatography (Fig. 7).

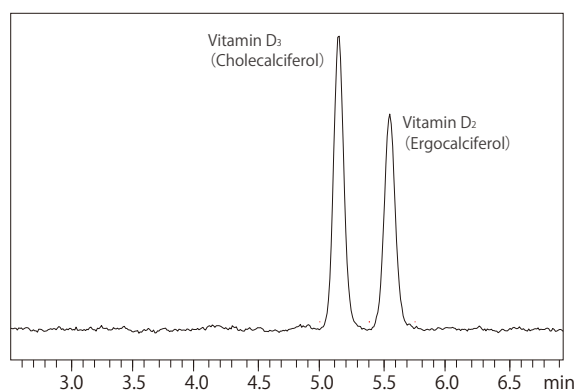


Fig. 6 Separating Vitamins D<sub>2</sub> and D<sub>3</sub> (UC-PyE Column)

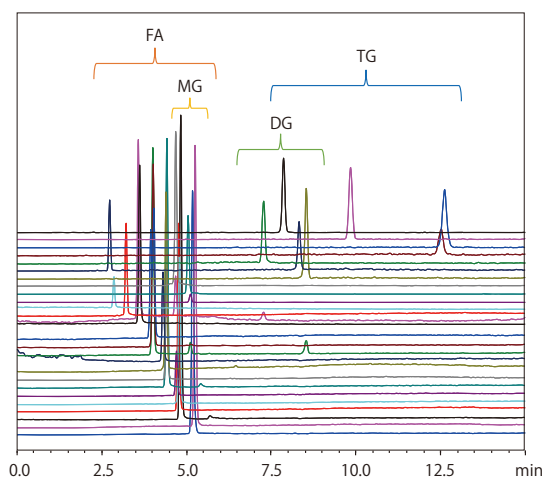


Fig. 7 Example of Simultaneously Analyzing Fatty Acids and Triglycerides by SFC

### 1-3-3 Low Viscosity — Achieving Ultra-Fast Analysis

Supercritical carbon dioxide has both low viscosity and a high diffusion coefficient. That enables analysis with a high linear velocity (flow speed with respect to the column cross-section). Fig. 8 shows a comparison of linear velocity and column pressure plots for HPLC and SFC.

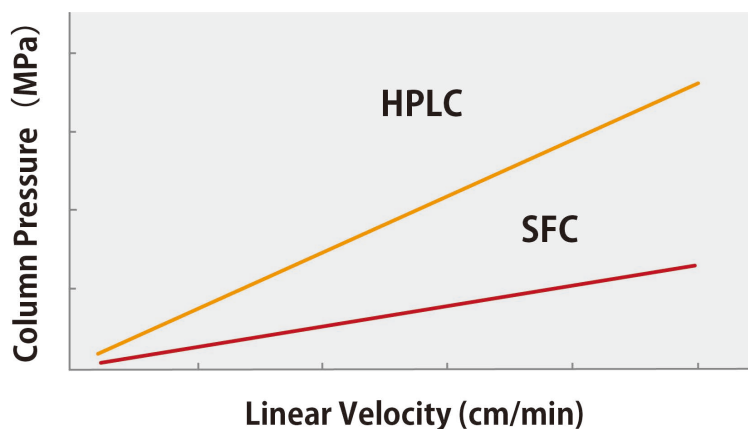


Fig. 8 Linear Velocity vs. Column Back Pressure

Though higher linear velocity can shorten HPLC analysis time, it results in increased column pressure. The upper range of column pressure is limited by instrument design and column capability. CO<sub>2</sub> in the supercritical state has lower viscosity than water and other organic solvents typically used for HPLC. That means shorter analysis times are possible for SFC than HPLC, without loss of resolution and separation pattern.

Supercritical carbon dioxide also has a higher coefficient of diffusion than solvents used for LC. The diffusion coefficient contributes to determining how quickly compounds reach equilibrium between stationary and mobile phases during chromatography and affects resolution. The van Deemter equation below indicates the relationship between physical parameters that determine the height equivalent to a theoretical plate (HETP) value, which is used to indicate column resolution.

$$\text{HETP} = A \cdot dp + \frac{B}{V} + C \cdot dp^2 \cdot V$$

- A: Diffusion coefficient for multiple flow channels
- B: Diffusion coefficient in the axial direction
- C: Resistance to mass transfer coefficient
- V: Linear velocity
- dp: Packing material particle diameter

Plots of the van Deemter equation for HPLC and SFC are shown in Fig. 9.<sup>\*3</sup> In the formula at the bottom of page 8, Term A is based on the multi-flow diffusion generated by the packing material in the column and is proportional to the packing material particle size. Term B is determined by the molecular diffusion in the axial direction of the column. This diffusion phenomenon occurs while compounds stagnate in the stationary or mobile phases and is inversely proportional to linear velocity, regardless of particle size. Term C is derived from how easily compounds can penetrate pores in the packing material and the mass transfer diffusion that contributes to linear velocity. The mass transfer diffusion value is proportional to linear velocity and the square of the particle diameter, which results in higher HETP values for small particles with low linear velocities. In the case of HPLC, mobile phases with a high linear velocity will result in larger mass transfer diffusion values, but supercritical carbon dioxide has a higher diffusion coefficient than HPLC solvents and can penetrate packing material pores more easily. That results in lower mass transfer diffusion values even for high linear velocities.

That makes SFC an analytical technique that can achieve ultra-fast analysis without sacrificing resolution even at high linear velocities.

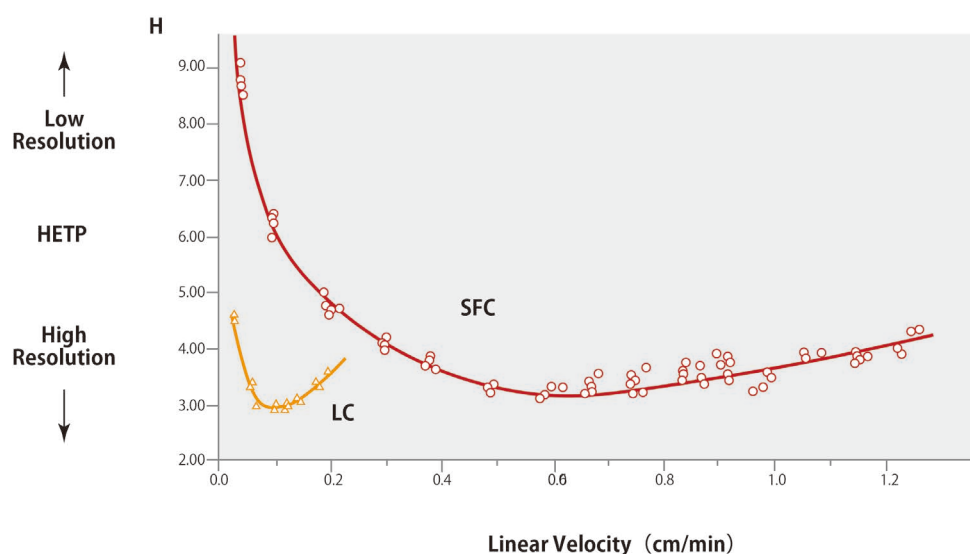


Fig. 9 van Deemter Plot Comparison for LC and SFC

These plots indicate the regions where the optimal linear velocity for SFC is higher than for HPLC. An example comparing the analysis of tocopherols by HPLC and SFC is shown in Fig. 10.<sup>\*4</sup>

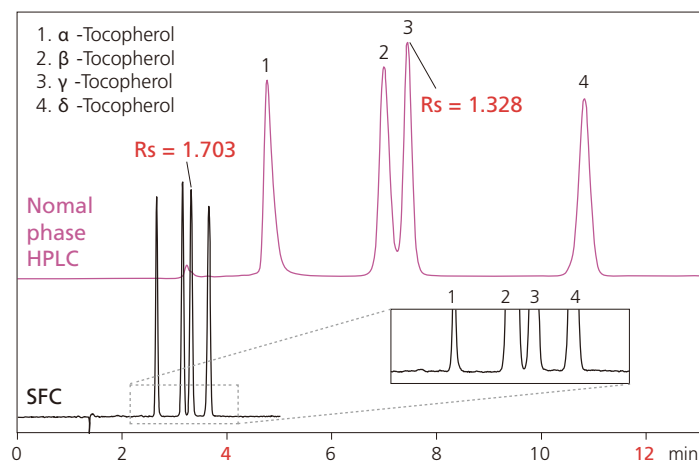


Fig. 10 Comparison Example of Analyzing Tocopherols by HPLC and SFC

### 1-3-4 Vaporization — Improved LC-MS Detection Sensitivity

Under certain temperature and pressure conditions, CO<sub>2</sub> will become a solvent, which is used as a mobile phase for SFC. If an MS is used as a detector, then a back pressure regulator is installed just before it so that only the modifier or make-up solvent reaches the MS. Consequently, unlike conventional LC-MS analysis, this reduces the diluting effects of the mobile phase, potentially increasing the sensitivity while using SFC-MS. This can be seen in the example below where 400 components were analyzed with 90% of these showing higher sensitivity with SFC-MS (Fig 11). Fig. 12 shows a comparison of results from using LC-MS and SFC-MS to measure samples prepared to 10 ppb.

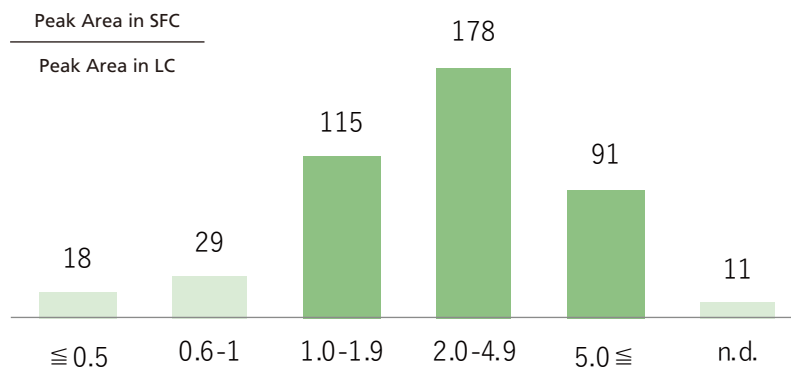


Fig. 11 Area Comparison for LC-MS and SFC-MS

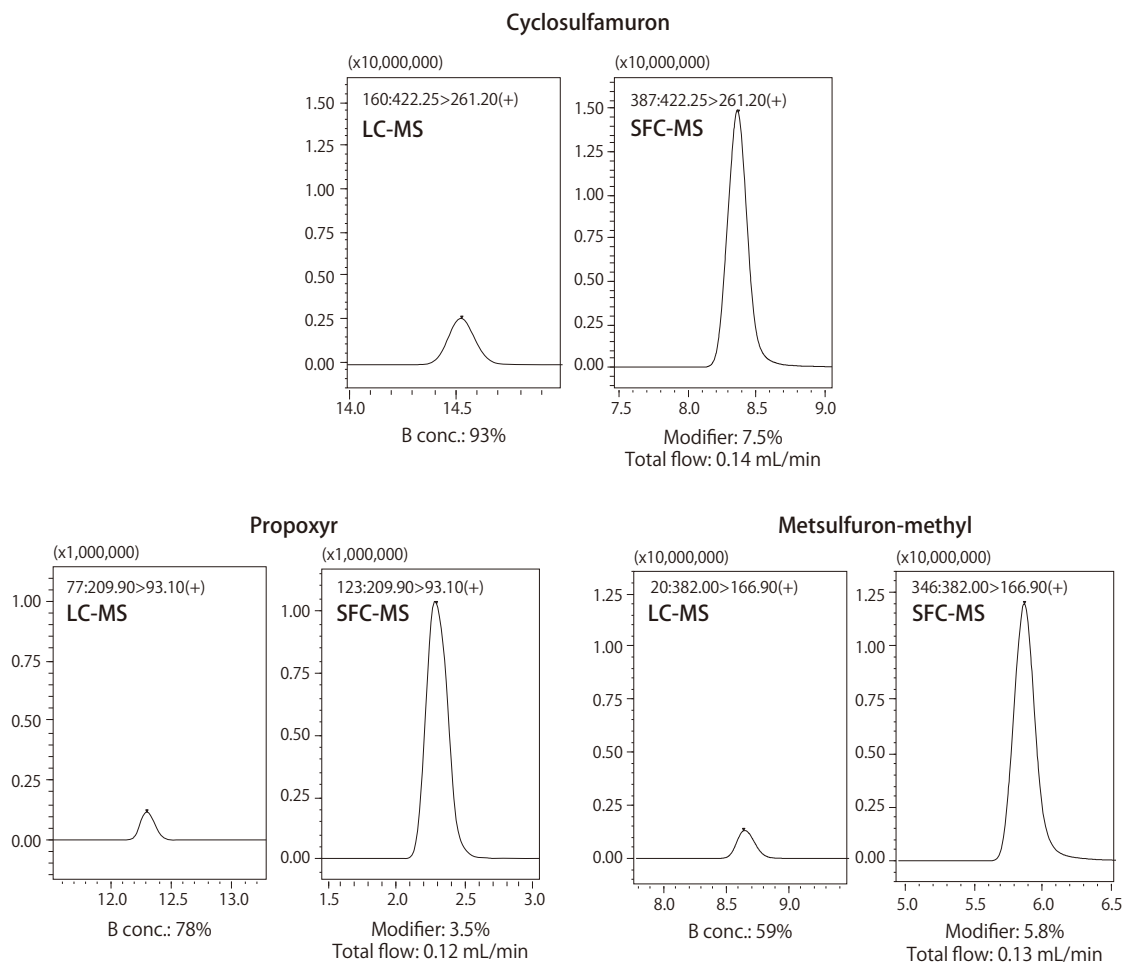


Fig. 12 Comparison of Mass Chromatograms of Standard Pesticide Mixture Solution Measured by LC-MS and SFC-MS

## 1-4 SFC at Preparative Range

GC and HPLC are commonly used to pre-separate samples, but SFC can be used in the same way as well. Preparative separation by SFC offers the following advantages.

### (1) Posttreatment after Preparative Separation

Supercritical carbon dioxide evaporates at ambient temperature and pressure conditions, which eliminates the need for posttreatment.

### (2) Solvent Cost for Large Preparative Separation Quantities

Using inexpensive and environmentally-friendly CO<sub>2</sub> can reduce the cost of purchasing and disposing of solvents.

### (3) Preparative Separation Recovery Rates

After preparative separation, HPLC requires posttreatment steps, such as solvent evaporation or concentration. In contrast, SFC solvents evaporate easily, which minimizes the need for posttreatment and prevents component fragmentation or decomposition during posttreatment.

If SFC is used for preparative separation, an organic solvent that can dissolve target components is sometimes added after column separation to prevent precipitation of target components within flow channels, and gas liquid separators can be used to remove these solvents. A variety of gas-liquid separators are currently being developed by various companies.

Just as with preparative separation by HPLC, SFC involves reviewing results from the analysis scale and then switching to a preparative size column, changing the mobile phase flow rate, and changing the sample injection volume. Equivalent separation can be achieved by increasing the flow rate and injection volume in proportion to the column cross-sectional area. Fig. 13 shows an example of changing from an analysis scale to a preparative scale. By using a column with the same stationary phase, separation equivalent to the analysis scale can be achieved.

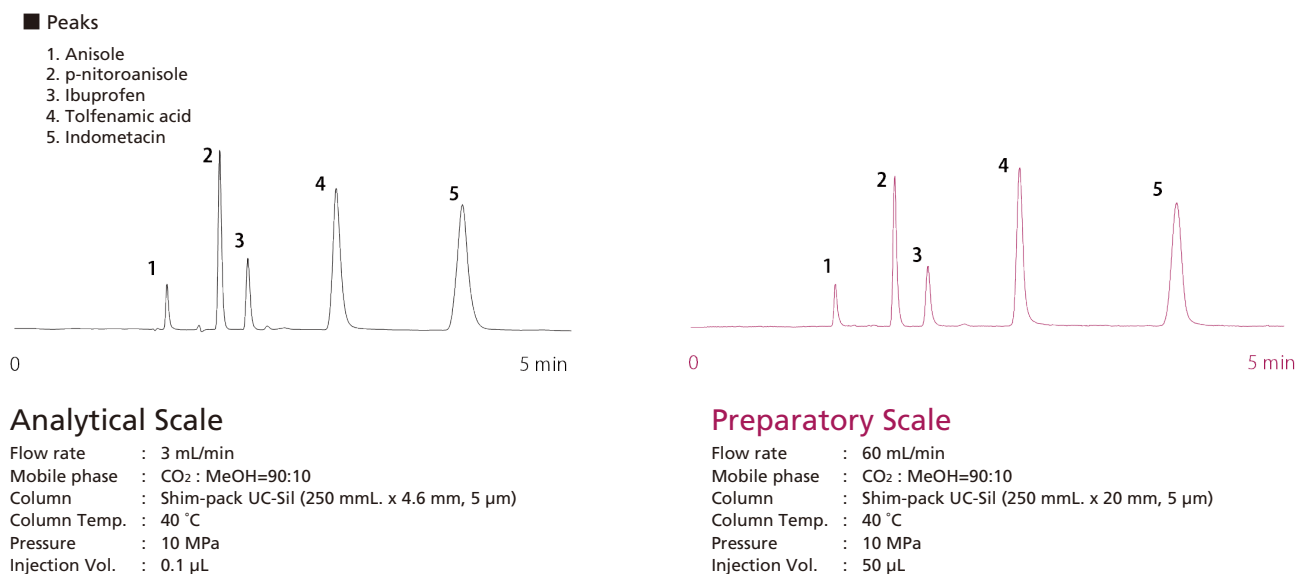


Fig. 13 Migrating from Analysis Scale to Preparative Scale with SFC

A representative example of using supercritical carbon dioxide for supercritical fluid chromatography is from manufacturing decaffeinated coffee. Previously, a solvent was used to remove the caffeine from coffee beans, but there were concerns about the health risks from the residual solvent. Since supercritical carbon dioxide is a safe substance that does not involve any health or environmental concerns, it is now increasingly used for manufacturing decaffeinated coffee.

The following lists some of the advantages of using supercritical carbon dioxide for manufacturing decaffeinated coffee.

- **Advantage 1: High penetration**

With higher penetration ability than water, it can easily pass through coffee beans. Supercritical carbon dioxide can penetrate the interior of coffee beans without having to grind the beans to increase their surface area.

- **Advantage 2: High solubility**

Substances that contact supercritical carbon dioxide are quickly dissolved and extracted from the coffee bean.

- **Advantage 3: Evaporation at ambient temperatures/pressures**


CO<sub>2</sub> residues are unlikely to remain in coffee beans.

- **Advantage 4: Low toxicity**

Even if residues remained in coffee beans, they would not be toxic.

In this way, the properties of supercritical carbon dioxide are used beneficially for manufacturing decaffeinated coffee.

#### Decaffeinated coffee manufacturing process using supercritical fluid CO<sub>2</sub>

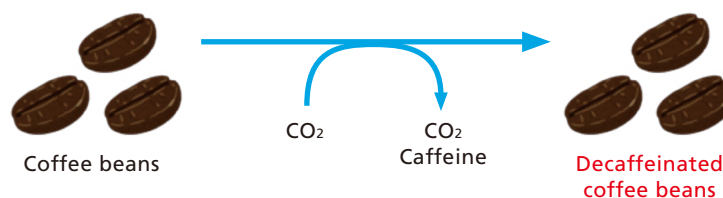

#### The properties of supercritical fluid CO<sub>2</sub> from the process of caffeine removal

##### 1. Penetration



##### Low viscosity

Penetrates into coffee beans.

##### 2. Volatilization




##### Naturally vaporized, no residue


Even if it remains, it is CO<sub>2</sub>, so there is almost no harm.

The method of using supercritical carbon dioxide to remove caffeine from coffee beans is based on supercritical fluid extraction technology used to extract target components from the interior of samples. Previous technologies required pretreatment that involved tissue homogenization, dissolution, and extraction (liquid-liquid extraction with a solid phase cartridge), but supercritical fluid extraction can extract target compounds by simply placing the sample in a chamber and filling it with supercritical carbon dioxide. Furthermore, compounds extracted by supercritical carbon dioxide can be injected directly into columns. That enables the automation of all steps from pretreatment and separation to analysis.

#### Conventional (QuEChERS) Pretreatment Method

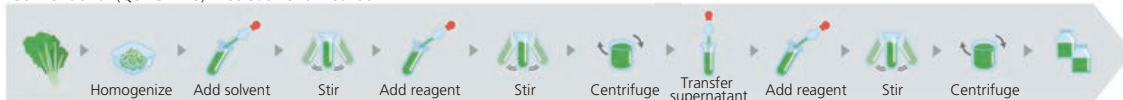

#### Nexera UC Pretreatment

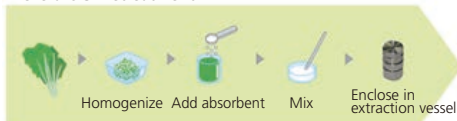

## Chapter 2 Nexera UC Supercritical Fluid Chromatograph

The Nexera UC supercritical fluid chromatograph product line includes SFC systems, SFC/UHPLC switching systems, preparative SFC systems, online SFE-SFC systems, and offline SFE pretreatment systems. Various Nexera UC systems can be configured by adding a newly developed SFC CO<sub>2</sub> solvent delivery unit (LC-30AD<sub>SF</sub>), a back pressure regulator (SFC-30A), and/or an SFE extraction unit (SFE-30A) to a base Nexera ultra high performance liquid chromatograph.

This chapter describes the instruments used in Nexera UC systems and as system extension technologies.

### 2-1 Components of SFC System

The SFC system product line includes a variety of systems, such as SFC-UV systems that use UV (or PDA) detection, UFMS (SFC-MS) systems that include a mass spectrometer capable of SFC speeds, and screening systems that can automatically switch between multiple columns and modifiers for considering various analytical conditions.

The appearance of a Nexera UC UFMS system is shown in Fig. 14.

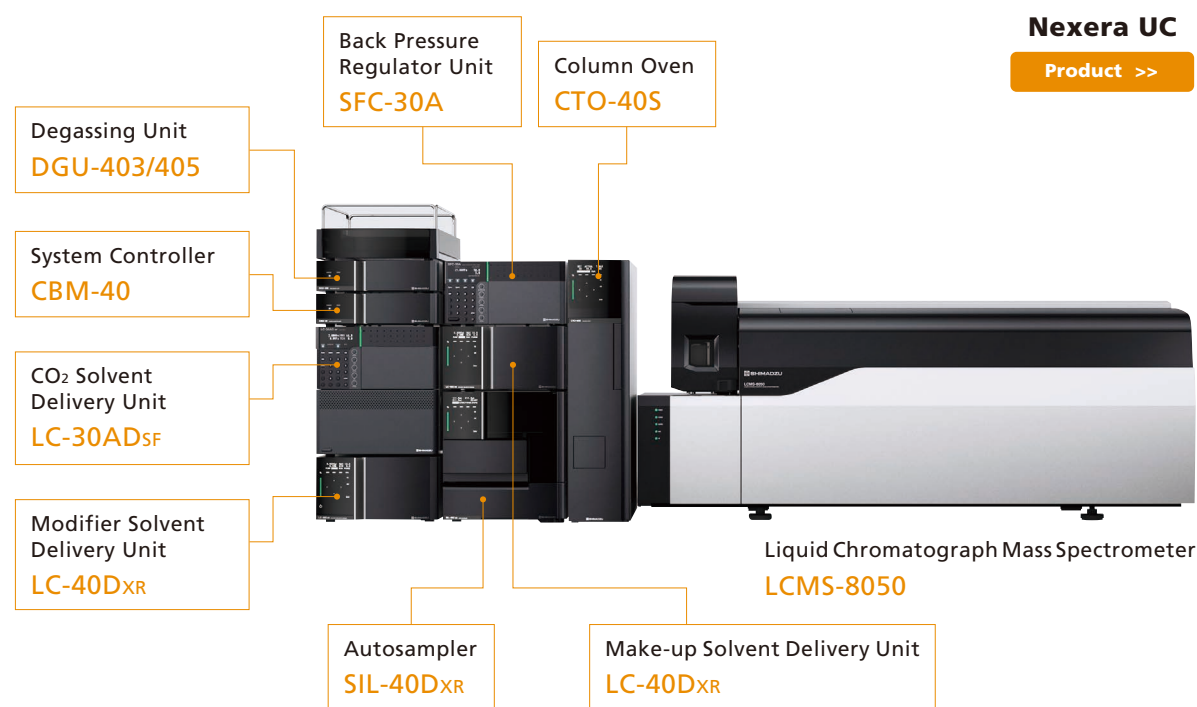


Fig. 14 Nexera UC UFMS System

## 2-2 Nexera UC/s: SFC/UHPLC Switching System

If both UHPLC and SFC are available during analytical condition development, optimal separation parameters can be selected from a wider range of potential separation parameters. With the Nexera UC/s system, the same system can be used for both UHPLC and SFC analysis modes. Fig. 15 shows the flow diagram of the Nexera UC/s system, which integrates both UHPLC and SFC capability in a single system. The system is configured with LC-30AD<sub>SF</sub> and SFC-30A units added to a standard UHPLC system. Both UHPLC and SFC analyses can be performed by switching between solvent delivery units and switching ON or OFF the pressure control mode of the back pressure control valve. That can reduce installation space and initial cost and increase the instrument utilization rate. By using the mobile phase solvent switching valve in combination with the column switching valve, mobile phase conditions can be changed automatically and continuously for up to twelve columns, which enables comprehensive measurements that improve method development efficiency.

For more information about Nexera UC/s systems, refer to Technical Report (C190-E212) "Improving Separation and Method Development Efficiency Using the Nexera UC/s UHPLC/SFC Switching System"

A chiral screening system designed specifically for separating chiral compounds is also available. It offers automatic modifier blending functionality for up to 12 columns using only the SFC mode, which eliminates the need to switch between UHPLC and SFC modes. A specific example of using this system to determine separation conditions for chiral compounds is described in Chapter 4.

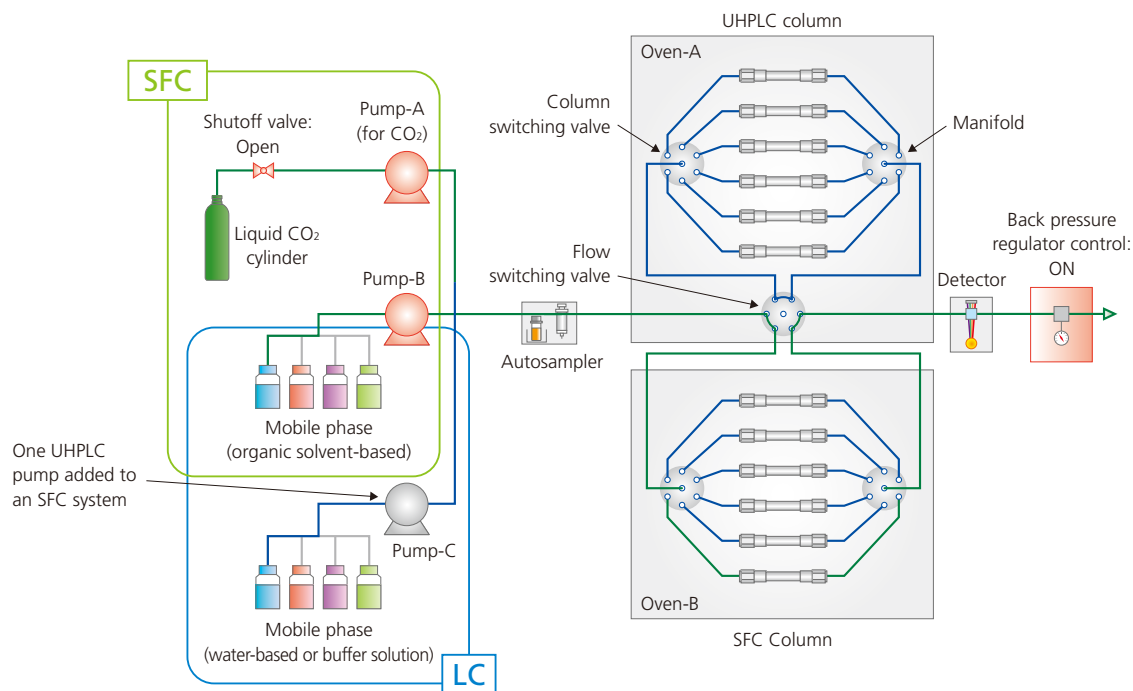


Fig. 15 Flow Diagram for Nexera UC/s: SFC/UHPLC Switching System (Green Line Indicates Section Used during SFC.)

### Technical Report

#### Improving Separation and Method Development Efficiency Using the Nexera UC/s UHPLC/SFC Switching System

Due to the difference in separation selectivity between ultra high performance liquid chromatography (UHPLC) and supercritical fluid chromatography (SFC) using supercritical carbon dioxide and organic solvents, SFC has garnered significant interest in recent years as a new separation method. By combining both UHPLC and SFC, the Nexera UC/s UHPLC/SFC switching system enables both UHPLC and SFC analysis using a single system. Since SFC has different separation characteristics than UHPLC, it can improve separation between isomers. By using two separation methods, UHPLC and SFC, for screening during the method development process, the system can result in determining superior analytical conditions more quickly. This article describes using a Nexera UC/s UHPLC/SFC switching system to evaluate analytical conditions for separating optical isomers.

[Technical Report >>](#)



## 2-3 Preparative SFC Systems

Preparative LC is a technique used to purify samples for specific target components. It is used in a wide range of fields, including chemicals, pharmaceuticals, and foods. Preparative LC serves as a powerful tool for achieving higher purity and recovery rate levels of target components, but it requires drying and powderizing steps. SFC can improve preparative workflow efficiency by significantly reducing the amount of work involved in the powderization process after preparative purification. The Nexera UC product line includes three systems—a stacked fraction system intended for large-volume fractionation, a multi-fraction system for separating multiple peaks, and an analytical fraction system intended for analysis-scale fractionation.

These products also feature Shimadzu's unique gas-liquid separator (patented) that inhibits sample dispersion and carryover to obtain high recovery rates. Therefore, excellent recovery rates can be achieved even for highly volatile compounds, such as the fragrance linalool, regardless of the flow rate or modifier concentration. The following sections describe the three product lines and the gas-liquid separator.

### Nexera UC Prep

[Product >>](#)

### 2-3-1 Stacked Fraction System

This system is optimized for large volume fractionation that involves repeatedly injecting compounds with up to several components (Fig. 16). The FRS-40 injection and collection unit can function as both an injector and fraction collector for repeatedly injecting the same sample to collect gram-level preparative fractions. It can inject volumes up to 20 mL\* and collect ten fractions. It supports flow rates from 10 to 150 mL/min and connecting columns with an internal diameter from 10 to 30 mm. For more information about stacked fraction systems, refer to Technical Report (C190-E247) "Improved Efficiency of Isomer Preparative Operations by Supercritical Fluid Chromatography with Stacked Injection"

\* When using optional parts

### 2-3-2 Multi-Fraction System

This system is suitable for applications involving fractionation of multiple components in samples with many peaks detected, such as for impurities in pharmaceuticals (Fig. 17). Volumes up to 2 mL\* can be injected using an autosampler that holds up to 162 samples (when using 1.5 mL vials). The FRC-40 SF fraction collector can collect up to 540 fractions (using 10 mL vials). It supports flow rates from 10 to 150 mL/min and connecting columns with an internal diameter from 10 to 30 mm.

\* When using optional parts



Fig. 16 Nexera UC Prep Stacked Fraction System



Fig. 17 Nexera UC Prep Multi-Fraction System

### 2-3-3 Analytical Fraction System

This system enables analysis-scale fractionation for applications that only require fraction quantities up to several milligrams, such as for checking synthesis (Fig. 18). By connecting an FRC-40 SF fraction collector to the Nexera UC system, the same system can be used for applications ranging from determining analytical conditions by method scouting to preparative separation of about a few milligrams. With a maximum flow rate of 5 mL/min, it supports using analytical size columns.



Fig. 18 Nexera UC Analytical Fraction System

### 2-3-4 LotusStream™ Gas-Liquid Separator

When carbon dioxide transitions from the supercritical fluid state to the gas state during preparative SFC, its volume immediately expands by about 500 times, which can cause eluate to splatter from the column, a factor leading to lower recovery rates.

The newly developed LotusStream (patented) gas-liquid separator uses multiple flow channels to limit the flow rate without increasing the tubing diameter. As a result, the carbon dioxide is discharged externally and the liquid travels along the column and then drips directly below without the eluate splattering (Fig. 19). For more information about LotusStream technology, refer to Technical Report (C190-E250) "Evaluating the Performance of the LotusStream Gas-Liquid Separator for Preparative Supercritical Fluid Chromatography"

[Movie >>](#)

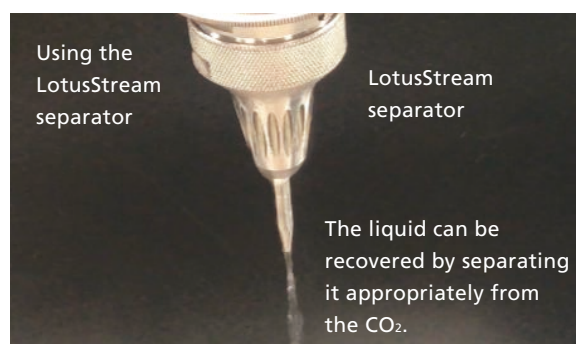
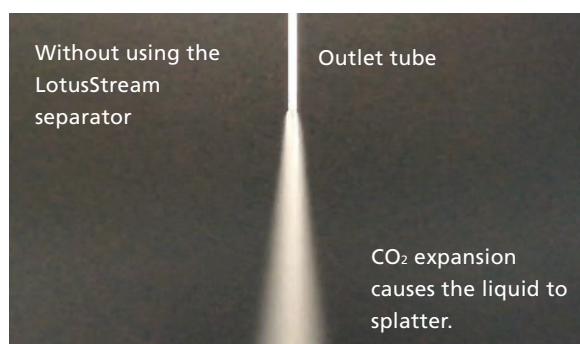


Fig. 19 Gas-Liquid Separation with and without the LotusStream Separator

#### Technical Report

#### Evaluating the Performance of the LotusStream Gas-Liquid Separator for Preparative Supercritical Fluid Chromatography

Preparative supercritical fluid chromatography (SFC) is one of the most common technologies in preparative purification. Unlike preparative liquid chromatography, preparative SFC is prone to low recovery rates because the liquid splatters if it is not appropriately separated from the CO<sub>2</sub> during recovery. This report provides an evaluation of recovery rates, carryover, and contamination using the newly developed LotusStream gas-liquid separator in the Nexera UC Prep preparative supercritical fluid chromatograph system.

[Technical Report >>](#)

## 2-4 Online SFE-SFC System

Online SFE-SFC is extraction and separation technology used to extract components by supercritical fluid extraction (SFE) and then inject the extract directly into a column for separation by SFC. The system flow diagram and the process flow from extraction to analysis are shown in Fig. 20. After using a valve inside the SFE unit to switch flow channels so that supercritical carbon dioxide flows into the extraction vessel, components can be extracted by letting the solution stand without solution flowing through the vessel (static extraction) or with solution flowing through the vessel (dynamic extraction). In the case of online SFE-SFC, the sample is sent to the analytical column during dynamic extraction. For more information about the operating principle of online SFE-SFC systems, refer to Technical Report (C190-E205) "Online Supercritical Fluid Extraction-Supercritical Fluid Chromatography (Online SFE-SFC)"

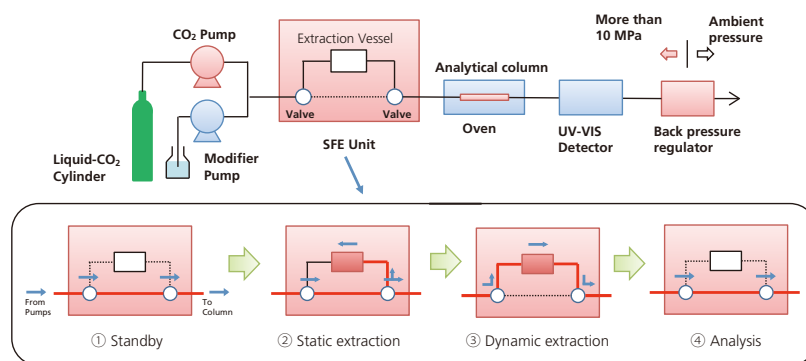


Fig. 20 Nexera UC Online SFE-SFC-MS System Flow Diagram

The entire online SFE-SFC process, from extraction to separation and detection, can be completed within a single system, which eliminates the need for any complicated pretreatment processes and enables automation. This can significantly reduce the time and effort required for the various operations involved in analysis. It also means that the entire process, from extraction to separation and detection, can be performed without exposure to light, without oxidation, and in a moisture-free environment. Therefore, the method is extremely useful for analyzing unstable compounds, such as compounds with components easily decomposed by light, easily oxidized, or easily hydrolyzed. Unlike offline SFE, online SFE-SFC does not involve preparing sample solutions, which means sensitivity can be increased easily because target components are not diluted by the sample solvent. Fig. 21 shows a comparison of analyzing reduced coenzyme Q10 by online SFE-SFC versus solvent extraction-SFC. For more information about this analysis, refer to Application News No. L496 "Analysis of Unstable Compounds Using Online SFE-SFC"

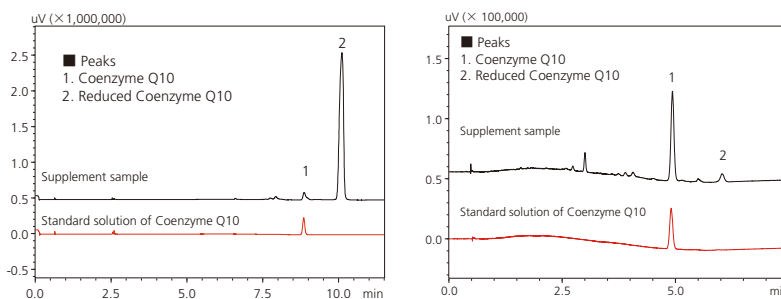


Fig. 21 Example of Analyzing Reduced Coenzyme Q10 (Left: Online SFE-SFC; Right: Solvent Extraction-SFC)

### Technical Report

#### Online Supercritical Fluid Extraction Supercritical Fluid Chromatography

An online supercritical fluid extraction-supercritical fluid chromatography (online SFE-SFC) system directly connects the supercritical fluid extraction (SFE) with supercritical fluid chromatography (SFC).

[Technical Report >>](#)

### Application News

#### Analysis of Unstable Compounds Using Online SFE-SFC

The reduced coenzyme Q10 is easily oxidized to form oxidized coenzyme Q10. This article compares the results from using both solvent extraction-SFC and online SFE-SFC analysis to analyze the reduced coenzyme Q10 contained in a vitamin supplement.

[Application News >>](#)

Conventional pretreatment methods by solvent extraction involve adding a solvent to the sample to extract components, which dilutes the extract solution and results in lower concentrations of compounds in the extract solution. Furthermore, for analysis-scale flow rates, only a few microliters of extract can be injected due to the small internal diameter of columns specified. Therefore, if a sample with a low extract concentration is analyzed, extraction, concentration, and other operations must be performed repeatedly due to inadequate sensitivity. In contrast, the Nexera UC online SFE-SFC system can inject extracts directly into the column, which results in injecting a much larger quantity of the sample into the column than when using solvent extraction. A comparison example\* is shown in Fig. 22. Direct injection with a Nexera UC system enables high-sensitivity analysis because it injects 10,000 times more of the sample into the column than when using 10 mL of solvent to inject 1  $\mu$ L by solvent extraction.

\* Assumes equivalent extraction efficiency and a 1 % concentration of the target compound in the sample.

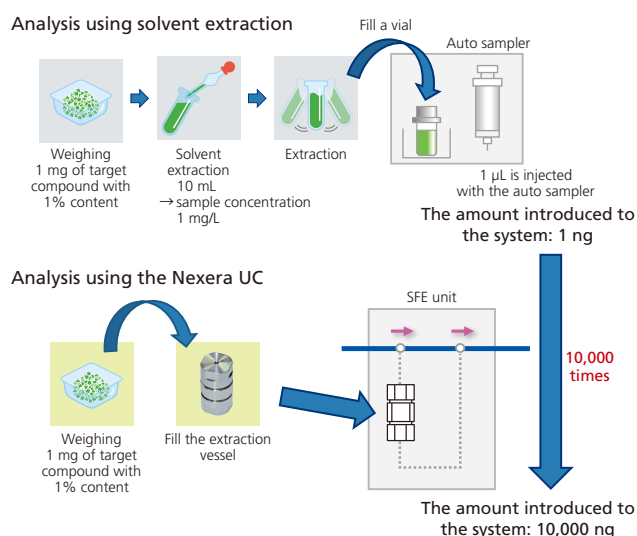


Fig. 22 High-Sensitivity Extract Analysis Using an Online SFE-SFC-MS System

Techniques for analyzing residual pesticides in foods must be quick and easy, including for pretreatment, due to the large number of regulated pesticides involved and to shorten inspection times. Using an online SFE-SFC system, the only operations involved are pretreating the sample by mixing it with a dehydrating agent and then filling the extraction vessel. Results from pesticide analysis using the online SFE-SFC-MS system are shown in Table 3. The system can be used for pesticides with a wide range of polarities, from hydrophilic compounds with a -1.4 LogPow value to hydrophobic compounds with a 6.9 LogPow value.

Table 3 Repeatability and Linearity of Pesticide Analysis Using an Online SFE-SFC-MS System

Compounds	LogPow	Repeatability (%RSD,n=5)	Range (ng/g)	R <sup>2</sup>	Compounds	LogPow	Repeatability (%RSD,n=5)	Range (ng/g)	R <sup>2</sup>
Ethofenprox	6.9	6.1	1-100	0.9991	Imazosulfuron	1.6	6.2	1-100	0.9998
Hexaflumuron	5.68	6.8	1-100	0.9992	Bensulfuron methyl	0.79	8.1	1-100	0.9996
Benzofenep	4.69	1.4	2-200	0.9990	Primisulfuron methyl	0.2	5.5	1-100	0.9994
Mepronil	3.66	4.6	1-100	0.9993	Halosulfuron methyl	-0.02	5.5	1-100	0.9996
Prometryn	3.34	2.7	1-100	0.9994	Azimsulfuron	-1.4	4.2	1-100	0.9998
Fenamidone	2.8	3.0	2-200	0.9991					
Ethlychlozate	2.5	3.0	1-100	0.9996					

## 2-5 Offline SFE-SFC System

This is an offline SFE pretreatment system designed specifically for extracting target components from solid samples. The extraction process flow is shown in Fig. 20. The extraction vessel is packed with the sample and placed in the SFE unit (Fig. 23A). The extraction vessel is filled with supercritical carbon dioxide for static extraction (Fig. 23B). Then supercritical carbon dioxide is pumped through the extraction vessel for dynamic extraction (Fig. 23C). The extract is trapped in a trap column and the eluate solution that contains the target components is collected in the fraction collector (Fig. 23D) for measurement with the specified analytical instruments, such as LC-MS, GC-MS, or FTIR. For more information about the operating principle of offline SFE pretreatment systems, refer to Technical Report (C190-E199) "Improved Sample Pretreatment Using Offline Supercritical Fluid Extraction"

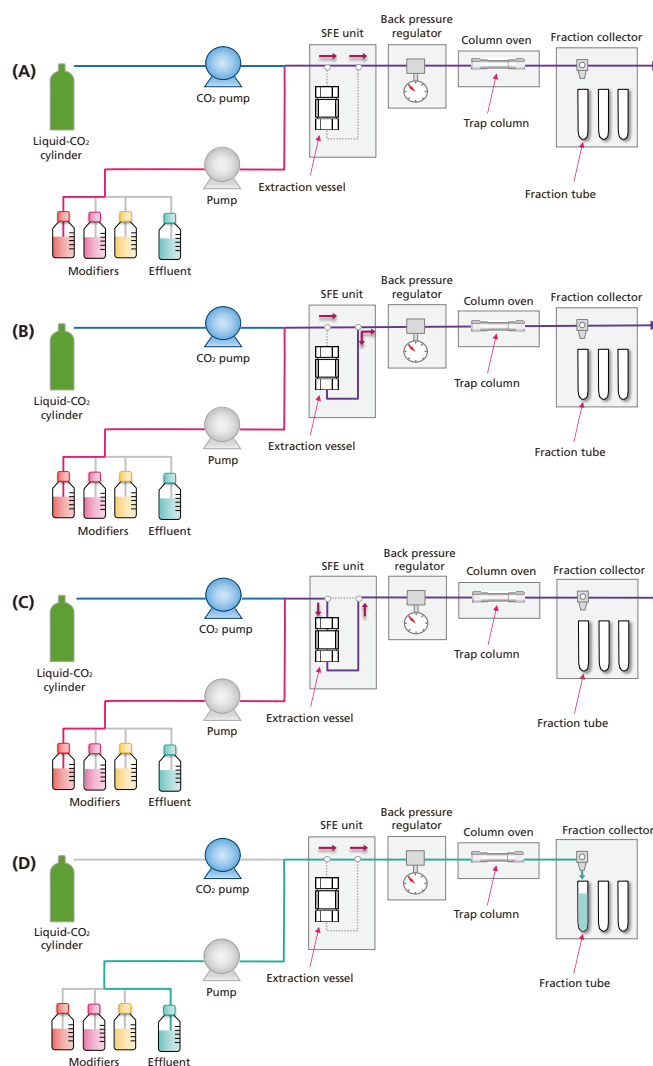


Fig. 23 Extraction Process Flow in Nexera UC Offline SFE Pretreatment System

### Technical Report

#### Improved Sample Pretreatment Using Offline Supercritical Fluid Extraction

Typically, components are extracted from solid samples using methods such as dissolution or solid-liquid extraction, but those methods are not well suited for processing multiple analytes, because of the long time required for extraction and the tedious steps involved in processes other than extraction, such as cleaning/preparing extraction equipment. However, due to the superior ability of supercritical fluids to dissolve and penetrate other substances, extraction methods using supercritical fluids can extract target components from solids with high efficiency. They also enable the automation of extraction steps. This article describes how to improve the efficiency of pretreatment for analysis using a Nexera UC SFE pretreatment system configured specifically for pretreatment for offline supercritical fluid extraction.

Technical Report >>

## Chapter 3 SFC Analytical Condition Parameters

The approach for the SFC analytical condition development process is basically the same as for determining analytical conditions for HPLC. HPLC can generally target any component in compounds that can be dissolved in the mobile phase. In contrast, SFC can be used to analyze any compound that is compatible with supercritical carbon dioxide and can be dissolved in an organic solvent.

If chromatography is used for quantitative analysis, then optimal analytical conditions should achieve resolution of 1.5 or more. The basic equation for calculating resolution  $R_s$  is indicated below.

$$R_s = \frac{1}{4} (\alpha - 1) \cdot \frac{k'}{1+k'} \cdot \sqrt{N}$$

$\alpha$ : Separation factor

$k'$ : Retention factor

$N$ : Number of theoretical plates

The  $\alpha$ ,  $k$ , and  $N$  values that affect  $R_s$  are all independent factors that vary depending on the column and mobile phase used, where  $N$  represents separation efficiency and  $\alpha$  selectivity. In particular, separation can be improved by increasing  $N$  and  $\alpha$  values.

Because the supercritical carbon dioxide used for SFC offers different properties than mobile phase solvents used for HPLC, when a column used for HPLC is used for SFC, it can sometimes result in different elution behavior (Fig. 24). This chapter describes the characteristics of columns and mobile phases involved in determining analytical conditions for SFC.

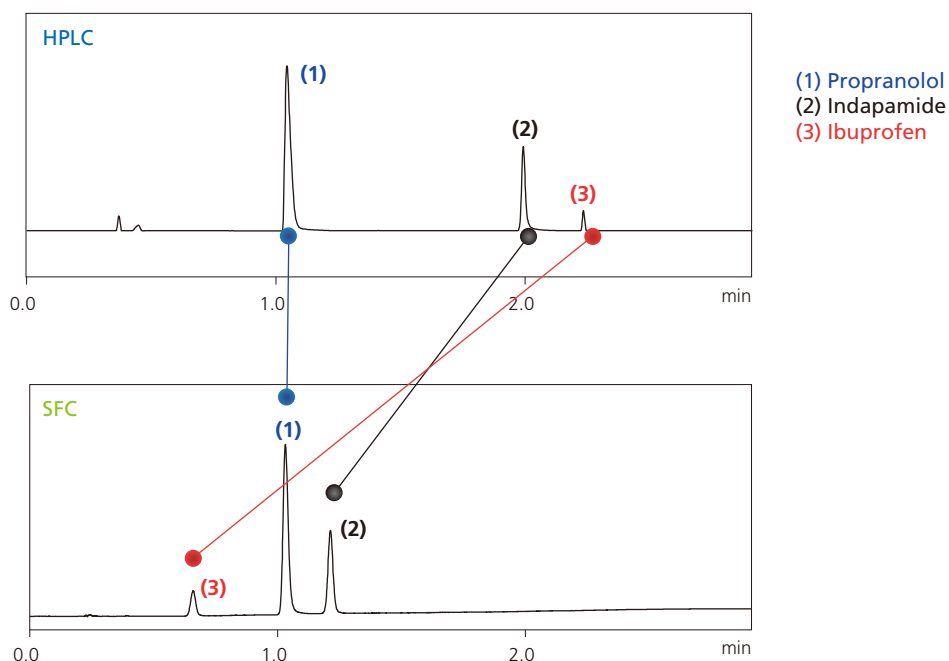


Fig. 24 Separation of Three Types of Pharmaceutical Ingredients (Upper: HPLC; Lower: SFC)

### Technical Report

#### Study of Effective Parameters for Separation of Isomers of Fluorescent Substance by Supercritical Fluid Chromatography

In order to detect the sites of trace amounts of cancers within the bodies of human patients in a short time, research and development is being performed with respect to fluorescent probes, molecules that emit light only after the reactions with biomarker enzymes selective for cancer sites. With rhodamine and fluorescein, which are commonly used as the nuclei for fluorescent probes, when molecules with substituents introduced by typical methods are synthesized, isomers are created that differ with respect to the site of substitution on the benzene ring region. In order to be used as the nuclei for fluorescent probes, it is important that these isomers be separated. Separation of isomers of the fluorescent substance carboxytetramethylrhodamine was optimized using the Nexera UC supercritical fluid chromatograph. This article will discuss the factors with an impact on the separation of the isomers.

Technical Report >>

### 3-1 SFC Columns

In principle, the various columns used for HPLC can also be used for SFC, but in order to ensure the CO<sub>2</sub> is kept in a supercritical fluid state, make sure the column has the pressure capacity necessary for maintaining the required pressure level. Polymer-based columns require particular care because the pressure increases during use as the stationary phase expands.

Since supercritical carbon dioxide has about the same polarity as hexane, the main SFC separation modes will presumably behave similar to normal-phase chromatography, but unlike HPLC, columns with properties of both normal-phase columns (such as silica, diol, and CN) and reversed-phase columns (such as C18, cholesteryl, and phenyl) can be used. Furthermore, the supercritical carbon dioxide used for SFC has high mobile phase diffusivity and low molecular density, making it more prone to secondary interactions than HPLC and resulting in major changes in separation behavior depending on the column stationary phase.<sup>\*5</sup> Therefore, to determine optimal conditions, it is important to consider analytical conditions using a wide variety of stationary phase types.

Though supercritical carbon dioxide has low polarity similar to hexane, as mentioned above, it can be mixed with highly polar solvents. That means when optimizing the mobile phase and stationary phase for SFC, normal-phase columns can be replaced with a reversed-phase column without changing the type or composition of mobile phase or modifier. Fig. 25 shows a comparison of the elution order for HPLC and SFC. For HPLC, mobile phases are matched to either the reversed-phase or normal-phase separation mode based on mobile phase properties. That involves purging the solvent from inside the system each time the separation mode is changed. For SFC, the same mobile phase conditions can be used for either the reversed-phase or normal-phase separation modes, which means separation results can be checked by simply replacing the column.

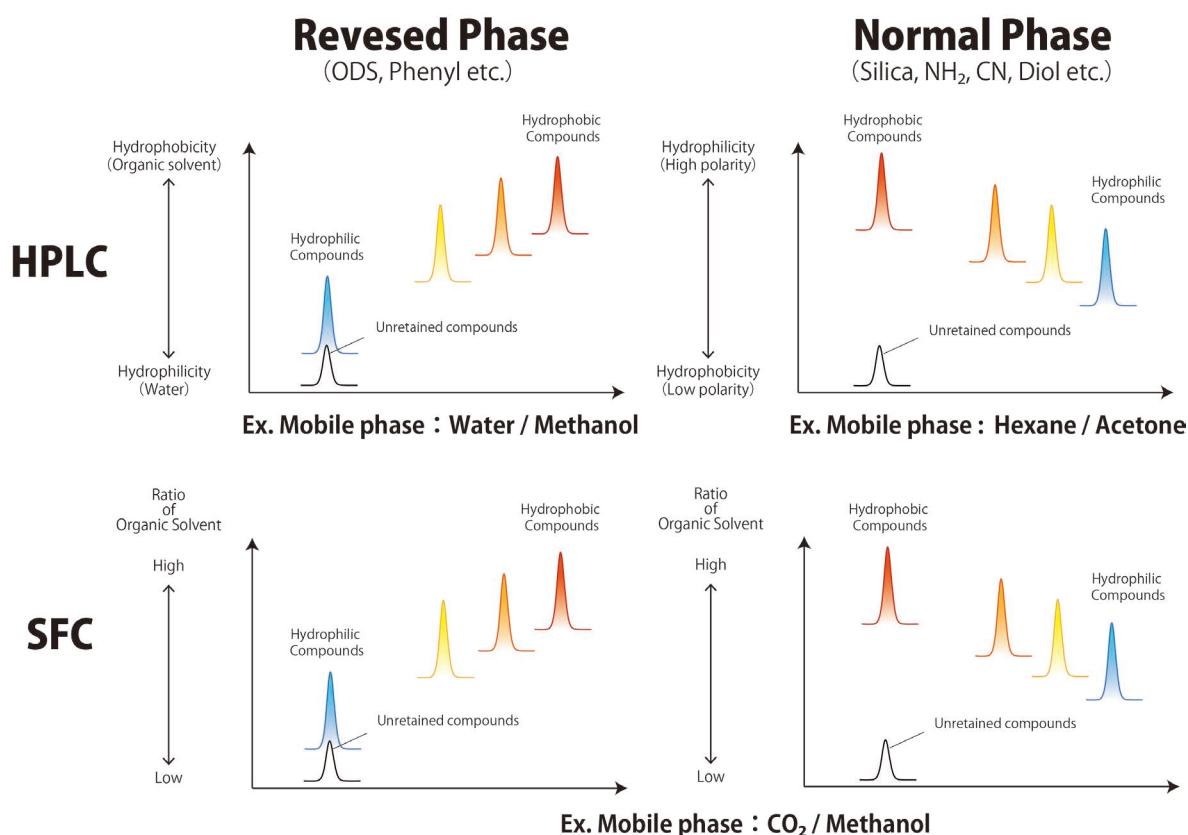


Fig. 25 Comparison of Elution Order for HPLC and SFC

### 3-1-1 Column Selection Guide

Given that normal-phase separation is the main separation mode used for SFC, normal-phase UC-Diol and UC-Diol II columns are widely used as the first choice for SFC systems. The next most widely used columns are UC-Py columns, due to the similarity in separation behavior to ethylpyridine-based columns. Whereas HPLC uses aqueous or non-aqueous mobile phases with significantly different compositions for reverse-phase and normal-phase analysis respectively, SFC uses the same mixture of supercritical fluid carbon dioxide mobile phase and a modifier (a water-soluble organic solvent such as methanol) regardless of the stationary phase. That means the same mobile phase composition can be used for serial analysis involving all types of columns. Since supercritical carbon dioxide can penetrate tiny pores more easily, it promotes higher interaction with the stationary phase. That means SFC offers the possibility of separating isomers or other compounds that are difficult to be separated by HPLC. In particular, columns with unique or multiple interactive effects can help improve separation using SFC.

Fig. 26 shows an index of elution order and retention behavior for a variety of columns used to analyze ibuprofen (acidic compound), indapamide (neutral compound), and propranolol (basic compound). Fig. 27 shows chromatograms obtained using diol group, hydroxyphenyl group, and cholesteryl group columns.

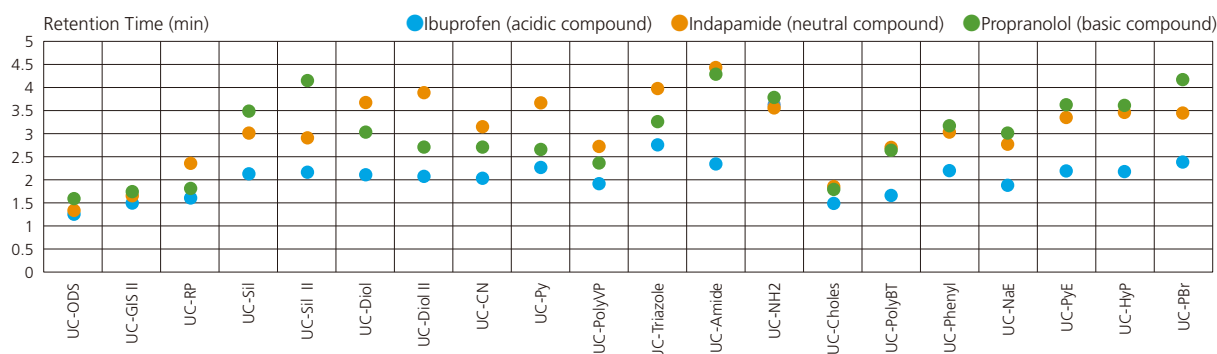


Fig. 26 Retention Behavior of Respective Shim-pack™ UC Series Products

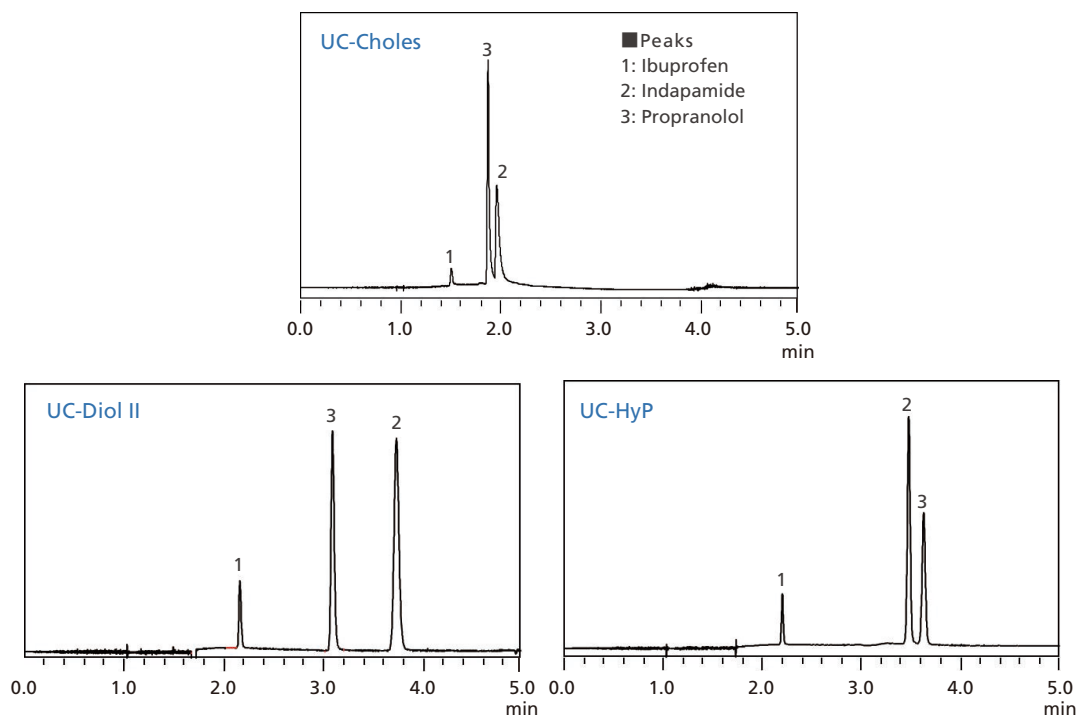


Fig. 27 Differences in Retention Behavior for Respective Stationary Phases (normal-phase for Diol, normal-phase + static electric interaction for HyP, and hydrophobic interaction for Choles)



### 3-1-2 Shim-pack™ UC Series Columns

Shimadzu offers the Shim-pack UC series of columns designed specifically for SFC analysis. The line of Shim-pack UC series columns is listed in Table 4.

For more details about Shim-pack UC series columns, refer to the corresponding brochure (C190-E251).

The set of six columns with different separation selectivity characteristics is especially useful when scouting for columns during method development (Fig. 28)

#### Shim-pack UC Series

[Product >>](#)

Table 4 Shim-pack UC Series Columns

Column	Carbon Loading	Surface Area	Pore Size	Bonded Phase
Shim-pack UC RP	9%	450 m <sup>2</sup> /g	10 nm	(Polar Embedded) Octadecyl group
Shim-pack UC GIS II	11%			Octadecyl group
Shim-pack UC Phenyl	9.5%			Phenyl group
Shim-pack UC CN	14%			Cyanopropyl group
Shim-pack UC Diol	20%			Diol group
Shim-pack UC Sil	-			(Unmodified silica)
Shim-pack UC Amide	18%			Carbamoyl group
Shim-pack UC NH <sub>2</sub>	8%			Aminopropyl group
Shim-pack UC-ODS	16%			3 um : 340 m <sup>2</sup> /g 5 um : 300 m <sup>2</sup> /g
Shim-pack UC Diol II	N.D.	Diol group		
Shim-pack UC Sil II	-	(Unmodified silica)		
Shim-pack UC Choles	20%	Cholesteryl Group		
Shim-pack UC PBr	8%	Pentabromobenzyl Group		
Shim-pack UC PyE	18%	Pyrenylethyl Group		
Shim-pack UC HyP	N.D.	3-Hydroxyphenyl Group		
Shim-pack UC Py		Pyridinyl Group		
Shim-pack UC Triazole		Triazolyl group		
Shim-pack UC-NaE	11%	Naphtylethyl group		
Shim-pack UC-PolyBT	N.D.	N.D.	N.D.	Polybutylene terephthalate (coated on silica gel)
Shim-pack UC-PolyVP				Poly(4-vinylpyridine) group



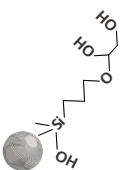
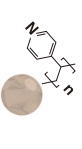
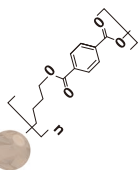
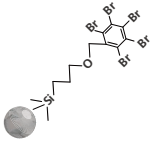
	6 columns set					
	UC-ODS	UC-Sil II	UC-Diol II	UC-PolyVP	UC-PolyBT	UC-PBr
Chemistry						
Features	The separation mode is reverse phase. Retention is provided through hydrophobic interaction.	This is excellent for retention of basic compounds and recognition of their tertiary structures.	The separation mode is normal phase. This inhibits non-specific interactions.	A favorable peak shape is obtained even without acid-base additives.	This is excellent for resolving aromatic compounds through π-π interactions.	With ODS, separation of poorly retained compounds is improved.

Fig. 28 Set of Six Shim-pack UC Series Columns

## 3-2 SFC Mobile Phases

When determining the mobile phase for HPLC, not only the solvent but also the acid/base level, the addition of ion-pair reagents or other modifiers, and other factors that help separation must be considered. For SFC, separation is affected by the quantity and type of modifier added, the acid or other additives added to the modifier, and the pressure and temperature applied to maintain the CO<sub>2</sub> in a fluid state. CO<sub>2</sub> is typically used as the mobile phase for SFC, but there are also cases of supercritical fluoroform being used.<sup>6</sup>

### 3-2-1 Modifiers

Due to the nonpolarity of supercritical carbon dioxide, a modifier such as methanol, ethanol, or acetonitrile can be added to provide polarity to the mobile phase. Increasing the modifier quantity could break down the supercritical fluid state of the CO<sub>2</sub> and potentially change it to a subcritical or liquid state. Nevertheless, it will continue to function as a mobile phase regardless of the state and is not a problem for chromatography. Adding a modifier increases the critical temperature and critical pressure necessary for maintaining supercritical carbon dioxide; hence, adding high percentages of a modifier is not desirable. For HPLC, reversed-phase chromatography and normal-phase chromatography have symmetric properties that can counteract the ability of the mobile phase solvent to elute components. For SFC, the supercritical carbon dioxide always has a low elution capability and the modifier has a high elution capability for both reserved-phase and normal-phase chromatography. Table 5 shows examples of solvents used as modifiers. Fig. 28 shows a comparison of chromatograms separated using different types of modifiers.

Table 5 Examples of Solvents Used as Modifiers

Organic Solvent	Log P <sub>ow</sub>	Organic Solvent	Log P <sub>ow</sub>
Hexane	3.9	Acetone	- 0.24
Toluen	2.69	Dioxane	- 0.27
Tetrahydrofuran	0.47	Ethanol	- 0.32
1-propanol	0.329	Acetonitrile	- 0.34
2-propanol	0.05	Methanol	- 0.82

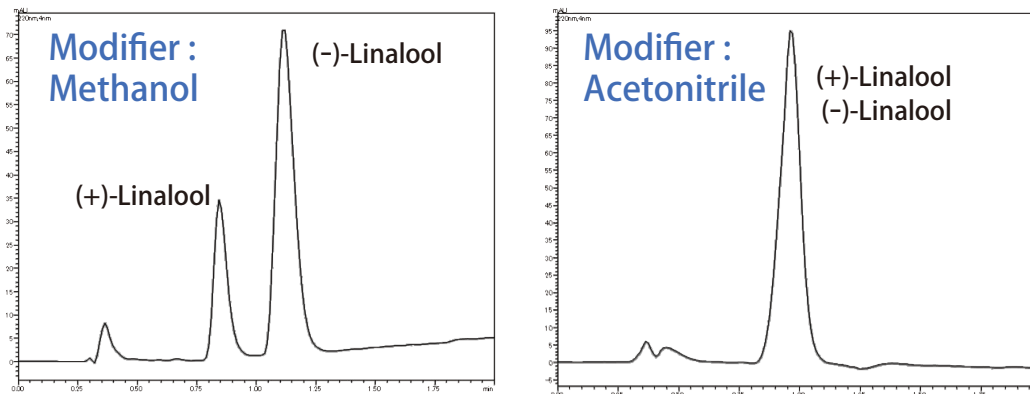


Fig. 29 Comparison of Modifiers





### 3-2-2 Using Additives to Adjust Separation and/or Improve Peak Shape

For HPLC analysis, buffer solutions, ion-pair reagents, or other additives are used to change separation selectivity or improve peak shape. Additives are used for the same reasons in SFC as well.

Additives used for SFC are added to the modifier rather than the supercritical carbon dioxide. The main additives used are listed in Table 6. Adding an acid such as acetic acid or a base such as an amine can improve peak shape by inhibiting ionization of target components or by masking secondary functional groups in the stationary phase (Fig.30). For LC-MS detection, a volatile salt such as ammonium formate is often used because ion suppression can occur. It has been reported that adding a small amount of water ionizes part of the supercritical carbon dioxide so that it behaves like an acidic additive.\*7


Table 6 Examples of Additives and Applicable Compounds

<b>Acid</b>	Acetic acid, Formic acid, Trifluoro acetic acid, Citric acid	<b>Ion-pair reagent</b>	Alkylbenzoic sulfonate salts Tetrabutyl ammonium salts
<b>Base</b>	Triethylamine	<b>Volatile salt</b>	Ammonium formate Ammonium acetate
<b>Other</b>	Water		

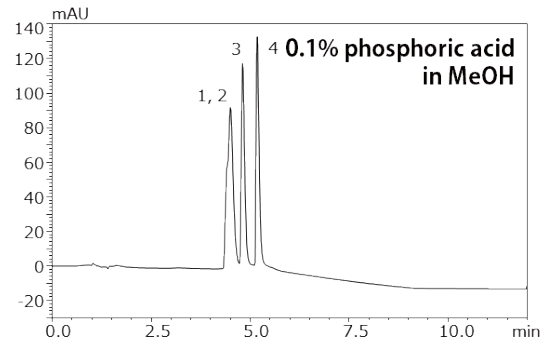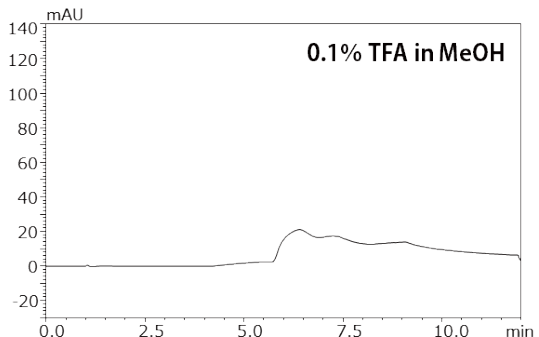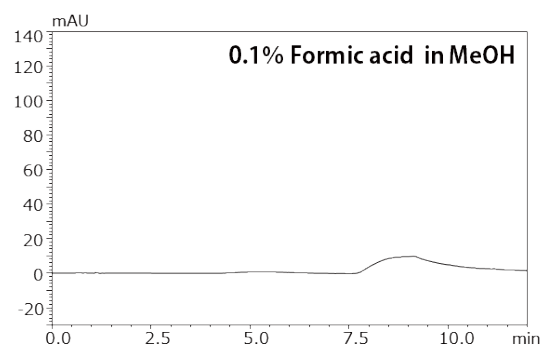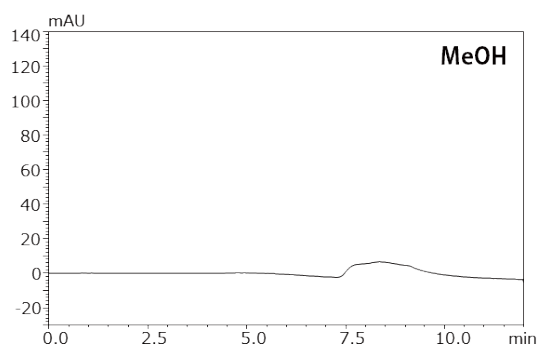


#### ■ Peaks


1. Hydrolyzed s-BPDA

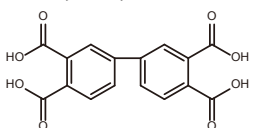

2. Hydrolyzed ODPA-C

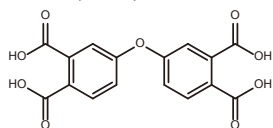

3. Hydrolyzed a-BPDA

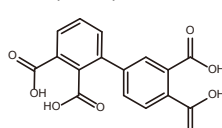

4. Hydrolyzed a-ODPA

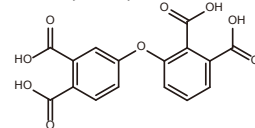


Fig. 30 Effect on Peak Shape of Salt or Acid Type Added to Modifier

### 3-2-3 Other Factors Affecting Separation

For HPLC, the column temperature affects separation only slightly, but for SFC, the temperature and pressure can affect the change in CO<sub>2</sub> state, which changes the density, diffusion coefficient, and other properties of the CO<sub>2</sub>. Therefore, the temperature and pressure setting values are factors that can potentially affect separation (Fig. 31).

As explained in 1-3-3 regarding the van Deemter equation, the large diffusion coefficient makes it easier for supercritical carbon dioxide to penetrate other substances, which results in a lower C value that is related to the mass transfer diffusion and enables higher separation efficiency. Increasing the temperature increases the diffusion coefficient and decreases the viscosity, which enables using a longer column and measuring samples at higher flow rates. Temperature and pressure mainly contribute to the number of theoretical plates, but the number of theoretical plates contributes to resolution by a factor equivalent to its square root. That means it does not have a large effect compared to modifiers, stationary phases, or other parameters that contribute to the selectivity of separation.

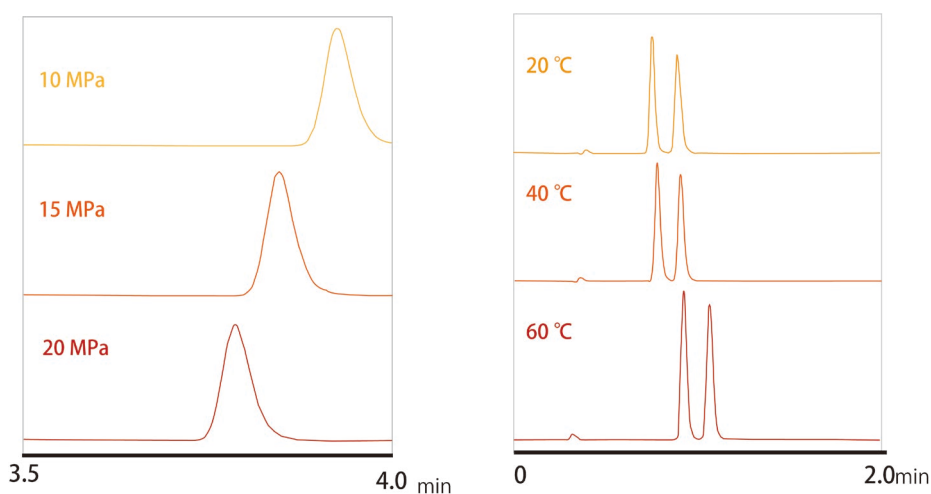


Fig. 31 Effects of Pressure and Temperature on SFC Analysis

### 3-3 Support for SFC Method Development: Method Scouting Systems

For typical HPLC or SFC analysis, scouting (searching) for optimal separation parameters requires a time-consuming analytical condition optimization process of manually adjusting solutions, replacing columns, and specifying gradient profile settings. In particular, because predicting the retention mode for SFC is more difficult than for regular HPLC, separation parameter scouting must be performed comprehensively.

The separation parameter optimization stage requires a process of checking retention behavior as the type and concentration of solvents and additives are varied. Due to the time and trouble involved in the modifier preparation process, it often results in discarding a greater quantity of unwanted solution.

The Nexera UC/s system described on page 14 includes an internal valve for automatically switching between columns and mobile phases/modifiers (Fig. 15) and includes functionality for automatically blending mobile phases to vary the types of solvents, and the types and concentrations of additives used. By using the valve to switch between columns with multiple separation modes, it can comprehensively scout separation parameters. Fig. 32 shows the mobile phase blending process used to automatically prepare modifiers.

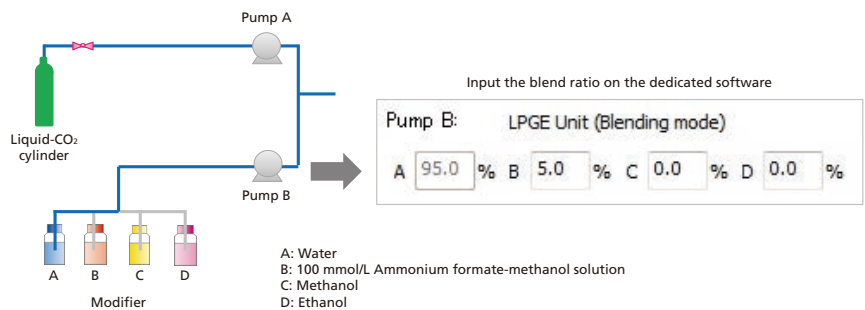


Fig. 32 Effects of Pressure and Temperature on SFC Analysis

Using dedicated Method Scouting Solution software, separation parameters can be scouted comprehensively for multiple analytical conditions (Fig. 33). Based on the chromatogram obtained, optional software is used to evaluate separation and rank the optimized parameters (Fig. 34). The software automatically identifies all chromatograms with resolution greater than a specified criterion (1.5 in this case) and then ranks the resolution in those chromatograms.

For more information about using a Nexera UC/s system to optimize separation parameters, refer to Technical Report (C190-E212) "Improving Separation and Method Development Efficiency Using the Nexera UC/s UHPLC/SFC Switching System"

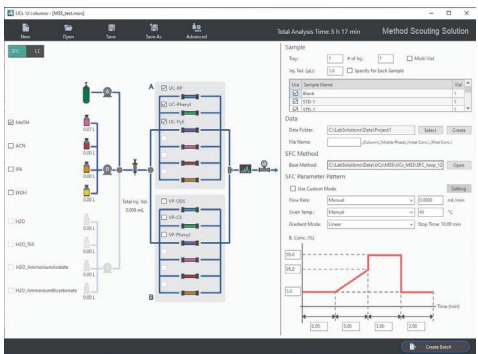


Fig. 33 Method Scouting Solution Operating Window for Nexera UC/s

Ranking	Run No.	Analytical Condition	Resolution	Separation factor	Symmetry factor		Retention factor		Area%		Peak number
					Peak1	Peak2	Peak1	Peak2	Peak1	Peak2	
1	32	Omeprazole_OZ-3_MeOH_20_40	7.965	1.921	1.16	1.159	6.583	12.644	49.829	50.171	2
2	17	Omeprazole_IC-3_MeOH_20_40	5.587	1.602	1.387	1.274	8.078	12.937	49.971	50.029	2
3	16	Omeprazole_IC-3_EtOH_20_40	5.382	1.639	1.915	1.661	8.617	14.124	49.984	50.016	2
4	31	Omeprazole_OZ-3_EtOH_20_40	5.377	1.599	1.169	1.162	7.229	11.561	49.778	50.222	2
5	1	Omeprazole_AD-3_EtOH_20_40	3.996	1.509	1.257	1.404	8.779	13.25	50.054	49.946	2
6	8	Omeprazole_AY-3_MeOH_20_40	3.55	2.08	1.178	1.145	3.652	7.597	49.974	50.026	2
7	11	Omeprazole_IA-3_MeOH_20_40	3.428	1.523	1.464	1.312	7.435	11.327	49.973	50.027	2
8	4	Omeprazole_AS-3_EtOH_20_40	2.515	1.673	1.657	1.518	1.244	2.081	49.754	50.246	2
9	10	Omeprazole_IA-3_EtOH_20_40	1.586	1.157	1.322	1.279	7.115	8.234	49.347	50.653	2

Fig. 34 Using Software to Rank Separation Parameters

## Chapter 4 Examples of SFC Solutions

Due to the properties of supercritical carbon dioxide, SFC not only offers a solution for separating substances that are difficult to separate by HPLC, but also can be used for pretreatment, as described above. This chapter describes examples of using SFC and the SFE pretreatment method.

- 4-1 Analyzing Triglycerides by SFC
- 4-2 Separating Chiral Compounds by SFC
- 4-3 Using Offline SFE to Quickly Extract Compounds from Natural Products
- 4-4 Fast Analysis of Ciclesonide Inhalant by Online SFE-SFC-QTOFMS
- 4-5 Preparative Purification of Volatile Aroma Components by SFC

### 4-1 Analyzing Triglycerides by SFC

Triglycerides are a type of neutral fat decomposed from fat and other tissue to supply energy to cells when needed. Triglycerides are a family of structurally similar compounds with low polarity, which makes them considered difficult to analyze by conventional HPLC. In this example, SFC was used to analyze the triglycerides contained in vegetable oil. For more details about this example, refer to Application News No. C189 "Analysis of Triglycerides Using the Nexera UC Supercritical Fluid Chromatograph"

The composition of triglycerides targeted in this analysis is shown in Table 7. For separation modes that use a retention mechanism based on hydrophobic interactions, such as used in ODS columns, a longer carbon chain length in the acyl group will result in greater retention, whereas more double bonds, given the identical carbon chain length, will result in less retention. Fig. 35 shows a chromatogram from analyzing fish oil diluted with hexane by 10,000 times. The food-grade oil contains many different triglycerides, which comprise different types of acyls depending on their source. In reality, an extremely large number of molecular types exist.

Table 7 Target Components and MRM

Compounds	Acyl composition		MRM	
Trilaurin	TG 36:0	C12:0/C12:0/C12:0	ESI(positive)	656.60>439.40
Trimyristin	TG 42:0	C14:0/C14:0/C14:0	ESI(positive)	740.70>495.45
Tripalmitolein	TG 48:3	C16:1/C16:1/C16:1	ESI(positive)	824.75>551.50
Tripalmitin	TG 48:0	C16:0/C16:0/C16:0	ESI(positive)	818.70>547.45
Trilinolenin	TG 54:9	C18:3/C18:3/C18:3	ESI(positive)	908.85>607.55
Trilinolein	TG 54:6	C18:2/C18:2/C18:2	ESI(positive)	902.80>603.55
Triolein	TG 54:3	C18:1/C18:1/C18:1	ESI(positive)	896.75>599.50
Tristearin	TG 54:0	C18:0/C18:0/C18:0	ESI(positive)	890.70>595.45
Triecosenoin	TG 60:3	C20:1/C20:1/C20:1	ESI(positive)	986.90>659.60
Triarachidin	TG 60:0	C20:0/C20:0/C20:0	ESI(positive)	992.95>663.65
Trierucin	TG 66:3	C22:1/C22:1/C22:1	ESI(positive)	1071.00>715.65
Tribehenin	TG 66:0	C22:0/C22:0/C22:0	ESI(positive)	1077.05>719.70

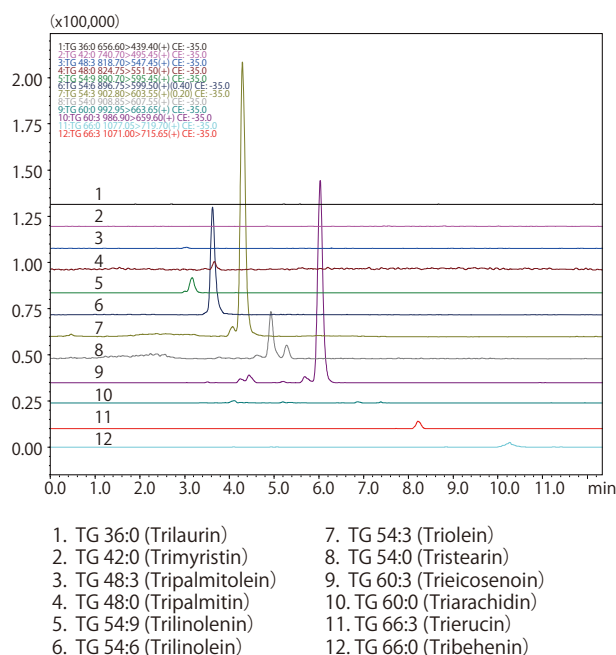


Fig. 35 Chromatogram of Fish Oil

#### Application News

Analysis of Triglycerides Using the Nexera UC Supercritical Fluid Chromatograph

[Application News >>](#)

## 4-2 Separating a Chiral API Molecule using SFC

Optically active (chiral) compounds contain asymmetric carbon structures in their molecules that cannot be superimposed on their mirror images. Such isomers have identical chemical and physical properties, but often have been confirmed to have activity differences, especially in biological organisms. Consequently, there is a growing need to analyze them during drug synthesis for drug discovery in the pharmaceutical industry.

HPLC is the main method used to separate chiral compounds by chromatography, but there are concerns about the cost and environmental impact of the organic solvents used for separation by normal-phase chromatography. The supercritical carbon dioxide used in SFC is nonpolar and has similar properties to n-hexane, which is typically used as a mobile phase in normal-phase chromatography. That makes it possible to migrate chiral separation methods from conventional normal-phase chromatography to SFC. However, that requires selecting from a wide variety of columns and modifiers during optimization of chiral compound separation parameters for SFC. Therefore, the following is an example of using a Nexera UC chiral screening system to determine the optimal chiral separation parameters for omeprazole. For more information about the analysis, refer to Application News No. L495 “Automated Optimization of Chiral Separation Parameters Using Nexera UC Chiral Screening System”

The Nexera UC chiral screening system can automatically switch between up to twelve columns, four types of modifiers, and various mixture ratios to automatically optimize from a large number of possible separation parameters. Each parameter setting is specified using the Method Scouting Solution software shown in Fig. 33. The columns used are listed in Table 8.

Table 8 Chiral Columns Used for Separation Parameter Optimization

Column	Stationary Phase
CHIRALPAK® IA-3/SFC (IA)	Amylose tris (3, 5-dimethylphenylcarbamate)
CHIRALPAK IB-3/SFC (IB)	Cellulose tris (3, 5-dimethylphenylcarbamate)
CHIRALPAK IC-3/SFC (IC)	Cellulose tris (3, 5-dichlorophenylcarbamate)
CHIRALPAK ID-3/SFC (ID)	Amylose tris (3-chlorophenylcarbamate)
CHIRALPAK IE-3/SFC (IE)	Amylose tris (3, 5-dichlorophenylcarbamate)
CHIRALPAK IF-3/SFC (IF)	Amylose tris (3-chloro-4-methylphenylcarbamate)
CHIRALPAK AD-3/SFC (AD)	Amylose tris (3, 5-dimethylphenylcarbamate)
CHIRALPAK AS-3/SFC (AS)	Amylose tris [(S)- <i>α</i> -methylbenzylcarbamate]
CHIRALPAK AY-3/SFC (AY)	Amylose tris (5-chloro-2-methylphenylcarbamate)
CHIRALCEL® OD-3/SFC (OD)	Cellulose tris (3,5-dimethylphenylcarbamate)
CHIRALCEL OJ-3/SFC (OJ)	Cellulose tris (4-methylbenzoate)
CHIRALCEL OZ-3/SFC (OZ)	Cellulose tris (3-chloro-4-methylphenylcarbamate)

Fig. 36 shows the results of separation patterns from a total of 36 possible combinations of 12 chiral columns and 3 types of modifiers (methanol, ethanol, and acetonitrile/ethanol mixture). Optional software was then used to rank the results, with chromatograms for the top three shown in Fig. 37. As a result, the parameters were successfully optimized while reducing the amount of work involved in the tedious process of optimizing separation parameters for chiral analysis by SFC.



Fig. 36 Chiral Screening Separation Results  
Modifiers: (a) Methanol, (b) Ethanol, and (c) Acetonitrile/ethanol = 75/25 (v/v)

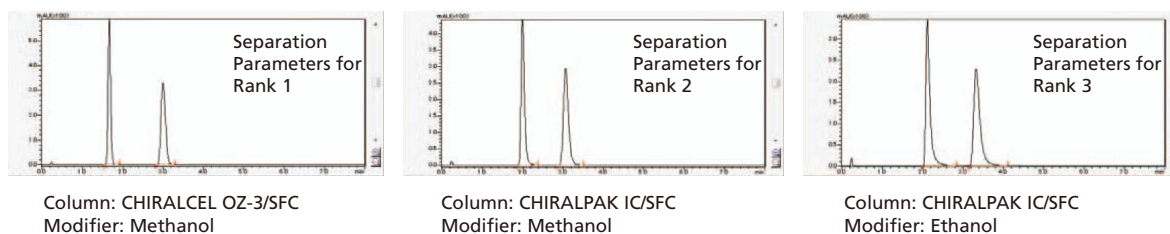


Fig. 37 Chromatograms for the Top Three Separation Parameters

## Application News

### Automated Optimization of Chiral Separation Parameters Using Nexera UC Chiral Screening System

Optically active (chiral) compounds contain asymmetric carbon structures in their molecules that cannot be superimposed on their mirror images. HPLC has been the main method used to chromatographically separate such chiral compounds, but in recent years, the use of supercritical fluid chromatography (SFC) has gained attention. This article describes an example of evaluating conditions for separating and analyzing the chiral compound omeprazole. In this example, a Nexera UC chiral screening system with method scouting functionality was used with 12 chiral analytical columns and three types of modifier solutions (methanol, ethanol, and a mixture of acetonitrile and ethanol).

[Application News >>](#)



## 4-3 Using Offline SFE to Quickly Extract Compounds from Natural Products

Natural products contain components that offer functional benefits used in pharmaceutical ingredients, essential oils, and supplements. Conventionally, such components have been extracted using solvent-based methods, but the extraction accuracy varies depending on the proficiency of the person performing the extraction. In this example, compounds were extracted from natural products by supercritical fluid extraction using a Nexera UC offline SFE system and corresponding fractions were analyzed by HPLC. For more information about this analysis, refer to Application News "Rapid Extraction of Various Compounds from Natural Products."

The samples included tea leaves, ginger, and nutmeg. The samples were homogenized, mixed with a dehydrating agent, packed in an extraction vessel, and loaded into the extraction system (Fig. 38). The system used is shown in Fig. 23 in 2-5 Offline SFE-SFC System. The extraction status was monitored with a PDA detector installed after the extraction system. The resulting chromatograms are shown in Fig. 39. Fig. 40 shows the results from HPLC analysis of the fractions indicated shaded red in Fig. 39. The results indicate that a wide variety of compounds were detected in fractions obtained from respective natural products by SFE.

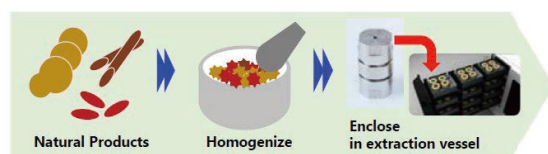


Fig. 38 Pretreatment Process Flow

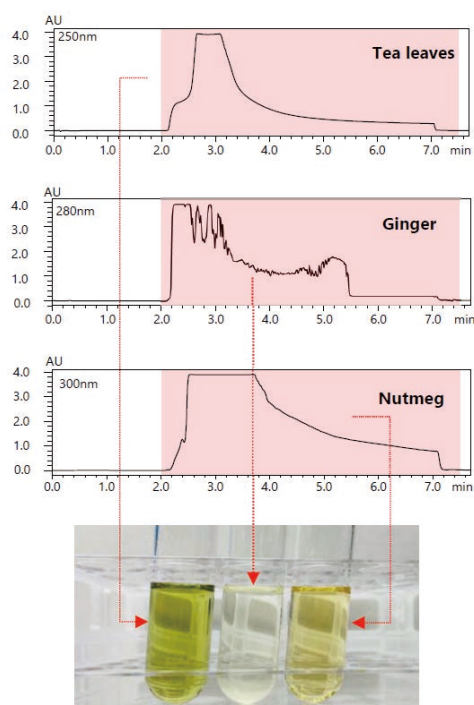


Fig. 39 Extracting Compounds from Natural Products

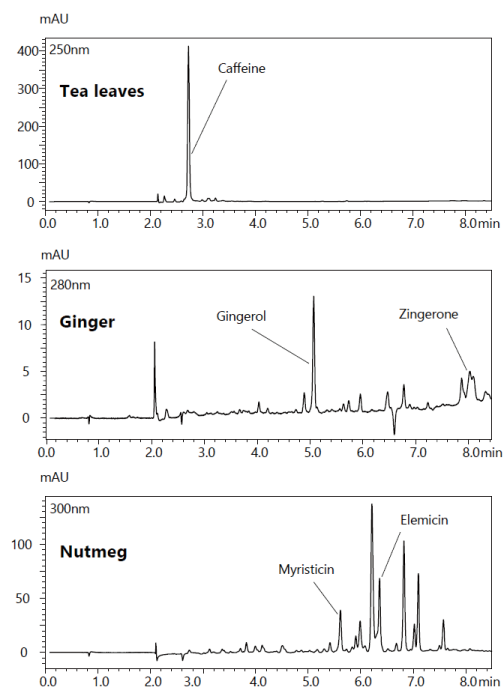


Fig. 40 HPLC Chromatograms of Fractions Extracted from Natural Products

### Application News

#### Rapid Extraction of Various Compounds from Natural Products

[Application News >>](#)

## 4-4 Fast Analysis of Ciclesonide Inhalant by Online SFE-SFC-QTOFMS

Because checking the content of active ingredients and impurities is essential for maintaining pharmaceutical quality, establishing a quick and convenient screening method would help shorten development periods and promote the use of the pharmaceuticals in actual practice.

In this example, screening was performed for active ingredients and impurities contained in the drug Alvesco, a ciclesonide inhalant used to treat bronchial asthma. For more information about this analysis, refer to Application News “Fast Analysis of Ciclesonide Inhalant by Online SFE-SFC-QTOFMS.”

In this analysis, components targeted include ciclesonide and related impurities A, B, and C listed in the European Pharmacopoeia (EP), related impurity IF1 indicated in the Alvesco pharmaceutical interview form (IF), and possible impurities benzothiazole (BT) and 2-mercaptobenzothiazole (MBT) that might be contained in vials.

Samples prepared by spraying Alvesco once onto a glass disc were sealed inside the extraction vessel and measured with an online SFE-SFC-QTOFMS system. Fig. 41 shows extracted ion chromatograms for Alvesco-related compounds, with accurate masses of detected peaks listed in Table 9. Ciclesonide, impurities A, B, and C, BT, and MBT were identified based on retention times and accurate mass values measured from the respective standard mixture samples. Impurity IF2, impurity IF3 (with two peaks presumably from diastereomer), and CH-CIC, which were not contained in the standard mixture sample, were identified based on accurate mass values. Impurities IF1 and IF4 were not detected.

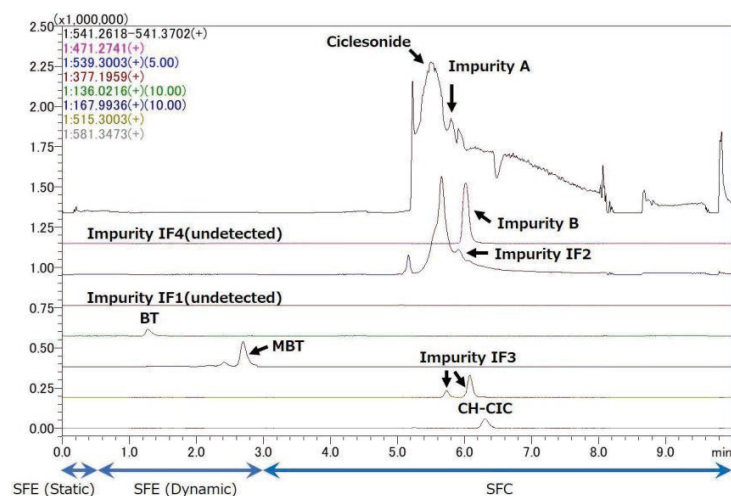


Fig. 41 Extracted Ion Chromatograms from Alvesco

Table 9 Compounds Detected from Alvesco by Online SFE-SFC-QTOFMS

Compound	Formula and Theoretical $m/z$ ( $[M+H]^+$ )	Observed $m/z$	Error (mDa)	RT (min)	Listed in	
					EP	IF
Ciclesonide	$C_{32}H_{44}O_7$ , $m/z$ 541.3160	541.3159	-0.1	5.48	✓	✓
Impurity A	$C_{32}H_{44}O_7$ , $m/z$ 541.3160	541.3158	-0.2	5.87	✓	✓
Impurity B	$C_{28}H_{38}O_6$ , $m/z$ 471.2741	471.2735	-0.6	6.06	✓	✓
Impurity C	$C_{32}H_{42}O_7$ , $m/z$ 539.3003	539.2993	-1.0	5.54	✓	
Impurity IF2	$C_{32}H_{42}O_7$ , $m/z$ 539.3003	539.2995	-0.8	5.96		✓
Impurity IF3 (2 diastereomers)	$C_{30}H_{42}O_7$ , $m/z$ 515.3003	515.2998	-0.5	5.74		✓
		515.3000	-0.3	6.08		✓
CH-CIC	$C_{35}H_{48}O_7$ , $m/z$ 581.3473	581.3463	-1.0	6.33		
BT	$C_7H_5NS$ , $m/z$ 136.0216	136.0211	-0.5	1.27		
MBT	$C_7H_5NS_2$ , $m/z$ 167.9936	167.9932	-0.4	2.70		

### Application News

Fast Analysis of Ciclesonide Inhalant by Online SFE-SFC-QTOFMS

[Application News >>](#)

## 4-5 Preparative Purification of Volatile Aroma Components by SFC

Conventionally, gas chromatography (GC) has been used to fractionate aroma components and other volatile compounds. Using GC for fractionation provides high resolution, but it involves issues with low sample loads per analysis and long analysis times. Fractionation by supercritical fluid chromatography (SFC) enables sample loads equivalent to liquid chromatography (LC) and faster analysis times. The liquified carbon dioxide gas that is used as the mobile phase for SFC evaporates at ambient temperatures and pressures, leaving only a small quantity of organic solvent in fractions after fractionation. That makes it easy to concentrate target components. This example describes using a Nexera UC supercritical fluid chromatograph analytical fraction system (refer to 2-3-3 Analytical Fraction System) for preparative separation of linalool as the target component, which is the key aroma component in lavender essential oil. For more information about this analysis, refer to Application News “Preparative Purification of Aroma Components Using a Supercritical Fluid Chromatograph.”

The analytical fraction system was used to separate the fractions containing a commercial linalool product (Fig. 42). A sample solution containing about 20 g/L of the target components was used. Then the (+)-linalool and (-)-linalool fractions obtained were used to determine the recovery rate of each component. That resulted in a high recovery rate of 97 % or more for both target components (Table 10).

The analytical fraction system was used to separate the fractions containing linalool in lavender essential oil (Fig. 43). The commercial lavender essential oil was diluted with methanol by 10 fold and filtered through a 0.2  $\mu$ m membrane filter before use. Then the (+)-linalool and (-)-linalool fractions obtained were reanalyzed with the analytical fraction system to confirm their purity. The resulting chromatograms from reanalysis are shown in Fig. 44 and purity confirmation results are indicated in Table 11. Fractionation using the analytical fraction system achieved high-purity preparative purification with over 99 % purity.

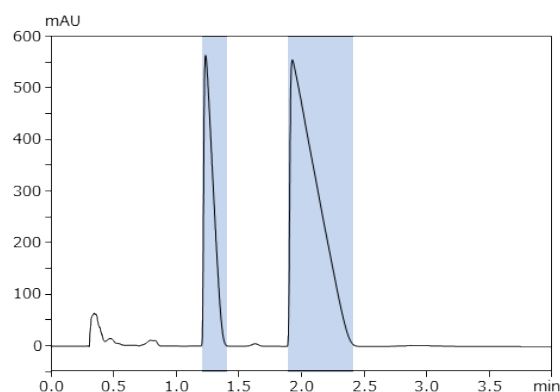


Fig. 42 Preparative Chromatogram of Linalool

Table 10 Recovery Rate of Each Fraction

	Recovery rate (%)
(+)-Linalool	97.6
(-)-Linalool	99.3

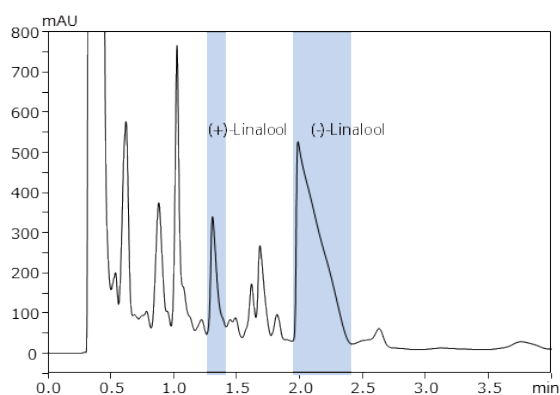


Fig. 43 Preparative Chromatogram of Lavender Essential Oil

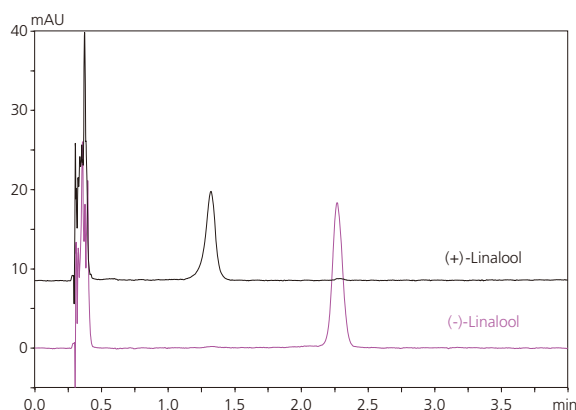


Fig. 44 Chromatograms from Reanalysis of Collected Fractions

Table 11 Purity of Target Components Contained in Fractions (Area Percent of Peaks in Detection Range 0.5 to 4.0 min)

	Area (%)
(+)-Linalool	100%
(-)-Linalool	99.4%

### Application News

#### Preparative Purification of Aroma Components Using a Supercritical Fluid Chromatograph

[Application News >>](#)

## Chapter 5 Frequently Asked Questions about SFC

### 5-1 SFC Installation

**Q1: Are there any special precautions for indoor use?**

Be careful to discharge exhaust gases into a fume hood to prevent CO<sub>2</sub> concentrations from becoming too high within the laboratory.

Installing a CO<sub>2</sub> meter is also recommended for the purpose of detecting CO<sub>2</sub> gas leaks within the room.

**Q2: Is a temperature-controlled (and humidity-controlled) room environment preferable for maintaining the carbon dioxide gas supply temperature?**

As long as the instruments are used within the manufacturer's recommended conditions, controlling the room temperature is probably not a significant concern because temperature is controlled in the column oven, but gas cylinders are affected by the ambient temperature where they are installed so care should be given to that temperature.

**Q3: Can an LC unit currently in use (Prominence™ LC-20A series) be used in combination with SFC?**

Yes, if the current Prominence unit is operated in combination with a modifier pump. For more details, contact Shimadzu.

**Q4: Can CO<sub>2</sub> cylinders typically distributed be used?**

Cylinders of 99.9 % CO<sub>2</sub> gas that are typically sold are more than acceptable for use. However, quality can vary so be sure to contact the gas vendor to confirm detailed specifications.

**Q5: About how much CO<sub>2</sub> gas is used per analysis?**

Assuming ten minutes per analysis and pumping 100 % CO<sub>2</sub> at a flow rate of 3 mL/min, about 30 g is used per analysis. Use that value to recalculate usage based on appropriate changes to the analysis time, modifier ratio (which reduces the proportion of CO<sub>2</sub>), and flow rate.

## 5-2 Samples and Target Components

### Q1: Do samples analyzed by SFC need to be dissolved in a liquid like HPLC?

In general, samples must be dissolved in a solution before injection, just like for HPLC. However, by connecting it to a supercritical fluid extraction (SFE) system, even components in solid samples can be analyzed.

### Q2: For SFC, how is sample solubility in the mobile phase checked?

Due to the polarity of supercritical carbon dioxide, the ideal mobile phase should be n-hexane or a mixture of n-heptane and 2-propanol, but typically samples are dissolved for use in 100 % methanol or another modifier.

### Q3: Can samples containing water be injected?

Yes, they can. However, for larger injection volumes, the sample solvent could affect peak shape or retention times.

### Q4: Does supercritical carbon dioxide cause any sample decomposition?

Supercritical carbon dioxide is known to cause slight acidity in mobile phases, but there is minimal risk of decomposing or otherwise altering compounds.

### Q5: Are there any examples of analyzing blood, urine, or other biological samples by injecting them directly into an SFC system?

We are not aware of any examples of analyzing blood by direct injection. However, there are examples of analyzing blood soaked in filter paper by SFE-SFC and examples of directly injecting urine samples for SFC-MS analysis.<sup>\*8-11</sup>

### Q6: Are there any samples that can be analyzed by HPLC but not SFC?

SFC analysis of amino acids, soluble vitamins, polyphenols, and other highly polar compounds is possible by increasing the modifier concentration and adding an acid. However, the resulting peak shapes are not as good as when using HPLC. For preparative applications, SFC offers the advantage of not requiring drying or purification steps after fractionation.

### Q7: What type of modifier solvent is typically selected for analyzing compounds that decompose in alcohol or water?

A non-protic solvent is often used in cases where SFC is useful for analyzing samples that decompose or change properties when exposed to moisture. Candidate non-protic solvents include acetonitrile and THF.

## 5-3 Mobile Phase

### Q1: What waste liquids remain after analysis? Are they the solvent with dissolved measurement target substances, modifiers, and make-up solution?

Yes, waste liquids include the solvent containing dissolved measurement target substances, modifiers, and make-up solutions. However, liquid waste bottles will also contain carbon dioxide gas and some solvent consumption occurs due to evaporation.

### Q2: What additives can be added to modifiers to improve peak broadening problems when measuring hydrophilic peptides?

There is a report of analysis by using TFA as a method for improving peak shape in peptide analysis by SFC.<sup>\*12</sup>

## 5-4 Columns

### Q1: Can HPLC columns also be used for SFC?

SFC columns are generally compatible with fittings used in Shimadzu HPLC systems so it is possible to connect them. However, confirm the pressure capacity and other requirements for actual use with respective column manufacturers.

### Q2: Can columns be used for both HPLC and SFC?

Yes, they can. However, it requires particular care regarding purging column liquids after analyses are finished. After using the column for SFC analysis, reduce the pressure inside the column to vaporize any CO<sub>2</sub> in the column or purge it with methanol or ethanol before disconnecting the column. Disconnecting the column with only supercritical carbon dioxide dries out the packing material due to evaporation of liquid in the column. Then if that column is used for HPLC, the high viscosity of the solvent used for the LC mobile phase, such as methanol or water, will cause surface tension around pores in the stationary phase that prevents penetration into the pores and interaction with functional groups. Because the supercritical fluid used for SFC analysis has lower viscosity than the mobile phase solvents used for HPLC, they can penetrate inside pores and interact with functional groups without being affected by surface tension, even if not filled with solvent. If a column used for HPLC analysis is used for SFC, any water in the column must be thoroughly purged with methanol or ethanol, because the water will not mix with the supercritical carbon dioxide.

### Q3: Is there a difference in the service life of columns used for SFC and HPLC?

Due to the lower viscosity of supercritical fluid, it applies less force on the packing material. In other words, it applies less back pressure on the column, which tends to result in longer column life than for HPLC.

### Q4: Since we cannot predict the elution order for separation by SFC, we do not know which columns are to be used or how they should be used. When screening columns for SFC, what separation parameters should be considered?

For SFC, the same mobile phase parameters can be used for everything from normal-phase columns to reversed-phase columns. Therefore, we recommend first fixing the mobile phase parameters and then checking separation using multiple columns. For advice about selecting columns, refer to 3-1-2 Column Selection. Given the multiple columns, solvent delivery parameters, and separation parameters such as temperature and pressure, there is a huge number of possible parameter combinations. To reduce the amount of work involved in considering separation parameters, you are recommended to use a method scouting system because it automatically sets and executes respective parameter combinations.

## 5-5 Miscellaneous

### Q1: Are retention time and peak intensity repeatability levels comparable to HPLC levels in terms of precision and accuracy?

Repeatability comparisons for HPLC and SFC using identical conditions have resulted in comparable repeatability and quantitative accuracy levels.

### Q2: We have been struggling with matrix and ghost peaking problems in our LC/MS analysis. Will the same phenomena occur for SFC as well?

Ionization can be inhibited by the matrix even for SFC/MS. However, since LC and SFC can result in significantly different separation behavior, sometimes SFC/MS can resolve such LC/MS matrix effect problems.

### Q3: Can SFE be used for samples that contain non-volatile salts?

Yes, as long as the non-volatile salts are not dissolved by the extraction solvent. However, beware of contamination and precipitation if using a mass spectrometer for detection.

## References

\*1

Simultaneous analysis for water- and fat-soluble vitamins by a novel single chromatography technique unifying supercritical fluid chromatography and liquid chromatography

J Chromatogr.A 2014 Oct3; 1362:270-7

\*2

Overview of the retention in subcritical fluid chromatography with varied polarity stationary phases

J Sep Sci. 2008 May; 31(8):1238-51

\*3

The dependence of reduced plate height on reduced velocity in carbon dioxide supercritical fluid chromatography with packed columns

Chromatographia 1987 Jan; 21-25

\*4

Technical Report "Supercritical Fluid Chromatography" (C190-E200)

\*5

Comparison of Retention Behavior between Supercritical Fluid Chromatography and Normal-Phase High-Performance Liquid Chromatography with Various Stationary Phases

Molecules 2019 Jul; 24(13): 2425

\*6

Comparison of Reversed-Phase HPLC Separation Using Carbon Dioxide and Fluoroform for Enhanced-Fluidity Liquid Mobile Phases Anal.Chem.1998,70,1595-1603

\*7

Packed column supercritical fluid chromatography of hydrophilic analytes via water-rich modifiers

Journal of Chromatography A, Volume 1250, 10 Aug 2012, P196-204

\*8

Carotenoids and apocarotenoids determination in intact human blood samples by online supercritical fluid extraction-supercritical fluid chromatography-tandem mass spectrometry

Analytica Chimica Acta , Volume 1032, 22 Nov 2018, P40-47

\*9

High-throughput phospholipid profiling system based on supercritical fluid extraction-supercritical fluid chromatography/mass spectrometry for dried plasma spot analysis Journal of Chromatography A, Volume 1250, 16 Jun 2012, P69-75

\*10

High-throughput phospholipid profiling system based on supercritical fluid extraction-supercritical fluid chromatography/mass spectrometry for dried plasma spot analysis Journal of Chromatography B, Volume 1076, 10 Aug 2012, P196-204

\*11

Use of on-line supercritical fluid extraction-supercritical fluid chromatography/tandem mass spectrometry to analyze disease biomarkers in dried serum spots compared with serum analysis using liquid chromatography/tandem mass spectrometry

Rapid Communications in Mass Spectrometry, Volume 31, 30 May 2017, P886-894

\*12

Advantageous use of SFC for separation of crude therapeutic peptides and peptide libraries

J Pharm Biomed Anal. ,Volume 185, 5 Jun 2020

Nexera, LotusStream, Shim-pack and Prominence are trademarks of Shimadzu Corporation or its affiliated companies in Japan and/or other countries.  
CHIRALPAK and CHIRALCEL are registered trademarks of Daicel Corporation.



Shimadzu Corporation  
[www.shimadzu.com/an/](http://www.shimadzu.com/an/)

**For Research Use Only. Not for use in diagnostic procedures.**

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

Company names, products/service names and logos used in this publication are trademarks and trade names of Shimadzu Corporation, its subsidiaries or its affiliates, whether or not they are used with trademark symbol "TM" or "®".

Third-party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®".

Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

The contents of this publication are provided to you "as is" without warranty of any kind, and are subject to change without notice. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication.