

# Sensitive detection and quantification of hydrogen sulfide as a gasotransmitter by combining monobromobimane-based derivatization to triple quadrupole LC/MS/MS

**ASMS 2012** WP25-512

Takako Hishiki <sup>1</sup>, Tsuyoshi Nakanishi <sup>2</sup> and Makoto Suematsu <sup>1,3</sup>

<sup>1</sup> Department of Biochemistry, School of Medicine, Keio University, Tokyo, Japan

<sup>2</sup> MS Business Unit, Shimadzu Corporation, Kyoto, Japan

<sup>3</sup> Japan Science and Technology Agency, Exploratory Research for Advanced Technology, Suematsu Gas Biology Project, Tokyo, Japan

# Sensitive detection and quantification of hydrogen sulfide as a gasotransmitter by combining monobromobimane-based derivatization to triple quadrupole LC/MS/MS

## Overview

Hydrogen sulfide (H<sub>2</sub>S) regulates various cellular functions as a gasotransmitter and reliable and accurate quantification of H<sub>2</sub>S concentration in biological samples is required. Here the quantitative analysis of endogenous H<sub>2</sub>S level in tissue samples was performed by coupling the monobromobimane-based assay to triple quadrupole

LC/MS/MS. Furthermore, this method has an advantage of the measurement at neutral pH to quantify biologically free-gaseous H<sub>2</sub>S not acid-labile H<sub>2</sub>S and this approach is expected to be a useful method to quantify free H<sub>2</sub>S level in tissue samples.


## Introduction

Hydrogen sulfide (H<sub>2</sub>S) is known to be a third gaseous mediator and also to have important roles on a vasodilation in both the peripheral and cerebral circulation. To elucidate various signal transduction pathways elicited by H<sub>2</sub>S, reliable and accurate measurement of H<sub>2</sub>S concentration in biological samples becomes increasingly significant. Until now, head space gas analysis using gas chromatography and measurement by polarographic sensor have been

reported to quantify H<sub>2</sub>S concentration in biological samples. However these approaches often yielded a higher value resulting from acid-labile H<sub>2</sub>S not from biologically free-gaseous H<sub>2</sub>S. Here we developed sensitive quantitative analysis of H<sub>2</sub>S in biological samples, mouse brain tissues by combining a monobromobimane-based derivatization at a neutral pH to SRM (selected reaction monitoring) using triple quadrupole LC/MS/MS.

## Methods

Thiol-specific derivatization agent monobromobimane was used to quantify the endogenous H<sub>2</sub>S in brain tissues. Frozen brains from 12-d-old mice were placed in two volumes of ice-cold 5 mM monobromobimane (mBBr) in 10 mM Tris-HCl (pH 7.5) and then homogenized. Methanol was added to precipitate proteins. The mixture was vortexed and centrifuged. A supernatant was desalted by a solid-phase extraction and filtrated through a 5-kDa cutoff filter. An aliquot of 10 µL was analyzed to quantify sulfide-dibimane (SDB), monobromobimane-based derivatives of H<sub>2</sub>S by LC/MS/MS instrument, Nexera UHPLC system and LCMS-8030 triple quadrupole mass spectrometer. SDB was eluted from an ODS column (150 mmL. X 2.1 mmI.D. 1.7 µm particle size) with a gradient of acetonitrile, detected in negative mode ESI / SRM mode.

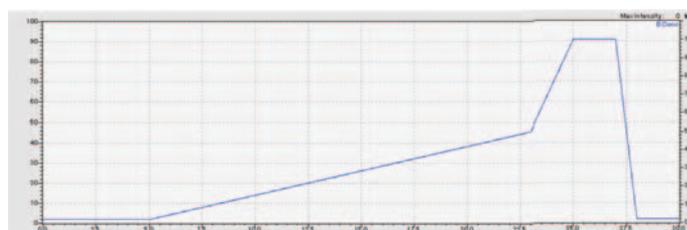

Column	: commercially available ODS column (150 mmL. X 2.1 mmI.D. 1.7 mm particle size)
Flow rate	: 0.2 mL/min
Column temp.	: 40°C
Mobile phase	: A) water containing 0.1% formic acid : B) acetonitrile containing 0.1% formic acid
Time prog.	: 2%B (0 → 5min) → 50%B (23 min) → 100%B (25 → 27 min) → 2%B (28 → 30 min)
Ionization mode	: ESI negative
DL temp.	: 250°C
HB temp.	: 400°C
Nebulizing gas flow	: 3.0 L/min
Drying gas flow	: 15 L/min
Transition	: <i>m/z</i> 413.10 → <i>m/z</i> 191.15 (Q1 : 30 V, CE : 20 V, Q3 : 20 V)

Fig. 1 Analytical condition of SRM to quantify endogenous hydrogen sulfide

# Sensitive detection and quantification of hydrogen sulfide as a gasotransmitter by combining monobromobimane-based derivatization to triple quadrupole LC/MS/MS

## Results

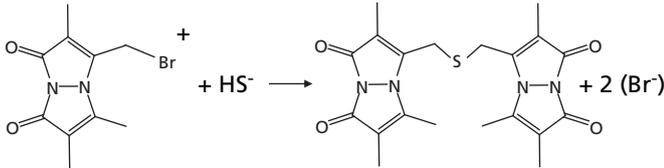


Fig. 2 Scheme of the derivatization of hydrogen sulfide with monobromobimane

Hydrogen sulfide is derivatized with monobromobimane reagent and forms sulfide-dibimane (SDB). Reacted SDB is stable at room temperature and is measured on SRM analysis. In this study, the transition of SDB was determined at  $m/z$  413  $\rightarrow$  191 at CE of 20 V in ESI/negative mode.

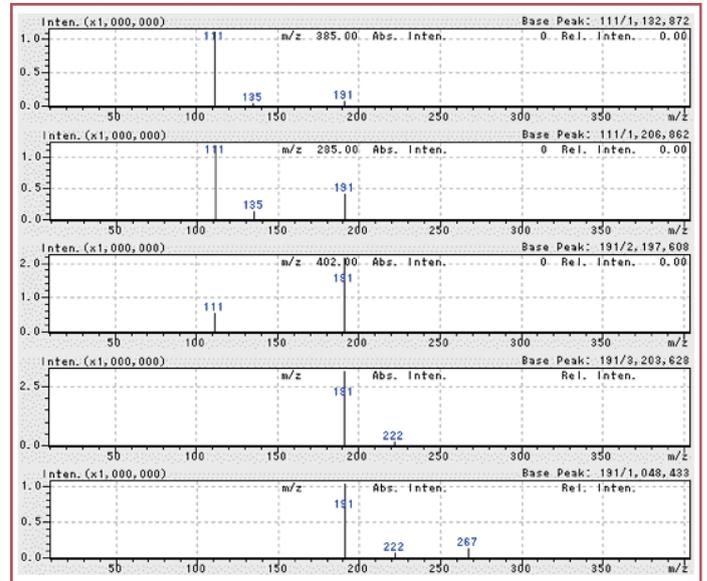


Fig. 3 MS/MS spectra of SDB standard on various collision energy (CE)

(A) A calibration curve of SDB standard.

A calibration curve was plotted by dilution series of SDB standard and were linear at a range of 1-5000 nM ( $R^2 > 0.99$ ).

(B) A calibration curve of NaHS solution forming H<sub>2</sub>S

A calibration curve was also plotted by using the derivatization with mBBr to H<sub>2</sub>S resulting from NaHS solution. Various concentrations of NaHS were reacted with excess monobromobimane (1 mM) for 10 min on ice.

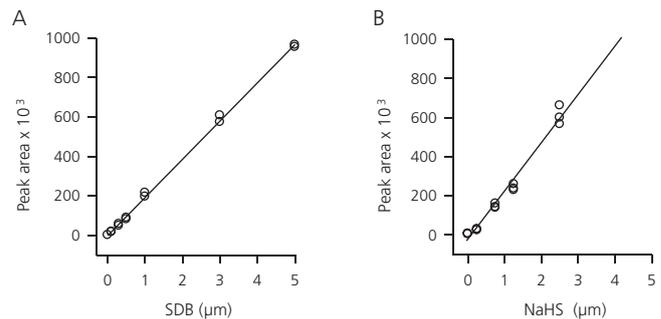


Fig. 4 Calibration curves of SDB and NaHS standard


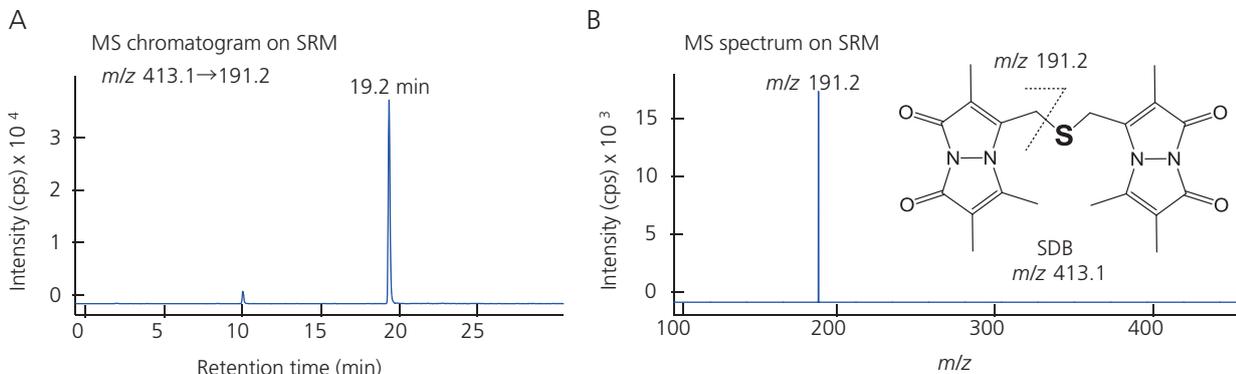


Fig. 5 MS chromatogram and MS spectrum on SRM of SDB derived from H<sub>2</sub>S in murine brain tissues

## Sensitive detection and quantification of hydrogen sulfide as a gasotransmitter by combining monobromobimane-based derivatization to triple quadrupole LC/MS/MS

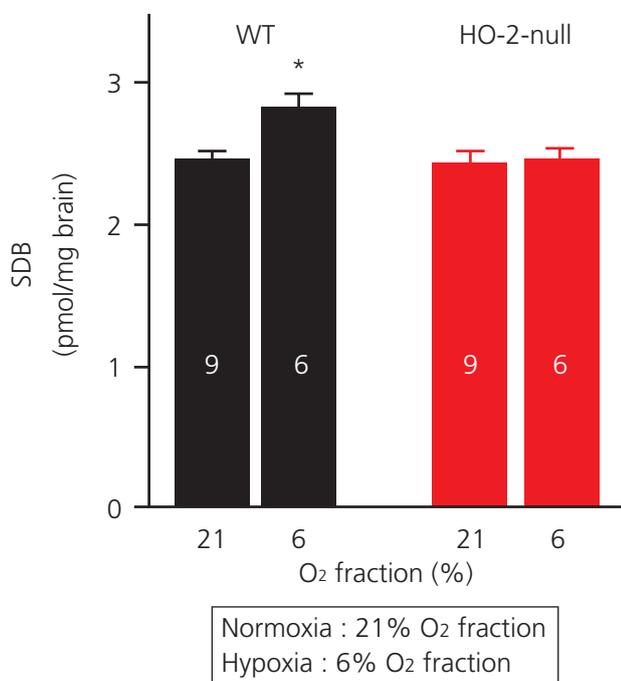

Endogenous H<sub>2</sub>S concentration in neonatal mouse brain was determined by SRM analysis using a monobromobimane-based assay. To quantify H<sub>2</sub>S level in brain tissues, a calibration curves was plotted by a SDB standard and applied as a external calibration. Fig. 6 shows H<sub>2</sub>S concentration in WT mouse increased under 6% O<sub>2</sub> fraction (hypoxia). On the contrary, the alteration of H<sub>2</sub>S concentration from heme oxygenase (HO)-2-null mouse brains could not be found between 21% and 6% fraction. Under normoxia, H<sub>2</sub>S levels were similar in WT and HO-2-null mice, presumably due to compensation from other sources of H<sub>2</sub>S in the transsulfuration pathway.

Error bar show a standard error (n=9 or 6).  
\*P<0.05 compared with WT at 21 % O<sub>2</sub>.

Fig. 6 Endogenous H<sub>2</sub>S concentration in neonatal brain of both WT and HO-2-null under normoxia/hypoxia

## Conclusions

Hydrogen sulfide is known to be as the gasotransmitter to control multiple biological functions. Therefore reliable and accurate quantitative analysis of H<sub>2</sub>S level in tissues is required. In this study we performed the quantitative analysis of hydrogen sulfide in mouse brain tissue by coupling the monobromobimane-based assay to triple quadrupole LC/MS/MS.

First analytical condition of SDB was determined on SRM by a SDB standard and then the linearity of calibration curves was validated by using a SDB standard and NaHS reagent forming H<sub>2</sub>S. Correlation coefficient of these calibration

curve was confirmed to be >0.99 and

The endogenous H<sub>2</sub>S level in mouse brain tissues was quantified as the derivatized SDB by triple quadrupole LC/MS/MS. This result shows the of the endogenous H<sub>2</sub>S is at the level of pmol/mg tissue. And the increase of the H<sub>2</sub>S level in WT mouse was confirmed under the hypoxia. On the other hand, the change of H<sub>2</sub>S level could not be found in the HO-2-null mouse. These results suggests this quantitative analysis by using triple quadrupole LC/MS/MS coupled with monobromobimane-based assay is helpful to quantify the endogenous H<sub>2</sub>S level in tissue samples.

## References

(1) Morikawa T, et al. (2012) Hypoxic regulation of the cerebral microcirculation is mediated by a carbon monoxide-sensitive hydrogen sulfide pathway. *Proc. Natl. Acad. Sci. U.S.A.* 109:1293-1298



Shimadzu Corporation  
[www.shimadzu.com/an/](http://www.shimadzu.com/an/)

For Research Use Only. Not for use in diagnostic procedures.  
The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.