

# Application News

## No. A591

### Spectrophotometric Analysis

## Evaluation of Fluorescence Anisotropy in Fluorescein and Rhodamine B Solutions

The evaluation of fluorescence anisotropy can be used for research into physical and chemical molecular properties, for example, as a method of investigating the electronic structure and solvent-solute interactions of small molecules as well as the conformation and dynamics of natural and synthetic macromolecules<sup>1)</sup>.

This article introduces the results of investigating the changes in fluorescence anisotropy of fluorescein solution and rhodamine B solution while varying the temperature of these solutions.

K. Maruyama, K. Sobue

### Fluorescence Anisotropy

Measurement of fluorescence anisotropy in solution is published in the referenced IUPAC Technical Report<sup>1)</sup>. The basic relations defined in this report are given below.

Fluorescence anisotropy ( $r$ ) is generally defined by the following equation.

$$r = \frac{(I_{VV} - I_{VH})}{(I_{VV} + 2I_{VH})} \quad (1)$$

where  $I_V$  is the intensity measured with the excitation polarizer and emission polarizer in the vertical position (V), and  $I_H$  is the intensity measured with the excitation polarizer in the vertical position (V) and the emission polarizer in the horizontal position (H). The denominator is the total fluorescence intensity and the numerator is the degree of polarization.

The fluorescence polarization ratio ( $p$ ) is generally defined by the following equation.

$$p = \frac{(I_{VV} - I_{VH})}{(I_{VV} + I_{VH})} \quad (2)$$

The relationship between fluorescence anisotropy and fluorescence polarization ratio is given by the following expression.

$$r = \frac{2p}{(3 - p)} \quad (3)$$

Generally, the fluorescence anisotropy is expressed using the angle  $\alpha$  between the absorption transition moment upon molecular excitation and the emission transition moment upon fluorescence emission, in the following way.

$$r = 0.2(3 \cos^2 \alpha - 1) \quad (4)$$

Equation (4) shows that the fluorescence anisotropy will be in the range of  $-0.2 \leq r \leq 0.4$  (and decreases due to molecular rotation in the time between absorption and emission).

Note that in actual measurement, sensitivity differences between the polarized components (vertical and horizontal directions) of the light emitted from the instrument must be corrected. The G-factor, which is defined by the following equation, is used for correction.

$$G = \frac{I_{HV}}{I_{HH}} \quad (5)$$

In actual measurement, the fluorescence anisotropy ( $r_{\text{exp}}$ ) is calculated using the following equation which incorporates the G-factor into equation (1).

$$r_{\text{exp}} = \frac{(I_{VV} - GI_{VH})}{(I_{VV} + 2GI_{VH})} \quad (6)$$


### Temperature Changes and Fluorescence Anisotropy

Fluorescein was dissolved into glycerol solution (glycerin and pure water prepared to a 95:5 ratio) to obtain a concentration of 0.2 mg/L. We determined the optimal excitation wavelength using 3D spectrum measurement, and performed measurement while varying the solution temperature under the conditions listed in Table 1. Figs. 1 and 2 and Table 2 show the measurement results.

As the temperature rises, energy dissipation, internal conversion, and intersystem crossing due to collisions between molecules occur more readily leading to a decrease in fluorescence intensity (temperature quenching)<sup>2)</sup>. Fig. 1 illustrates the decrease in fluorescence intensity as the temperature rises.

**Table 1 Measurement Conditions**


|                               |  |
|-------------------------------|--|
| Instrument used               | : RF-6000, polarized light measurement attachment, constant-temperature single cell holder |
| Excitation wavelength         | : 497 nm (fluorescein solution)<br>555 nm (rhodamine B solution)                           |
| Fluorescence wavelength range | : 505 nm to 650 nm (fluorescein solution)<br>565 nm to 700 nm (rhodamine B solution)       |
| Scan speed                    | : 600 nm/min   |
| Sampling interval             | : 1.0 nm   |
| Bandwidth                     | : Ex 3.0 nm/Em 5.0 nm  |
| Sensitivity                   | : Low  |

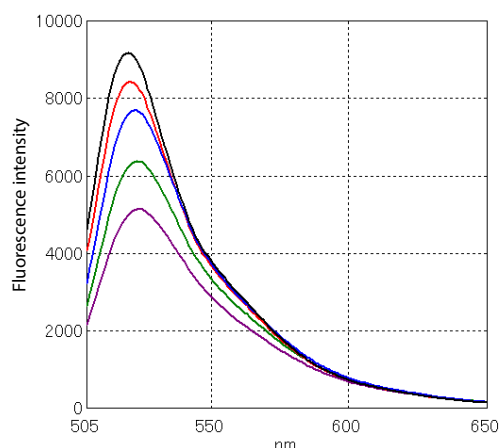

**Fig. 1 Fluorescence Spectra of Fluorescein Solution (Without Polarizer)**

**Black: 10 °C, Red: Room Temperature, Blue: 40 °C, Green: 60 °C, Purple: 78 °C**

As shown by Fig. 2, the intensities required to calculate the fluorescence anisotropy were measured through the combination of excitation and emission polarizers. We can see that the sensitivity of the horizontal component is higher on the instrument used.

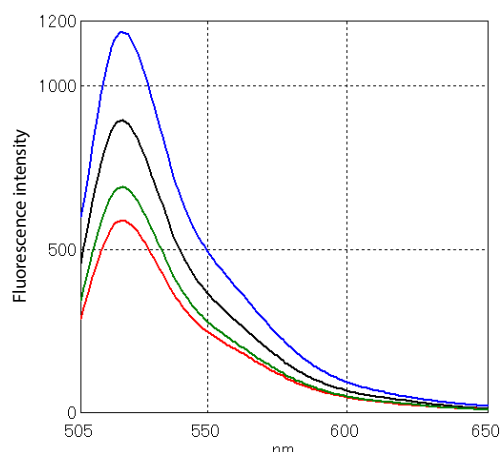


**Fig. 2 Polarization Spectra of Fluorescein Solution (10 °C)**  
Black:  $I_{VV}$ , Red:  $I_{VH}$ , Blue:  $I_{HH}$ , Green:  $I_{HV}$

**Table 2 Fluorescence Anisotropy Values and Fluorescence Polarization Ratios of Fluorescein Solution at Various Temperatures**


| Temperature      | Fluorescence wavelength | $I_{VV}$ | $I_{VH}$ | G    | r    | p     |
|------------------|-------------------------|----------|----------|------|------|-------|
| 10 °C            | 520 nm                  | 931.4    | 588.2    | 0.59 | 0.36 | 0.457 |
| Room temperature | 521 nm                  | 804.6    | 525.2    | 0.59 | 0.35 | 0.444 |
| 40 °C            | 522 nm                  | 644.8    | 501.4    | 0.59 | 0.28 | 0.371 |
| 60 °C            | 523 nm                  | 450.7    | 437.9    | 0.59 | 0.20 | 0.271 |
| 78 °C            | 524 nm                  | 304.1    | 359.3    | 0.59 | 0.13 | 0.178 |

Table 2 shows the decrease in fluorescence anisotropy as the temperature rises. The three causes of depolarization are generally considered to be intramolecular relaxation, rotational Brownian motion, and intermolecular excitation energy transfer<sup>3)</sup>. Intermolecular excitation energy transfer does not need to be taken into account for this sample because fluorescein is the only fluorescent molecule in the solution. Furthermore, since intramolecular relaxation is not really affected by temperature changes, we can infer that the decrease in fluorescence anisotropy is occurring due to the increased rotational Brownian motion as the temperature rises.

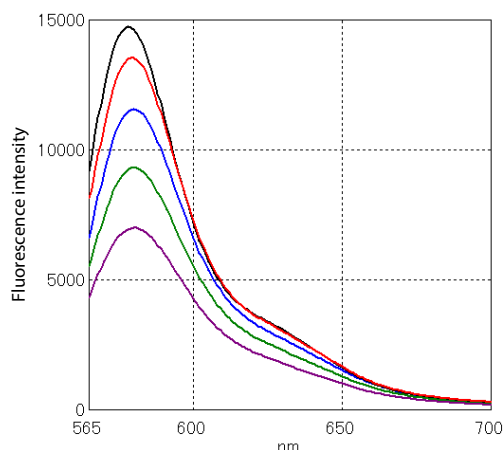

**Fig. 3 Fluorescence Spectra of Rhodamine B Solution (Without Polarizer)**  
Black: 10 °C, Red: Room Temperature, Blue: 40 °C, Green: 60 °C, Purple: 78 °C

Next, rhodamine B was dissolved into glycerol solution (glycerin and pure water containing 0.5 mM NaOH prepared to a 95:5 ratio) to obtain a concentration of 0.2 mg/L. We determined the optimal excitation wavelength using 3D spectrum measurement, and performed measurement while varying the solution temperature under the conditions listed in Table 1. Figs. 3 and 4 and Table 3 show the measurement results.

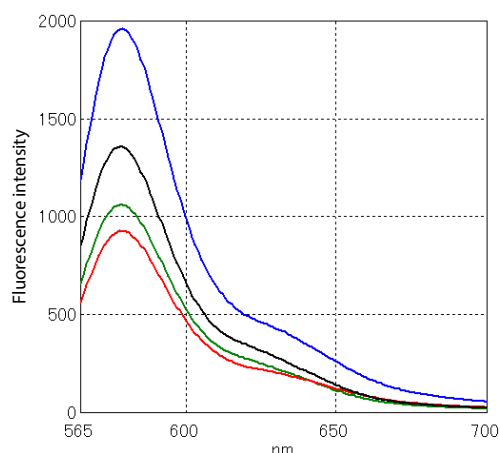

**Fig. 4 Polarization Spectra of Rhodamine B Solution (10 °C)**  
Black:  $I_{VV}$ , Red:  $I_{VH}$ , Blue:  $I_{HH}$ , Green:  $I_{HV}$

**Table 3 Fluorescence Anisotropy Values and Fluorescence Polarization Ratios of Rhodamine B Solution at Various Temperatures**

| Temperature      | Fluorescence wavelength | $I_{VV}$ | $I_{VH}$ | G    | r    | p     |
|------------------|-------------------------|----------|----------|------|------|-------|
| 10 °C            | 578 nm                  | 1359.8   | 928.2    | 0.54 | 0.36 | 0.461 |
| Room temperature | 579 nm                  | 1244.3   | 885.4    | 0.54 | 0.35 | 0.445 |
| 40 °C            | 580 nm                  | 987.1    | 764.5    | 0.54 | 0.32 | 0.410 |
| 60 °C            | 580 nm                  | 710.5    | 627.6    | 0.54 | 0.27 | 0.354 |
| 78 °C            | 580 nm                  | 489.3    | 481.2    | 0.54 | 0.23 | 0.306 |

As with the fluorescein solution, the rhodamine B solution also showed a decrease in fluorescence intensity and fluorescence anisotropy as the temperature rises. We can again infer that the decrease in fluorescence anisotropy is caused by the increased rotational Brownian motion as the temperature rises, which is consistent with the fluorescein solution.

## Conclusion

We measured the fluorescence anisotropy of fluorescein solution and rhodamine B solution while varying the temperature of these solutions. We confirmed that decreases in fluorescence intensity and fluorescence anisotropy coincide with temperature increases. The decrease in fluorescence anisotropy is most likely caused by the increased rotational Brownian motion as the temperature rises. This demonstrates that measuring fluorescence anisotropy helps to investigate information related to physical and chemical molecular properties.

## References

- 1) Marcel Ameloot, Martin vandeVen, A. Ulises Acuna, and Bernard Valeur. Fluorescence anisotropy measurements in solution: Methods and reference materials (IUPAC Technical Report). *Pure Appl. Chem.*, 2013, vol. 85, No. 3, p. 589-608
- 2) Yasuharu Nishikawa, Keizo Hiraki, "Fluorescence and Phosphorescence Spectrometry", Kyoritsu Shuppan (in Japanese)
- 3) Kazuhiko Kinoshita, Koushin Mihashi (Editor), "Fluorometry", Gakkai Publishing Center (in Japanese)

First Edition: Mar. 2019



For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. Shimadzu disclaims any proprietary interest in trademarks and trade names used in this publication other than its own. See <http://www.shimadzu.com/about/trademarks/index.html> for details.

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

Shimadzu Corporation

[www.shimadzu.com/an/](http://www.shimadzu.com/an/)