

Application News

No. M291

Gas Chromatography Mass Spectrometry

Application of Pyrolysis-GC/MS to Food Contamination Analysis

In recent years, cases of food contamination by foreign matter have increased, heightening the necessity of contamination analysis by food manufacturers. The Fourier transform infrared spectrometer (FTIR) and energy dispersive X-ray fluorescent spectrometer (EDX) are used in identification of contaminants by instrumental analysis. The pyrolysis-GC/MS method and thermal extraction-GC/MS method are employed in contamination analysis by GC/MS, and enable qualitative analysis of resin materials and additives contained in trace organic contaminants.

This article introduces the results of an analysis of the resins in a food packaging material using the pyrolysis-GC/MS method, assuming food contamination. A Shimadzu OPTIC-4 multimode inlet for GC/MS was used in the analysis by pyrolysis-GC/MS. Because the OPTIC-4 enables high-speed heating (60 °C/s) to a maximum temperature of 600 °C, diverse sample injection modes are available and simple pyrolysis is possible.

K. Kawamura, Y. Kawakita, Y. Kudo



Fig. 1 GCMS-QP™2020 NX + OPTIC-4

■ Sample and Analysis Conditions

Fig. 1 shows the appearance of the GCMS-QP™2020 NX and OPTIC-4 multimode inlet used in the analysis.

A commercially-available food packaging material was used as the real sample material. The sample material was cut with a cutter knife to obtain a sample weighing approximately 0.2 mg, which was inserted in the DMI microvial of the OPTIC-4 and then set in the DMI insert liner (Fig. 2).

Table 1 shows the instrument system and analysis conditions.

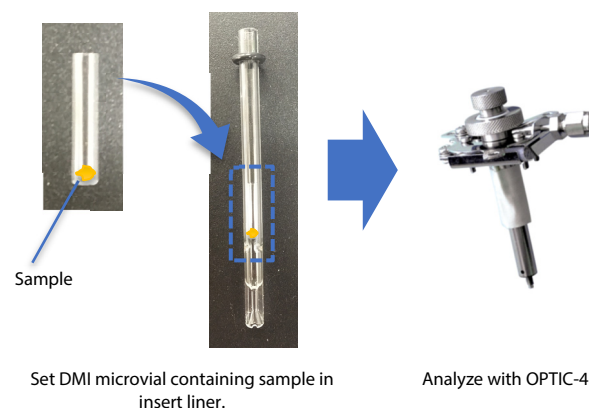


Fig. 2 Sampling Procedure for Pyrolytic Analysis Using OPTIC-4

Table 1 Analysis Conditions

Inlet	: OPTIC-4	
GC/MS	: GCMS-QP2020 NX	
Column	: UA-5 (MS/HT)-30M-0.25F (L30 m, 0.25 mm I.D., df=0.25 µm) (Frontier Laboratories Ltd.)	
OPTIC-4 conditions		
Vent time	: 30 s	
Equilibrium time	: 5 s	
End time	: 40 min	
Injection temperature	: 40 °C (10 s) → (60 °C/s) → 600 °C (3 min) → 320 °C (hold)	
Carrier gas	: Helium	
Carrier control mode	: Flow control	
Start column flow	: 1.0 mL/min	
End column flow	: 1.0 mL/min	
Split flow	: 100 mL/min	
Septum purge flow	: 3 mL/min	
GC conditions		
Column oven temp	: 40 °C (2 min) → 20 °C/min → 320 °C (16 min)	
MS conditions		
Interface temp.	: 280 °C	
Ion source temp.	: 230 °C	
Ionization method	: EI	
Measurement mode	: Scan (m/z 29 - 800)	
Event time	: 0.3 s	

Qualitative Analysis of Resin Material

Fig. 3 shows the obtained pyrogram (total ion chromatogram obtained by pyrolysis-GC/MS). From a reference⁽¹⁾ containing pyrolysis data on resins, this is a distinctive pyrogram of polyethylene (PE), in which

hydrocarbon species are arranged at equal intervals. Therefore, it could be inferred that the foreign matter in this experiment contained PE as the base material.


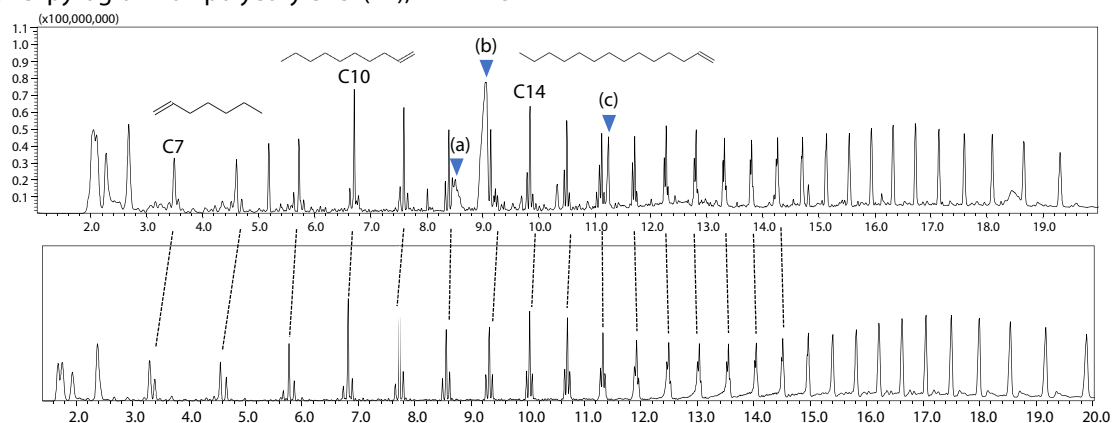



Fig. 3 Top: Pyrogram of Analyzed Foreign Matter, Bottom: Pyrogram of PE

In addition to the peaks seen in the pyrogram of PE, three distinctive peaks ((a) to (c)) were also detected in the pyrogram of the real sample. Compound identification of these peaks was carried out using the NIST Library and the above-mentioned reference⁽¹⁾. As a result, it was found that (b) is caprolactam, a compound characteristically seen as a pyrolysis product of polyamide (PA), and (a) and (c) were identified


respectively as 4-(vinylloxycarbonyl) benzoic acid and benzoic acid, which are compounds characteristically seen as pyrolysis products of polyethylene terephthalate (PET). Based on these results, the foreign matter measured in this experiment was estimated to be a composite resin containing polyamide (PA) and polyethylene terephthalate (PET) in addition to polyethylene (PE).

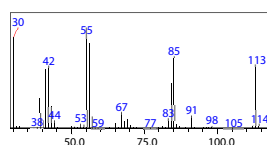



Peak (b): Caprolactam

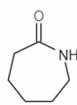


Characteristic component as pyrolysis product of PA

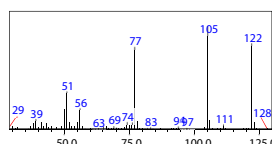

Peak (a): Benzoic acid

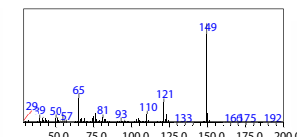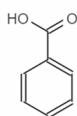

Peak (c): 4-(vinylloxycarbonyl)benzoic acid

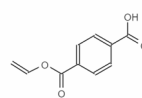

Characteristic components as pyrolysis product of PET

Fig. 4 Mass Spectra of Peaks (a) to (c) and Identified Compounds

Conclusion

As an analysis of contaminants in food products, the resins contained in an assumed foreign matter sample were analyzed by the pyrolysis-GC/MS method in an OPTIC-4 multimode inlet. As a result, qualitative analysis of the composite resin was possible from the pyrogram and pyrolysis products. Thus, this experiment demonstrated

the possibility of qualitative analysis of resin materials by using the pyrolysis-GC/MS method, including analysis of trace contaminants and contaminants in multilayer films, which are difficult to analyze by FTIR. By using this analysis technique, it is considered possible to identify the source of contamination and take appropriate countermeasures.

<Reference>

- (1) S. Tsuge, H. Ohtani, C. Watanabe: Pyrolysis-GC/MS Data Book of Synthetic Polymers – Pyrograms, Thermograms and MS of Pyrolyzers – , 1st Edition, Elsevier, 420 (2011)

GCMS-QP is a trademark of Shimadzu Corporation in Japan and/or other countries.

First Edition: Apr. 2020



Shimadzu Corporation

www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. Shimadzu disclaims any proprietary interest in trademarks and trade names used in this publication other than its own. See <http://www.shimadzu.com/about/trademarks/index.html> for details.

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

© Shimadzu Corporation, 2020

Related Products

Some products may be updated to newer models.



> GCMS-QP2020 NX

Gas Chromatograph Mass Spectrometer

Related Solutions

> Food Contamination

> Price Inquiry

> Product Inquiry

> Technical Service /
Support Inquiry

> Other Inquiry