

Analysis of Additives in Food Containers Using LCMS-9030 Quadrupole Time-of-Flight Liquid Chromatograph Mass Spectrometer

T. Goda

User Benefits

- ◆ Simple qualitative analysis of various types of additives contained in polymer materials is possible by using this system.
- ◆ Highly accurate quantitative analysis is realized by utilizing the high resolution of the LCMS-9030, even with complex samples.
- ◆ The LCMS-9030 demonstrates its effectiveness in evaluations of polymer materials and development/improvement of functional polymers.

Introduction

Additives such as antioxidants, ultraviolet absorbers, and flame retardants are added to polymers, namely, plastics and rubber. Because performance and durability can be dramatically enhanced by adding appropriate types and quantities of additives in the development and production of polymer materials, it is extremely important to obtain information on the additives contained in those materials.

Application News No. C79 introduced an example of quantitative analysis of polymer additives using a triple quadrupole (TQ) type high performance liquid chromatograph mass spectrometer (LC-MS). This Application News introduces an example of qualitative analysis and quantitative analysis of the polymer additives in food containers using the LCMS-9030, which is a quadrupole time-of-flight (QTOF) LC-MS.

Analysis of Food Containers

In this experiment, the polymer additives contained in five types of food containers were analyzed. Samples for the analysis were prepared by adding 1 mL of THF to 0.1 g of finely-chopped food container materials (pack, film), conducting sonication for 1 min, adding 1 mL of methanol to the sample, filtering the supernatant with a 0.2 µm filter, and then diluting the sample with methanol.

Fig. 1 shows the result of peak picking of the food A film using the "ANALYZE" function of LabSolutions Insight Explore™.


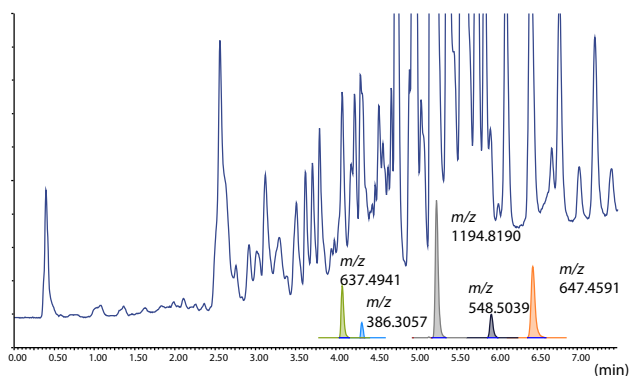


Fig. 1 Result of Peak Picking of Food A Film Using Insight Explore

Peaks were detected on the extracted ion chromatogram (EIC) at m/z 637.4941, 386.3057, 1194.8190, 548.5039, and 647.4591. These m/z correspond to the m/z of ions originating from Irganox® 1098, CYANOX® 425, Irganox® 1010, Irganox® 1076, and Irgafos® 168, respectively.

As an example, Fig. 2 shows the composition estimation results for the peak (peak X) on the EIC at m/z 637.4941. It was found that the compositional formula is $C_{40}H_{64}N_2O_4$.

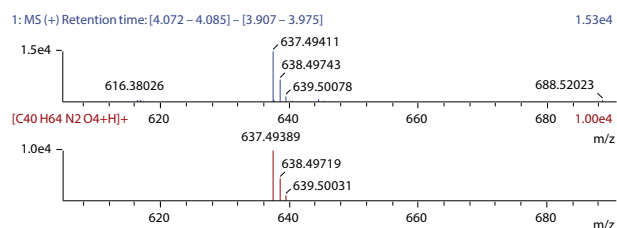

| # | Score | Pred. (M) | Pred. m/z | Meas. m/z | Diff. (mDa) | Formula (M) | Ion | Diff. (ppm) | Iso Score | DBE |
|---|-------|-----------|-----------|-----------|-------------|--|--------------------|-------------|-----------|------|
| 1 | 98.74 | 636.48661 | 637.49389 | 637.49411 | 0.22 | C ₄₀ H ₆₄ N ₂ O ₄ | [M+H] ⁺ | 0.345 | 98.74 | 10.0 |
| 2 | 88.67 | 636.48690 | 637.49418 | 637.49411 | -0.07 | C ₃₉ H ₇₀ O F P S | [M+H] ⁺ | -0.110 | 98.55 | 5.0 |
| 3 | 83.95 | 636.48639 | 637.49367 | 637.49411 | 0.44 | C ₃₈ H ₆₆ N ₂ O F ₂ S | [M+H] ⁺ | 0.690 | 93.61 | 6.0 |
| 4 | 71.89 | 636.48723 | 637.49451 | 637.49411 | -0.40 | C ₃₂ H ₇₀ N ₄ O ₄ P ₂ | [M+H] ⁺ | -0.627 | 69.07 | 1.0 |
| 5 | 64.90 | 636.48729 | 637.49457 | 637.49411 | -0.46 | C ₃₃ H ₆₄ N ₈ O ₂ S | [M+H] ⁺ | -0.722 | 72.46 | 6.0 |
| 6 | 60.34 | 636.48638 | 637.49366 | 637.49411 | 0.45 | C ₃₄ H ₇₅ N ₂ P ₃ S | [M+H] ⁺ | 0.706 | 67.38 | 0.0 |
| 7 | 56.03 | 636.48683 | 637.49410 | 637.49411 | 0.01 | C ₃₃ H ₇₂ N ₄ O F ₃ | [M+H] ⁺ | 0.016 | 73.36 | 0.0 |
| 8 | 54.55 | 636.48672 | 637.49400 | 637.49411 | 0.11 | C ₃₁ H ₆₆ N ₆ O ₄ F ₃ | [M+H] ⁺ | 0.173 | 60.67 | 2.0 |

Fig. 2 Result of Composition Estimation of Peak X (Top: Measured Spectrum, Middle: Theoretical Spectrum, Bottom: Candidates of Compositional Formula)

Table 1 Measurement Conditions

| [HPLC conditions] | |
|---------------------|--|
| Column | :Kinetex® 2.6u XB-C18 (75 mm x 2.1 mm I.D., 2.6 µm) |
| Mobile phase A | :Water containing 10 mmol/L Ammonium formate |
| Mobile phase B | :Methanol |
| Flow rate | :0.5 mL/min |
| Time program | :35 % (0 min) – 100 % (4-7.5 min) – 35 % (7.51-10 min) |
| Column Temp. | :40 °C |
| Injection vol. | :2 µL |
| [MS conditions] | |
| Ionization | :ESI positive / ESI negative |
| Mode | :Scan, MS/MS |
| Nebulizing gas flow | :2.0 L/min |
| Drying gas flow | :10.0 L/min |
| Heating gas flow | :10.0 L/min |
| DL temp. | :250 °C |
| BH temp. | :400 °C |
| Interface temp. | :300 °C |

Compound Search

In addition, an analysis was carried out using ACD/MS Structure ID Suite (Advanced Chemistry Development, Inc., ACD) to confirm the structural formula and compound name of this peak. This software makes it possible to list compounds from the PubChem offline database, which contains approximately 100 million items, based on precise mass and compositional formula information, and assign ranks of compounds from the degree of coincidence (assignment rate) between the product ion obtained by fragment prediction and the product ion observed in the measured MS/MS spectrum.

A database search by the compositional formula $C_{40}H_{64}N_2O_4$ found 71 candidate compounds.

Fig. 3 shows the result of the ranking of the candidate compounds found in the database search. The compound PubChem CID 90004 showed the highest assignment rate. As the result of an online search using PubChem, this compound was identified as Irganox® 1098 (Fig. 4).

Pubchem ID 90004 : Irganox® 1098

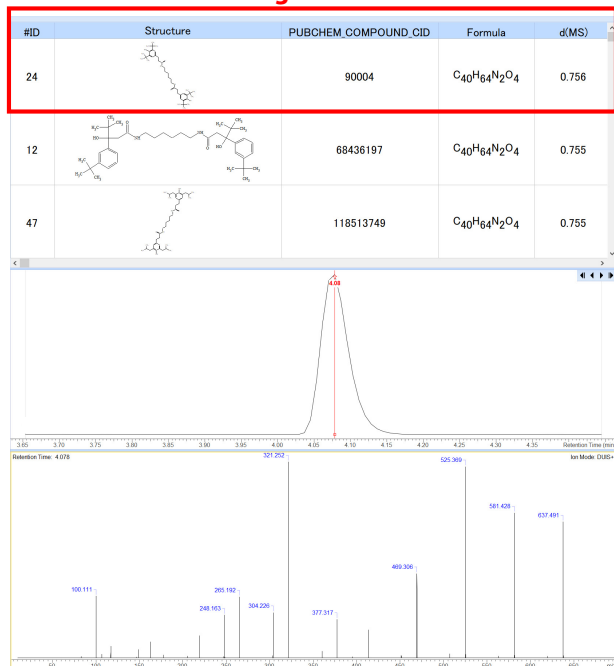


Fig. 3 Results of Database Search and Ranking of Assignment Rates Using ACD/MS Structure ID Suite

For a sure qualitative analysis, a standard sample of the additive was measured. The result of confirmation of the retention time of the above-mentioned Irganox® 1098 coincided with the retention time of peak X. The MS/MS spectra of Irganox® 1098 and peak X were also compared (Fig. 5). These results confirmed that peak X is Irganox® 1098.

As described above, it is possible to search the structural formula and compound name of target peaks by a process of ① LC-MS and LC-MS/MS measurement, ② peak picking using Insight Explore, ③ composition estimation of the target peak, ④ search for the structural formula using an offline database, ⑤ narrowing of the candidate compounds based on the assignment rates of fragments, ⑥ search for the structural formula and compound name using the online database, and ⑦ confirmation analysis using a standard sample using Shimadzu LCMS-9030 and LabSolutions Insight Explore and ACD/MS Structure ID Suite.

PubChem CID: 90004

Structure:

Find Similar Structures

Chemical Safety: [Laboratory Chemical Safety Summary \(LCSS\)](#)
[Datasheet](#)

Molecular Formula: $C_{40}H_{64}N_2O_4$

2.4 Synonyms

2.4.1 Depositor-Supplied Synonyms

23128-74-7
Antioxidant 1098
N,N'-(Hexane-1,6-diyl)bis(3-(3,5-di-tert-butyl-4-hydroxyphenyl)propanamide)
3,3'-bis(3,5-di-tert-butyl-4-hydroxyphenyl)-N,N'-hexamethylenedipropionamide
EINECS 245-442-7
Irganox 1098
UNII-918T54D300
n,n'-hexamethylenebis(3,5-di-tert-butyl-4-hydroxyhydrocinnamamide)
N,N'-Hexane-1,6-diylbis[3-(3,5-di-tert-butyl-4-hydroxyphenyl)propionamide]
Benzenepropanamide, N,N'-1,6-hexanediylbis(3,5-bis(1,1-dimethylethyl)-4-hydroxy-

Fig. 4 Result of Compound Search Using Online Database

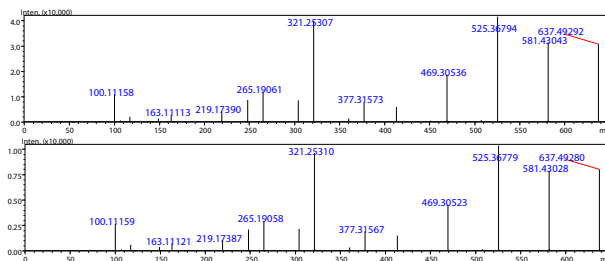


Fig. 5 MS/MS Spectra of (Top) Standard Sample of Irganox® 1098 and (Bottom) Peak X

Quantitative Analysis of Polymer Additives in Food Containers

MS/MS chromatograms were acquired for a quantitative analysis of the polymer additives contained in the food containers.

Fig. 6 shows the respective calibration curves, and Table 2 shows the range of the calibration curve and the coefficient of determination (R^2) of each compound.

Samples of the food containers were extracted by the method described above and were diluted from 10 times to 1000 times with methanol. As a result of the quantitative analysis/calculations, it was found that the concentration of Irgafos® 168 in the samples diluted 1000 times was from 1.85 to 40 ppb, and the concentrations in the packs and films were in the range of 37 to 800 mg/g. Table 3 shows the results of the quantitative analysis of each compound including Irgafos® 168, and Fig. 7 shows representative MS/MS chromatograms of the food A film.

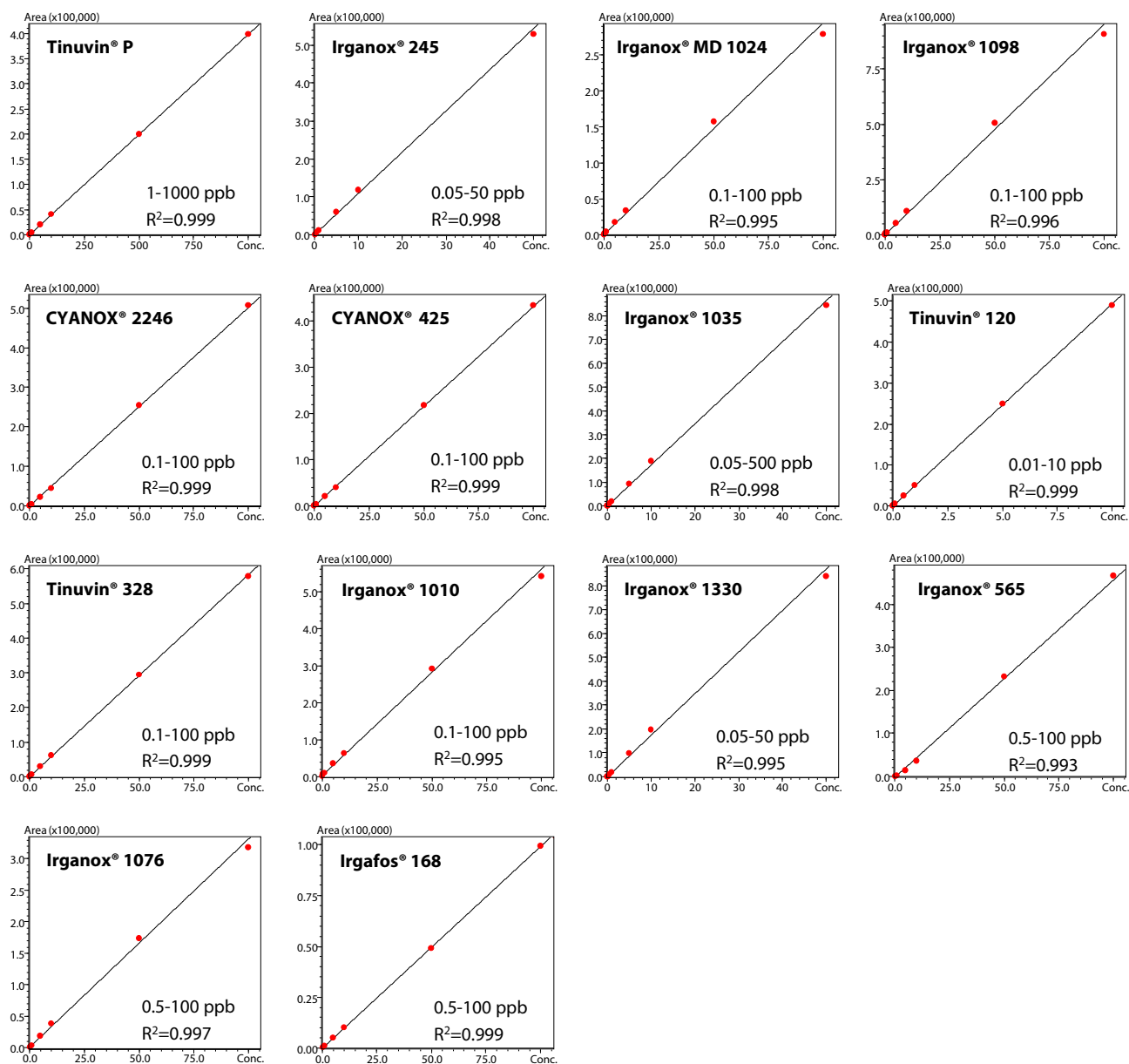


Fig. 6 Calibration Curves of 14 Polymer Additive Components

Table 2 Range of Calibration Curves and Coefficients of Determination of 14 Polymer Additive Components

| Compound name | Ionization method | Precursor ion | Monitor ion | Range of calibration curve (ppb) | Coefficient of determination (R^2) |
|------------------|-------------------|---------------|-------------|----------------------------------|--|
| Tinuvin® P | ESI positive | 226.0975 | 120.0556 | 1-1000 | 0.999 |
| Irganox® 245 | ESI positive | 604.3844 | 177.1279 | 0.05-50 | 0.998 |
| Irganox® MD 1024 | ESI positive | 570.4265 | 181.0972 | 0.1-100 | 0.995 |
| Irganox® 1098 | ESI positive | 637.4939 | 321.2537 | 0.1-100 | 0.996 |
| CYANOX® 2246 | ESI negative | 339.2330 | 163.1128 | 0.1-100 | 0.999 |
| CYANOX® 425 | ESI negative | 367.2643 | 367.2643 | 0.1-100 | 0.999 |
| Irganox® 1035 | ESI positive | 660.4292 | 249.1485 | 0.05-50 | 0.998 |
| Tinuvin® 120 | ESI positive | 439.3207 | 233.1531 | 0.01-10 | 0.999 |
| Tinuvin® 328 | ESI positive | 352.2383 | 282.1601 | 0.1-100 | 0.999 |
| Irganox® 1010 | ESI positive | 1194.8179 | 1194.8179 | 0.1-100 | 0.995 |
| Irganox® 1330 | ESI positive | 792.6289 | 219.1743 | 0.05-50 | 0.995 |
| Irganox® 565 | ESI positive | 589.3968 | 250.1009 | 0.5-100 | 0.993 |
| Irganox® 1076 | ESI positive | 548.5037 | 475.4146 | 0.5-100 | 0.997 |
| Irgafos® 168 | ESI positive | 647.4588 | 647.4588 | 0.5-100 | 0.999 |

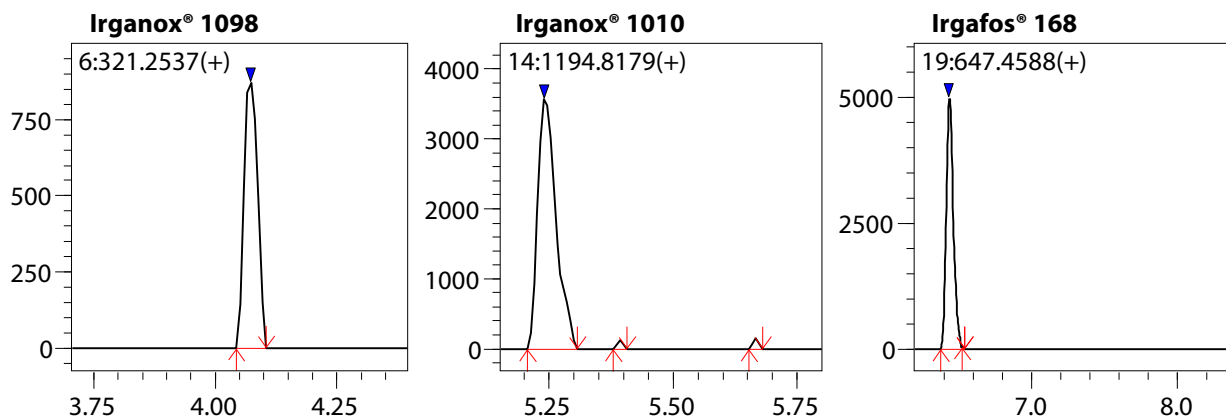


Fig. 7 Representative MS/MS Chromatograms of Food A Film

Table 3 Results of Quantitative Analysis of Polymer Additives in Food Containers

| Compound name | Concentration (mg/g) | | | | | | | | | |
|------------------|----------------------|-------------|-------------|-------------|-------------|-------------|-------------|-------------|-------------|-------------|
| | Food A pack | Food A film | Food B pack | Food B film | Food C pack | Food C film | Food D pack | Food D film | Food E pack | Food E film |
| Tinuvin® P | ---- | ---- | ---- | ---- | ---- | ---- | ---- | ---- | ---- | ---- |
| Irganox® 245 | ---- | ---- | ---- | 0.043 | ---- | ---- | ---- | ---- | ---- | ---- |
| Irganox® MD 1024 | 0.823 | 0.695 | 0.627 | 0.486 | 0.479 | 0.430 | 0.376 | 0.400 | 0.318 | 0.278 |
| Irganox® 1098 | ---- | 7.104 | ---- | 8.64 | ---- | ---- | ---- | ---- | ---- | ---- |
| CYANOX® 2246 | ---- | 0.021 | ---- | ---- | ---- | ---- | ---- | ---- | ---- | ---- |
| CYANOX® 425 | 0.130 | 3.132 | ---- | 0.069 | ---- | ---- | ---- | ---- | ---- | ---- |
| Irganox® 1035 | 0.011 | 0.012 | ---- | ---- | ---- | ---- | ---- | ---- | ---- | ---- |
| Tinuvin® 120 | ---- | 0.005 | ---- | ---- | ---- | ---- | ---- | ---- | ---- | ---- |
| Tinuvin® 328 | 0.024 | ---- | ---- | ---- | ---- | ---- | 0.268 | ---- | ---- | ---- |
| Irganox® 1010 | 9.544 | 51.094 | 1.698 | ---- | 14.054 | 76.426 | 6.260 | 58.466 | 15.218 | 113.920 |
| Irganox® 1330 | ---- | ---- | ---- | ---- | ---- | ---- | ---- | ---- | ---- | 0.004 |
| Irganox® 565 | ---- | 0.159 | ---- | ---- | ---- | ---- | ---- | ---- | ---- | 0.135 |
| Irganox® 1076 | 2.140 | 8.366 | ---- | 25.450 | 2.636 | 2.482 | 7.994 | 9.644 | 1.484 | 8.438 |
| Irgafos® 168 | 111.04 | 339.94 | 119.64 | 37.1 | 253.68 | 799.66 | 350.10 | 616.62 | 205.24 | 126.96 |

Conclusion

A workflow that includes the processes from detection to qualitative analysis and quantitative analysis of the functional additives contained in polymer materials was realized by using a Shimadzu LCMS-9030 quadrupole time-of-flight liquid chromatograph mass spectrometer and LabSolutions Insight Explore and ACD/MS Structure ID Suite.

This technique is expected to contribute to efficient development and improvement of better synthetic polymer materials.

LCMS and LabSolutions Insight Explore are trademarks of Shimadzu Corporation or its affiliated companies in Japan and/or other countries. Tinuvin, Irganox, and Irgafos are trademarks of BASF SE in Japan and/or other countries. CYANOX is a trademark of Cytec Technology Corp. in the United States and/or other countries. Kinetex is a trademark of Phenomenex Inc. in the United States and/or other countries.



Shimadzu Corporation
Analytical & Measuring Instruments Division
Global Application Development Center

www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country. The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. See <http://www.shimadzu.com/about/trademarks/index.html> for details.

Third party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®".

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

First Edition: Jul. 2021