

Application News

DPiMS™ QT Probe Electrospray Ionization Kit
LCMS™-9030 Quadrupole Time-of-Flight Liquid Chromatograph Mass Spectrometer

Detection of Isoflavones and Its Metabolites in Foods by DPiMS QT and LCMS-9030

E. Imoto, T. Murata

User Benefits

- ◆ Detection of target compounds in processed foods is possible with simple sample preparation.
- ◆ Compounds can be measured quickly without column separation regardless of its polarity.
- ◆ Samples can be measured with stable mass accuracy.

■ Introduction

With the heightened health-consciousness of recent years, soybean isoflavones have attracted attention as functional components which are useful for health. However, with the global distribution of food products, it is necessary to verify that the components actually contained in food products are the same as the indicated contents.

Here, the measurement targets were the soy isoflavones daidzein, genistein, and glycitein and their glycosides. An ODS column is generally used in simultaneous analyses of isoflavones and their glycosides by LC or LC/MS. It takes excessive time to separate these compounds completely, because isoflavones and their glycosides have different polarities. In addition, since this technique also includes a process for extraction of the isoflavones in foods, the development of an easy and fast analysis technique for all of the steps from sample preparation to measurement has been needed.

This Application News introduces a new analysis technique using a DPiMS QT which is a probe kit of electrospray ionization unit and an LCMS-9030 quadrupole time-of-flight mass spectrometer (Fig. 1). The DPiMS QT makes it possible to conduct direct analysis and minimize the time from sample preparation to analysis.



Fig. 1 Appearance of DPiMS™ QT and LCMS™-9030

■ Analytical Conditions and Sample Preparation Method

In the DPiMS QT, the attached probe repeatedly carries out sampling from the sample plate, and simultaneously ionizes the sample adhering to the probe surface by applying a voltage to the probe tip. Then, the ionized samples are introduced directly into the mass spectrometer. Table 1 shows the analytical conditions.


In this Application News, the daidzein, daidzin, genistein, genistin, and glycitein in tofu, soybeans boiled in water, soy milk, and a soybean processed chocolate snack were measured. Fig. 2 shows the structural formulae of these compounds.

The sample preparation procedure is as follows. First, 10 mg each of the tofu, boiled soybeans, soy milk, and soybean processed chocolate snack were weighed, and 1 mL of 50 % ethanol was added and mixed for 1 min. The mixtures were then centrifuged, and 10 μ L of the supernatant was dripped on the sample plate and measured.


Table 1 Analytical Conditions

Mass spectrometer


System	: DPiMS QT+LCMS-9030
Polarity	: Positive
DL temp	: 250 °C
Heat block temp	: 50 °C
Interface Voltage	: 3.5 kV
TOF-MS	: m/z 100-800
Measurement Time	: 0.5 min

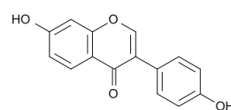


Daidzein
M.W. 254.0579

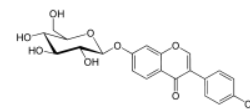


Daidzin
M.W. 416.1107

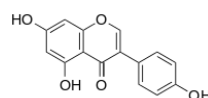

Genistein
M.W. 270.0528

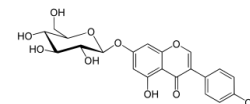

Genistin
M.W. 432.1056

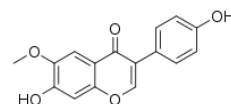

Glycitein
M.W. 284.0685

Fig. 2 Structural Formulae of Isoflavones and Their Glycosides

■ Analysis of Foods Containing Isoflavones

Fig. 3 show the MS spectra obtained from extracts of the tofu, boiled soybeans, soy milk, and soybean processed chocolate snack. With a simple sample preparation time of about 5 min and measurement time of 0.5 min, daidzein, daidzin, genistein, genistin, and glycitein were detected from the tofu, boiled soybeans, and soy milk, and daidzein, genistein, and glycitein were detected from the soybean processed chocolate snack. The mass error (ppm) calculated from the theoretical and measured m/z of each isoflavone was in the range of -2.4 ppm to 2.2 ppm, indicating that stable mass accuracy was realized, including in the matrix.

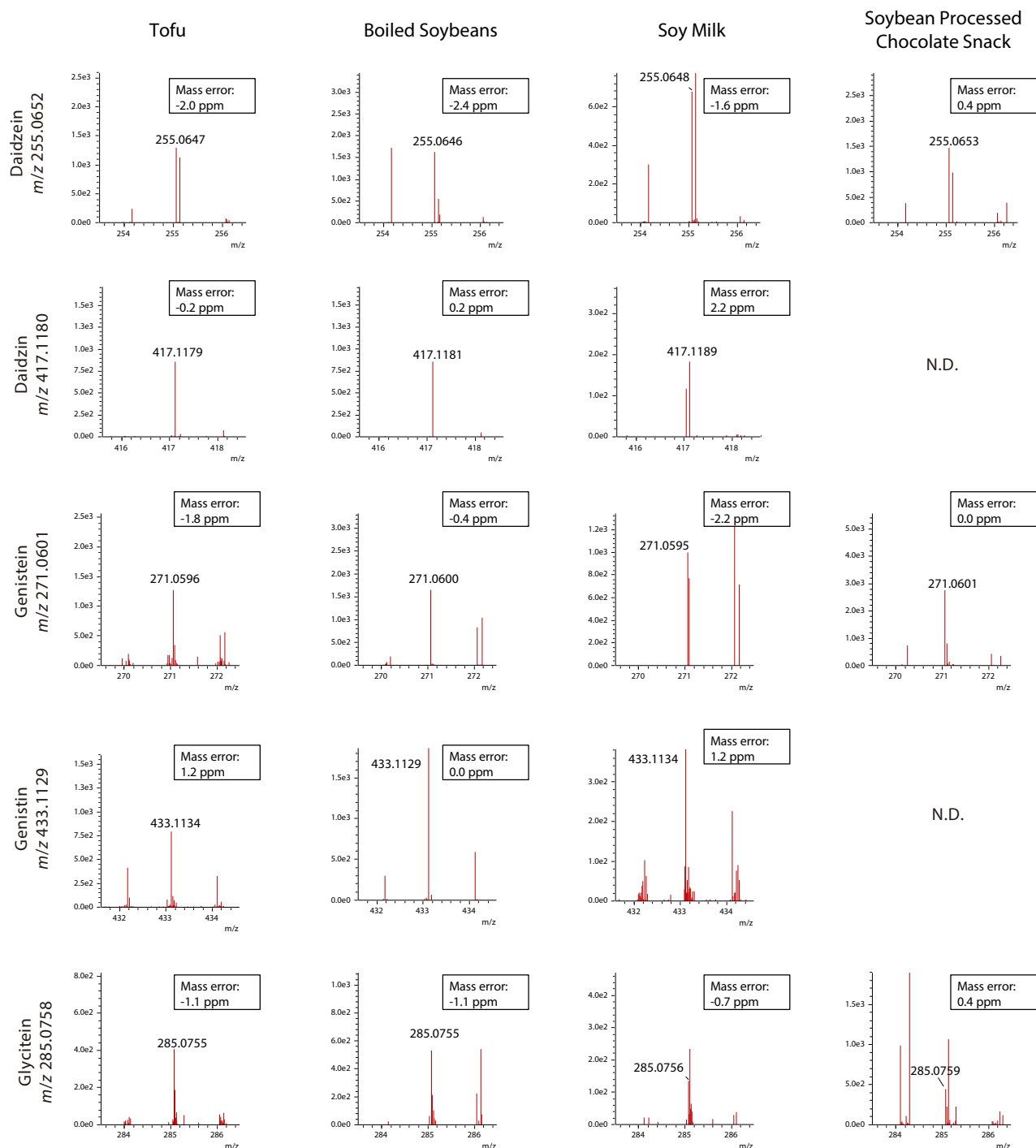


Fig. 3 MS Spectra of Extracts of Tofu, Boiled Soybeans, Soy Milk, and Soybean Processed Chocolate Snack and Mass Error (ppm) from Theoretical and Measured m/z

Conclusion

Isoflavones and their glycosides were detected in tofu, boiled soybeans, soy milk, and a soybean processed chocolate snack with stable mass accuracy by using a combination of a DPiMS QT and LCMS-9030.

The time required for the analysis was substantially reduced in comparison with the LC or LC/MS methods. The sample preparation time was approximately 5 min and the measurement time was 0.5 min.

LCMS and DPiMS are trademarks of Shimadzu Corporation or its affiliated companies in Japan and/or other countries.

For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.



Shimadzu Corporation
Analytical & Measuring Instruments Division
Global Application Development Center

www.shimadzu.com/an/

01-00169-EN

First Edition: Sep. 2021

The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. See <http://www.shimadzu.com/about/trademarks/index.html> for details.

Third party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®".

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

➤ Please fill out the survey

Related Products

Some products may be updated to newer models.



➤ LCMS-9030

Quadrupole Time-of-Flight Liquid
Chromatograph Mas...



➤ DPiMST™ QT

Kit for Direct Probe Ionization Mass
Spectrometer



➤ LCMS-9050

Quadrupole Time-of-Flight Liquid
Chromatograph Mas...

Related Solutions

➤ Food and Nutrition

➤ Price Inquiry

➤ Product Inquiry

➤ Technical Service /
Support Inquiry

➤ Other Inquiry