

Application Data Sheet

No. 4

Autograph Precision Universal Tester

Material Testing & Inspection

Tensile Tests of Rubber Dumb-bell Specimens

Standard No. ISO37: 2005 (JIS K6251: 2010)

Introduction

Rubber materials have characteristic mechanical properties including elasticity and flexibility, and are widely used for industrial parts, construction materials, and housewares. In particular, a diverse range of synthetic rubber materials with differing properties suited to match their application have been developed. Measuring these mechanical properties is extremely important to ensure quality control and for new materials development. This Data Sheet introduces an example of the evaluation of three synthetic rubber (main components: chloroprene [1]; urethane [2]) specimens (dumb-bell test pieces). Static tensile tests were performed, and the specimens were evaluated with respect to tensile strength, stress at given elongation, and elongation at break, which are aspects of their basic mechanical properties.

T. Murakami

Measurements and Jigs

In tensile tests of rubber dumb-bell specimens, the grips must tighten automatically. When tensile loads are applied to rubber materials, they elongate and their thickness decreases. For this reason, if there is no automatic tightening mechanism, the specimen will inadvertently break free of the grip before the maximum load is applied, making favorable measurements impossible. Accordingly, in rubber tensile tests, it is necessary to use pneumatic parallel grippers, pantograph grips, eccentric roller type grips, Henry Scott type grips, or other grips equipped with this feature.

Measurement Results



Fig. 1: Test Status

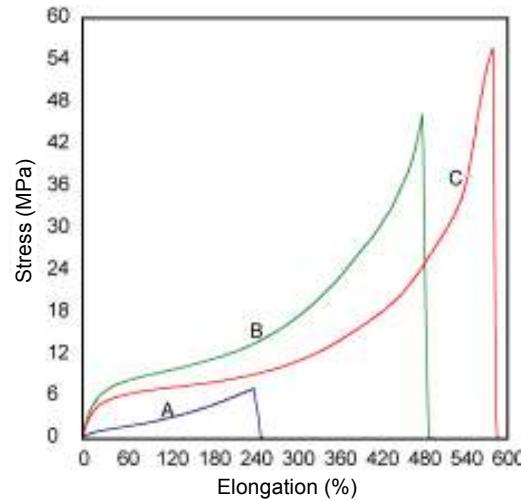


Fig. 2: Relationship Between Stress and Elongation

Table 1: Test Conditions

Item	Set Value
Test Speed	500 mm/min
Initial Distance between Grips	60 mm
Gauge Length	20 mm

Table 2: Test Results

Sample	Main Component	Tensile Strength (MPa)	Stress at 100 % Elongation (MPa)	Stress at 200 % Elongation (MPa)	Elongation at Break (%)
A	Chloroprene	7.1	2.4	5.4	242.3
B	Urethane	47.1	9.3	11.9	478.1
C	Urethane	55.5	7.0	8.4	572.6

The tensile test results for the three samples are shown in Table 2. A graph showing the stress-elongation relationship for the samples is shown in Fig. 2. Clear differences in mechanical properties such as tensile strength and elongation at break are apparent between the samples.

Rubber (Dumb-bells test pieces) Tensile Test System

Tester: AGS-X
Load Cell: 1 kN
Test Jig: 1 kN pneumatic flat grips (Single-side file teeth grip faces)
Extensometer: SES-1000 type extensometer for soft specimens
Software: TRAPEZIUM LITE X



AGS-X Table-Top Precision Universal Tester

Features

- A high-precision load cell is adopted. (The high-precision type is class 0.5; the standard-precision type is class 1.) Accuracy is guaranteed over a wide range, from 1/500 to 1/1 of the load cell capacity. This supports highly reliable test evaluations.
- Crosshead speed range
Tests can be performed over a wide range from 0.001 mm/min to 1,000 mm/min.
- High-speed sampling
High-speed sampling, as fast as 1 msec.
- TRAPEZIUMX LITE X operational software
This is simple, highly effective software.
- Jog controller (optional)
This allows hand-held control of the crosshead position. Fine position adjustment is possible using the jog dial.
- Optional Test Devices
A variety of tests can be conducted by switching between an abundance of jigs in the lineup.

First Edition: February 2013



Shimadzu Corporation

www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedures.
The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

© Shimadzu Corporation, 2013